

Paradigm[®] Veo[™]

User Guide

©2008 Medtronic MiniMed, Inc. All rights reserved.

Bolus Wizard®, CareLink®, Dual Wave®, Guardian®, MiniLink®, Paradigm®, Quick-serter®, Quick-set®, Silhouette®, Sof-set®, and Square Wave® are registered trademarks of Medtronic MiniMed, Inc.

Easy Bolus™, Enlite™ and Veo™ are trademarks of Medtronic MiniMed, Inc.

Energizer® is a registered trademark of Eveready Battery Company.

Glucagon Emergency Kit® is a registered trademark of Eli Lilly and Company.



Medtronic

Medtronic MiniMed
Northridge, CA 91325
USA
800 646 4633
818 576 5555

EC	REP
----	-----

Medtronic B.V.
Earl Bakkenstraat 10
6422 PJ Heerlen
The Netherlands
31 (0) 45 566 8000
www.medtronicdiabetes.com



6025644-021_a

REF MMT-554 MMT-754

Contacts:**Africa:**

Medtronic Africa (Pty) Ltd.
Tel: +27 (0) 11 677 4800

Australia:

Medtronic Australasia Pty. Ltd.
Tel: 1800 668 670 (product orders)
Tel: 1800 777 808 (customer help)

Bangladesh

Sonargaon Healthcare Pvt Ltd
Mobile: (+91)-9903995417
or (+880)-1714217131

België/Belgique:

N.V. Medtronic Belgium S.A.
Tel: 0800-90805

Brasil:

Medtronic Comercial Ltda.
Tel: +(11) 3707-3707

Canada:

Medtronic of Canada Ltd.
Tel: 1-800-284-4416 (toll free/sans-frais)

China:

Medtronic (Shanghai) Ltd.
Tel: +86 40 0820 1981
or 80 0820 1981

Croatia

Medtronic Adriatic
Tel: +385 1 488 11 20
Fax: +385 1 484 40 60

Danmark:

Medtronic Danmark A/S
Tel: +45 32 48 18 00

Deutschland:

Medtronic GmbH
Geschäftsbereich Diabetes
Telefon: +49 2159 8149-370
Telefax: +49 2159 8149-110
24-Std-Hotline: 0800 6464633

Eire:

Accu-Science LTD.
Tel: +353 45 433000

España:

Medtronic Ibérica S.A.
Tel: +34 91 625 05 42
Fax: +34 91 625 03 90
24 horas: +34 901 120 335

Europe:

Medtronic Europe S.A. Europe, Middle East and Africa Headquarters
Tel: +41 (0) 21-802-7000

France:

Medtronic France S.A.S.
Tel: +33 (0) 1 55 38 17 00

Hellas:

Medtronic Hellas S.A.
Tel: +30 210677-9099

Hong Kong:

Celki International Ltd.
Tel: +852 2332-3366

India:

India Medtronic Pvt. Ltd
Tel: (+91)-80-22112245 / 32972359
Mobile: (+91)-9611633007

Italia:

Medtronic Italia S.p.A.
Tel: +39 02 24137 261
Fax: +39 02 24138 210
Servizio assistenza tecnica: N° verde
24h: 800 712 712

Japan:

Medtronic Japan Co. Ltd.
Tel: +81-3-6430-2019

Latin America:

Medtronic, Inc.
Tel: 1(305) 500-9328
Fax: 1(786) 709-4244

Latvija:

Ravemma Ltd.
Tel: +371 7273780

Magyarország:

Medtronic Hungária Kft.
Tel: +36 1 889 0697

Malaysia:

Medtronic International Ltd.
Tel: +60-3 7953 4800

Middle East and North Africa:

Regional Office
Tel: +961-1-370 670

Nederland, Luxembourg:

Medtronic B.V.
Tel: +31 (0) 45-566-8291
Toll Free: 0800-3422338

New Zealand:

Medica Pacifica
Tel: +0800 106 100
24-hr After-Hours:
+0800 633 487

Norge:

Medtronic Norge A/S
Tel: +47 67 10 32 00
Fax: +47 67 10 32 10

РОССИЯ:

Medtronic B. V.
Tel: +7 495 580 73 77
24h: 8-800-200-76-36

Polska:

Medtronic Poland Sp. Z.o.o.
Tel: +48 22 465 6934

Portugal:

Medtronic Portugal Lda
Tel: +351 21 7245100
Fax: +351 21 7245199

Puerto Rico:

Medtronic Puerto Rico
Tel: 787-753-5270

Republic of Korea:

Medtronic Korea, Co., Ltd.
Tel: +82.2.3404.3600

Schweiz:

Medtronic (Schweiz) AG
Tel: +41 (0)31 868 0160
24-Stunden-Hotline: 0800 633333
Fax Allgemein: +41 (0)318680199

Serbia & Montenegro:

Epsilon
Tel: +381 11 311 8883

Singapore:

Medtronic International Ltd.
Tel: +65 6436-5097
or +65 6436-5090

Slovenija:

Zaloker & Zaloker d.o.o.
Tel: +386 1 542 51 11
Fax: +386 1 542 43 32

Slovenská republika:

Medtronic Slovakia O.Z.
Tel: +421 268 206 944
Fax: +421 268 206 999

Sri Lanka

Swiss Biogenics Ltd.
Mobile: (+91)-9003077499
or (+94)-777256760

Suomi:

Medtronic Finland Oy
Tel: +358 20 7281 232
Help line: +358 20 7281 200

Sverige:

Medtronic AB
Tel: +46 8 568 585 10
Fax: +46 8 568 585 11

Taiwan:

Medtronic-Taiwan Ltd.
Tel: +886.2.2183.6093
Fax: +886.2.2501.7680
Toll Free: 0800.082.999

Thailand:

Medtronic (Thailand) Ltd.
Tel: +66 2 232 7400 ext 1

Turkiye:

Medtronic Medikal Teknoloji
Ticaret Ltd. Sirketi.
Tel: +90 216 4694330

USA:

Medtronic Diabetes Global
Headquarters
Tel: +1-800-826-2099
24 Hour HelpLine: +1-818-576-5555
To order supplies: +1-800-843-6687

United Kingdom:

Medtronic Ltd.
Tel: +44 1923-205167

Yisrael:

Agentek
Tel: +972 3649 3111

Österreich:

Medtronic Österreich GmbH
Tel: +43 (0) 1 240 44-0
24 - Stunden - Hotline: 0820 820 190

Česká republika:

Medtronic Czechia s.r.o.
Tel: (+420-233) 059401

Important Safety Information Regarding Your Paradigm Insulin Pump (includes all models)

Avoid Immersing Your Pump In Water

Although it is unlikely that water damage will occur if your pump is splashed or briefly dunked, **you should avoid immersing your Paradigm insulin pump in water. To swim or participate in other water activities, always disconnect from your Paradigm pump and reconnect after you are out of the water.**

If you inadvertently submerge your pump in water, dry the pump quickly using a soft, clean towel and verify that it is working properly by selecting **Selftest** from the pump's Utilities Menu. If you believe that water has entered your pump or you observe any other possible pump malfunction, please check your blood glucose, treat high blood glucose (if necessary) with an injection and contact your local help line or representative for further assistance. Symptoms of high blood glucose include fatigue, excessive thirst and nausea. You should always contact your healthcare professional if you experience excessively high or low blood glucose levels, or if you have any questions about your care.

Electrostatic Discharge

Although your Paradigm pump is designed to be unaffected by typical levels of electrostatic discharge (ESD), very high levels of ESD can result in a reset of the pump's software with an associated pump error alarm. In most cases, exposure to high levels of ESD will trigger the pump's A-13 alarm although, under certain circumstances, high level ESD exposure can cause A-44, Bolus Stopped or Max Delivery alarms. High levels of ESD are more likely in situations where the relative humidity is very low, such as inside a heated building during the winter in areas where it is cold outside.

If your pump experiences an A-13 or other error alarm, press the **ESC** and **ACT** buttons to clear the alarm. If you are unable to clear the alarm by pressing **ESC** and **ACT**, you may need to remove and replace the pump's battery to clear the alarm. After clearing the alarm, you should always verify that your pump is set to the correct date and time and that all other settings (basal rate, max basal and bolus limits, etc.) are programmed to the desired values, since the software reset could erase your previously programmed settings. Please see the Alarms and Alerts section of this User Guide for more details regarding what to do if your pump displays an error alarm or other alert message.

Please contact your local helpline or representative to report any error alarms or other problems that occur with your pump.

Warranty

Medtronic Diabetes warrants the Medtronic MiniMed Insulin Pump against defects in materials and workmanship for a period of 4 years from the date of purchase.

During the warranty period, Medtronic Diabetes will, at its discretion, either repair or replace (with a new or recertified pump, at Medtronic Diabetes' discretion) any defective pump or motor, subject to the conditions and exclusions stated herein. In the event that a pump is repaired or replaced, the warranty period will not be extended.

This warranty is valid only if the Medtronic MiniMed Insulin Pump is used in accordance with the manufacturer's instructions. This warranty will not apply:

- If damage results from changes or modifications made to the pump by the user or third persons after the date of manufacture.
- If non-Medtronic infusion sets and/or reservoirs are being used with the pump.
- If damage results from service or repairs performed by any person or entity other than the manufacturer.
- If damage results from a *Force Majeure* or other event beyond the control of the manufacturer.
- If damage results from negligence or improper use, including but not limited to: improper storage, submersion in water or physical abuse, such as dropping or otherwise.

This warranty shall be personal to the original user. Any sale, rental or other transfer or use of the product covered by this warranty to or by a user other than the original user shall cause this warranty to immediately terminate. This warranty does not apply to batteries, infusion sets, reservoirs, and other accessories.

The remedies provided for in this warranty are the exclusive remedies available for any breach hereof. Neither Medtronic Diabetes nor its suppliers or distributors shall be liable for any incidental, consequential, or special damage of any nature or kind caused by or arising out of a defect in the product. Notwithstanding the foregoing, the terms of this warranty do not affect the statutory rights of consumers or seek to exclude any liability which is not permitted by law.

All other warranties, other than mandatory statutory warranties, expressed or implied, are excluded, including the warranties of merchantability and fitness for a particular purpose.

Contents

Chapter 1: Introduction	1
Assistance	1
Emergency kit	1
Consumables	2
Accessories	2
How to wear your pump	3
How to use this guide	3
User safety	4
Indications	4
Contraindications	5
Warnings	5
Pump	5
Reservoir and infusion sets	5
Sensor	5
Transmitter	5
Magnetic fields	6
X-rays, MRIs and CT scans	6
Precautions	6
Avoid extreme temperatures	6
Infusion sets and sites	7
Sensor	7
Adverse reactions	7
Notice	7
Insulin pump and RF accessories	7
RF interference from other devices	8
Chapter 2: Introduction to pump therapy	11
Theory of insulin pump therapy	11

Understanding pump therapy	12
Basal rate	12
Meal bolus	12
Gram counting	12
Exchange counting	12
BG targets	12
Insulin sensitivity	12
Active Insulin	13
Blood glucose and A1C testing	13
BG testing	13
A1C	13
Low blood glucose (hypoglycemia)	14
Hypoglycemia protocol: the rule of 15	15
High blood glucose (hyperglycemia)	16
Hyperglycemia protocol	16
Diabetes ketoacidosis (DKA)	17
Preventing DKA	17
Sick day management	18
Sick day protocol	18
Sick day supplies	18
Nutrition	19
Carbohydrate counting	19
Carbohydrate gram counting	19
Carbohydrate exchange system	19
Reading a food label	20
Fat and blood glucose	21
Protein and blood glucose	21
Other things that affect blood glucose	22
Healthy eating habits	22
Exercise	22
Chapter 3: The basics	25
Your pump	25
Install battery	26
Pump buttons	27
The pump screen	29
HOME screen	29

Screen icons	29
Battery	29
Time display	30
Reservoir volume	30
Alert and alarm icons	30
Sensor icons	30
Scroll bar	31
Screen backlight	31
Beep/vibrate	31
Operating modes	31
Normal mode	32
Special mode	32
Attention mode	32
Menus	33
MAIN MENU	33
BOLUS MENU	33
SUSPEND	33
SENSOR	33
CAPTURE EVENT	34
BASAL MENU	34
RESERVOIR + SET	34
UTILITIES MENU	34
STATUS screen	34
If you remove your pump	35
Chapter 4: Basic programming	37
Setting the time and date	37
Selecting the language	39
Bolus	39
Setting the Normal bolus	40
Normal meal bolus using the exchange system	41
Bolus practice:	41
Normal meal bolus practice using exchanges	42
Normal meal bolus using carbohydrate counting	42
Normal meal bolus practice using carbohydrates	43
Meal bolus, correction bolus and insulin sensitivity	43
Practice: Meal bolus	44

Review your bolus deliveries	44
Bolus details	45
Maximum bolus limit	47
Example 1: Max bolus	47
Example 2: Max bolus	47
Scroll rate	47
Example: Scroll rate	48
BG Reminder	48
Capturing Events	49
Turning on the Capture Event feature	49
Entering BG measurements	50
Entering insulin injection information	51
Entering carbohydrate information	52
Entering exercise information	52
Entering Other markers	53
Viewing Capture Event history	53
Basal	53
Start and stop times	54
Your basal settings	55
Basal programming and delivery	55
Current basal delivery	56
Daily basal rate(s)	56
Setting the Max basal rate	57
Example 1: Max basal	58
Example 2: Max basal	58
Stopping your pump	58
Resume pump delivery	59
Example: Suspend function	60
Practice: Suspend function	60
Practice: Resume basal delivery after a suspend	60
Chapter 5: Starting on insulin	61
Prepare your pump for use	61
Filling the reservoir	61
Filling the reservoir	62
Changing your infusion set	64
Removing the reservoir	64

Rewinding your pump	64
Inserting the reservoir in your pump	65
Filling the tubing	66
Inserting the infusion set	67
Quick-set infusion set (with Quick-serter®)	69
Filling the cannula	69
Disconnecting Quick-set	71
Reconnecting Quick-set	71
Record keeping for diabetes management	72
Determining your pump settings	72
Chapter 6: Using the Bolus Wizard feature	75
What is it?	75
Carbohydrate counting	75
Your blood glucose reading	75
Your personal Bolus Wizard feature settings	75
How the Bolus Wizard feature works	78
Bolus Wizard warnings	79
HIGH BG	79
LOW BG	79
MAX BOLUS EXCEEDED	79
How to program the Bolus Wizard feature	79
Turning on the Bolus Wizard feature	80
Selecting the Carb units	80
Setting the Carb/Exch ratios	81
Setting the BG units	82
Insulin sensitivity	82
Setting the BG Targets	83
About active insulin	84
Active insulin time	85
Review your Bolus Wizard feature settings	85
Meter option	86
Meter rules	86
Add, delete, review meter IDs	87
Normal bolus using Bolus Wizard feature	87
Bolus Wizard feature examples	89
Example 1: Blood glucose on target (normal blood glucose) and no active insulin	89

Example 2: Blood glucose above target (high blood glucose) and no active insulin	90
Example 3: Blood glucose below target (low blood glucose) and no active insulin	91
Example 4: Blood glucose above target (high blood glucose) with active insulin	92
Example 5: Blood glucose below target (low blood glucose) with active insulin	93

Chapter 7: Optimizing pump therapy 95

Square Wave and Dual Wave bolus	95
Turning on the Dual Wave/Square Wave option	96
Square Wave or Dual Wave bolus without Bolus Wizard feature	96
Square Wave bolus practice	98
Dual Wave bolus practice	99
Using the Bolus Wizard feature for a Square Wave or Dual Wave bolus	99
Easy bolus	101
Easy bolus setup	101
Step value setup	102
Delivering Easy Bolus	102
Example 1: Easy bolus	103
Your turn: Easy bolus practice	104
Missed Bolus Reminder	104
Add, delete, and review bolus reminders	104
Add	105
Delete	105
Review	105
Basal patterns	105
Turning on the patterns	106
Programming a pattern	106
Select a pattern	107
Example 1: Basal patterns	107
Example 2: Basal patterns	108
Your turn:	108
Temp basal rates	108
How does temp basal work?	108
Temp basal types	109
Insulin rate (U/H)	109
Percent of basal	110

Chapter 8: Sensor features	117
Entering your sensor settings	117
Turning on the sensor	117
Turning on the Glucose Alerts	118
Setting the Glucose Limits	118
Selecting the BG units	119
Setting the Glucose Limits start time	120
Setting the High Repeat	121
Setting the Low Repeat	122
Setting the Predictive glucose alert	122
Setting the Rate of Change Alerts	124
Setting the Low Glucose Suspend	125
Setting the Cal Repeat	125
Setting the Cal Reminder	126
Turning on the Auto Calibration	126
Entering the transmitter ID	126
Setting up the Weak Signal feature	127
Setting the sensor glucose graph timeout	128
Reviewing your settings	128
Using the Sensor Demo	129
The transmitter	129
Transmitter charger	129
Installing a battery in the gray charger	130
Installing a battery in the blue charger	130
Charging the transmitter	131
Starting the sensor	132
Inserting the Enlite™ sensor	132
Connecting the transmitter to the sensor	136
Preparing the sensor for communication	137
Calibrating the sensor	137
Entering meter BG	138
Chapter 9: Using your sensor	141
Status screens	141
Reading the sensor glucose graphs	142
Opening and viewing the graphs	143
The graphs	144

Examples of continuous sensor glucose graphs	145
3-hour graph	145
6-hour graph	146
12-hour graph	146
24-hour graph	147
How to check for rapid changes in sensor glucose	148
The rapid change arrows	148
How to make glucose alerts silent	149
Calibration history	150
Sensor alert history	150
Disconnecting the transmitter and removing the sensor	151
Disconnecting the sensor from the transmitter	151
Removing the sensor	151
Storage	151
Using your system in water	151

Chapter 10: Utilities 153

Alarm review	153
Setting your alert type	153
Auto-off	154
LOW RESV WARNING (Low reservoir warning)	154
Review daily insulin totals	155
Pump data management	156
Insulin delivery information	156
Sensor glucose information	157
Meter blood glucose information	158
AUC calculation	160
Personal reminders	162
Alarm clock	162
Remote control option	162
Turn on remote control option	163
Add, delete, review remote control IDs	163
Other Devices option	163
Block feature	164
Turning block on	164
Example 1: Block	164
Example 2: Block	164

Lock keypad feature	165
Locking the keypad	165
Unlocking the keypad	165
Selftest	165
User settings	166
Saving the settings	167
Restoring the settings	167
Clearing the settings	168
History	168
Chapter 11: Therapy management software	169
CareLink® Personal Software	169
Reports	169
Logbook	170
Chapter 12: Insulin pump therapy follow-up	171
Recommended follow-up	171
Everyday	171
Every month	171
Every three months	171
Laboratory tests	172
Every visit	172
Annually	172
Chapter 13: Troubleshooting, alerts and alarms	173
My pump has a NO DELIVERY alarm	173
What happens if I leave the pump battery out too long?	174
Why doesn't my pump battery last very long?	174
What is a CHECK SETTINGS alarm?	175
My screen appears distorted	175
The pump is asking me to rewind	175
My bolus stopped	175
My pump buttons are not acting right during a bolus	176
My pump will not display my blood glucose reading from my meter	176
I dropped my pump	176
I submerged my pump in water	177
I cannot get to the USER SETTINGS screen	177

Alerts	177
Example:	178
What to do	178
Pump alert conditions	179
LOW RESERVOIR	179
LOW BATTERY	179
Sensor alert conditions	179
WEAK SIGNAL	179
ALERT SILENCE	179
MISSED BOLUS	180
LOST SENSOR	180
LOW TRANSMTR	180
CHARGE TRANSMTR	181
CHANGE SENSOR	181
SENSOR END	181
CAL ERROR	181
METER BG NOW	182
METER BG BY	182
LOW X.X MMOL/L (XX MG/DL)	182
HIGH XX.X MMOL/L (XXX MG/DL)	182
LOW PREDICTED	182
HIGH PREDICTED	183
RISE RATE	183
FALL RATE	183
SENSOR ERROR	183
Alarms	183
What to do	184
Alarm conditions	185
A (ALARM)	185
AUTO OFF	185
BATT OUT LIMIT	185
BOLUS STOPPED	185
BUTTON ERROR	185
CHECK SETTINGS	186
E (ERROR)	186
EMPTY RESERVOIR	186
FAILED BATT TEST	186

FINISH LOADING	186
MAX FILL REACHED	186
MAX DELIVERY	187
MOTOR ERROR	187
NO DELIVERY	187
NO RESERVOIR	187
OFF NO POWER	187
RESET	187
WEAK BATTERY	187
LOW SUSPEND	188
Testing the transmitter	189
Connecting the tester	189
Disconnecting the tester	190
Troubleshooting sensor features	191
Reconnect old sensor	191
Find lost sensor	191
Understanding your transmitter, tester, and charger	191
Chapter 14: Maintenance	193
Battery	193
Storage	194
Cleaning your pump	194
Cleaning the transmitter	194
Cleaning the Serter	195
Chapter 15: Pump specifications	197
Alarms and error messages	197
Alarm history	197
Audio frequency	197
Backlight	198
Basal	198
BG Target	198
Bolus delivery	198
Maximum bolus	198
Bolus history	198
Bolus units	199
Bolus Wizard feature	199

Carb ratios	199
Carb units	199
Daily totals	199
Default screen	199
Delivery accuracy	199
Pump motor	200
Dual Wave bolus	200
Easy bolus	201
Filling the infusion set	201
Infusion pressure	201
Insulin sensitivity	201
Low resv (reservoir) warning	202
Meter value	202
Normal bolus	202
Occlusion detection	202
Percent temp basal	202
Power supply	203
Reservoir + Set history	203
Program safety checks	203
Pump size	203
Pump weight	203
Remote control	203
Reservoir	203
Square Wave bolus	204
Temporary (temp) basal rate	204
Time and date screen	204
Environmental conditions	204
Status screen	204
Bolus Wizard feature specifications	207
Bolus Wizard feature examples	210
Insulin delivery default settings	213
Bolus Wizard feature default settings	215
Sensor features default settings	215
Guidance and manufacturer's declaration	219
Icon table	225

Glossary 227

Index 235

Introduction

Thank you for choosing Medtronic Diabetes (legally known as Medtronic MiniMed, Inc.) as your partner in helping you gain better control of your diabetes. Whether you are beginning pump therapy for the first time or upgrading from a previous model, we believe that the combination of state-of-the-art technology and the simple, menu-driven programming of the pump will provide many benefits.

This user guide is designed to help you to understand pump therapy and the operation of your MiniMed Paradigm® Veo™ system. We strongly recommend that you work closely with your healthcare professional for a safe and complete pump start.

Assistance

Please contact your local help line or representative for assistance. Refer to the Medtronic Diabetes International Contacts list at the beginning of this user guide for the help line or representative in your area.

Emergency kit

Keep an emergency kit with you at all times to make sure that you always have necessary supplies. Inform a family member, co-worker, and/or friend where this emergency kit is kept. Please refer to the *User safety* section in this chapter for more information on pump safety. Your emergency kit should include these items:

- Fast-acting glucose tablets
- Blood glucose monitoring supplies
- Urine ketone monitoring supplies
- Extra Paradigm® compatible infusion set and Paradigm reservoir
- Insulin syringe and fast-acting insulin (with dosage instructions from your healthcare professional)
- Wallet card
- Dressing and adhesive
- Glucagon Emergency Kit®
- Extra AAA alkaline batteries (Energizer® brand is recommended)

CAUTION: If you give yourself insulin by using a syringe, the Bolus Wizard feature will not be able to correctly determine the active insulin in your system. Consult with your healthcare professional on how long you need to wait after a manual injection before you can rely on the active insulin calculation of your Bolus Wizard feature.

Consumables

The pumps use disposable reservoirs and infusion sets for insulin delivery. Installation instructions for Paradigm reservoir and Paradigm compatible infusion sets are provided in the *Starting on insulin* chapter.

- **Reservoirs** – The 554 pump is only used with the 176-unit Paradigm reservoir (MMT-326A). The 754 pump can be used with either the 300-unit Paradigm reservoir (MMT-332A) or the 176-unit reservoir, depending on your insulin needs.
- **Infusion sets** – Medtronic Diabetes provides a variety of Paradigm-compatible infusion sets to fit your needs. Contact your healthcare professional for help in choosing an infusion set that fits your needs. Change your infusion set every two to three days.

WARNING: For your protection the pump has undergone extensive testing to confirm appropriate operation when used with Paradigm reservoirs and Paradigm compatible infusion sets manufactured or distributed by Medtronic Diabetes. We recommend using Medtronic Diabetes infusion sets and reservoirs as we cannot guarantee appropriate operation if the pump is used with reservoirs or infusion sets offered by third-parties and therefore we are not responsible for any injury or malfunctioning of the pump that may occur in association with such use.

Accessories

- **Meter** – Your pump can be used with an optional blood glucose meter powered by MWT1 technology (where or if available). MWT1 is the wireless Radio Frequency (RF) technology that is used to transmit information from the meter to the pump. You can program your pump to automatically receive your BG reading from this meter. All meters referenced in this user guide are blood glucose meters supported by MWT1 technology.
- **Remote control** – The optional Paradigm remote control can be used with the pump to deliver normal boluses and suspend/resume the pump from a distant location. (This user guide provides programming instructions for the remote control. Refer to the remote control user guide for operating instructions.)
- **Medtronic MiniLink® Transmitter** – The transmitter (MMT-7703) is a small device that connects to the sensor. It comes with a tester and a charger. When connected to a sensor that is inserted in the body, the transmitter automatically initializes the sensor and begins to periodically send glucose data to the pump using a radio signal.

- **Sensor** — The sensor (MMT-7002/MMT-7003/MMT-7008) continuously converts tiny amounts of glucose from your fatty layer under the skin into an electronic signal. This signal is sent to the transmitter.
- **CareLink USB Upload Device** — The Medtronic Diabetes CareLink USB (MMT-7305) is used to download the Paradigm 554 or 754 pump data to the diabetes management software using a USB port on your computer.

Not all devices or accessories are available in all countries where the pump is approved.

To order supplies, contact your local representative. Refer to the enclosed contact card for the Medtronic Diabetes office in your area.

How to wear your pump



There are different ways to wear your pump. Medtronic Diabetes has optional accessories that can hide, protect, and add to the convenience of wearing a pump. Refer to the accessories catalog for more information.



- **Holster** — To wear the pump on your belt.
- **Pump clip** — To wear the pump underneath your clothing.
- **Activity guard** — If you are active in sports, or you are a child, use the guard to protect the pump from disconnecting.
- **Leather case** — Fine leather lined with nylon. Styling complements business and formal wear. Velcro flap provides easy access for programming. Wear it vertically with the built-in belt clip.

How to use this guide

NOTE: *This user guide shows sample screens only. Your pump screens may be slightly different.*

For step-by-step instructions, refer to the appropriate sections in this guide. Refer to the Glossary for definitions of terms and functions. The terms and symbols used in this guide are in the table below.

Conventions	Meaning
Press	to push and release the button
Hold	to push and keep pressure on the button
Select	to press  or  to highlight a screen item you want to select
Exit the menus	press ESC until the HOME screen appears

Conventions	Meaning
Pump buttons	always bold and uppercase; for example, ESC , ACT
Screen and menu names	always uppercase; for example, MAIN MENU, REWIND screen
Menu selections	always bold; for example, 24-Hour Setup , On , Off
Flashing (blinking) screen item	you can change the value for that item with the  or  buttons
NOTE and TIP	additional helpful information
CAUTION	warns of a potential hazard which, if not avoided, may result in minor or moderate injury to the equipment
WARNING	notifies you of a potential hazard which, if not avoided, could result in death or serious injury. It may also describe potential serious adverse reactions and safety hazards
Go to the...screen.	<p>when a step instructs you to go to a certain screen, the path to that screen is shown. For example: Go to the ALARM MENU.</p> <p>Main > Utilities > Alarm</p> <ol style="list-style-type: none"> 1 From the MAIN MENU, select Utilities and press ACT. 2 In the UTILITIES MENU, select Alarm and press ACT. 3 The ALARM MENU appears.

User safety

Indications

The Paradigm 554/754 pump systems are indicated for the continuous delivery of insulin, at set and variable rates, for the management of diabetes mellitus in persons requiring insulin. In addition, the pump system is indicated for continuous or periodic monitoring of glucose levels in the fluid under the skin, and possible low and high blood glucose episodes. The pump displays continuous glucose values and stores this data so that it can be analyzed to track patterns and improve diabetes management. Pump history can be downloaded to a computer for analysis of historical glucose values.

The continuous glucose values provided by the Paradigm 554/754 pump systems are not intended to be used directly for making therapy adjustments. Rather, they provide an indication that a confirmation fingerstick measurement may be required. All therapy adjustments should be based on measurements obtained using a home glucose monitor and not based on the value displayed by the pump.

Contraindications

Pump therapy is not recommended for people who are unwilling or unable to perform a minimum of four blood glucose tests per day and to maintain contact with their healthcare professional. Successful insulin pump therapy requires sufficient vision or hearing to allow recognition of the pump signals and alarms.

Warnings

Pump

The pump is not suitable for use in the presence of a flammable anaesthetic mixture with air, oxygen or nitrous oxide.

Reservoir and infusion sets

Standard Luer sets are not compatible with the Medtronic MiniMed Paradigm pump. Medtronic Diabetes Paradigm reservoirs and Paradigm-compatible infusion sets are specifically designed for use with the pump. Do not modify your Paradigm reservoir or Paradigm-compatible infusion set.

Do not put any other drugs/medications inside your reservoir to use with this pump. Only insulin that has been prescribed by your physician can be used in this pump.

Sensor

The sensor may create special needs regarding your medical conditions or medications. Please discuss these conditions and medications with your doctor before using the sensor.

Bleeding, swelling, irritation and/or infection at the insertion site are possible risks associated with inserting the sensor and sometimes result from improper insertion and maintenance of insertion site.

Transmitter

Product contains small parts and may pose a choking hazard for young children.

Optional occlusive dressing should be removed if irritation or reaction to this develops.

The transmitter should be disconnected from the sensor while traveling on an aircraft, or if it interferes with another transmitting device.

Magnetic fields

Do **not** use pump cases that have a magnetic clasp.

Do **not** expose your insulin pump to MRI equipment or other devices that generate very strong magnetic fields. The magnetic fields in the immediate vicinity of these devices can damage the part of the pump's motor that regulates insulin delivery, possibly resulting in over-delivery and severe hypoglycemia.

Your pump must be removed and kept outside the room during magnetic resonance imaging (MRI) procedures.

If your pump is inadvertently exposed to a strong magnetic field, discontinue use and contact your local help line or representative for further assistance.

X-rays, MRIs and CT scans

If you are going to have an X-ray, CT scan, MRI or other type of exposure to radiation, **take off your pump, sensor, transmitter, meter and remote control** before entering a room containing any of these equipment.

The Paradigm pump is designed to withstand common electromagnetic interference, including airport security systems. Be sure to carry the Emergency Card provided, when you are traveling. This card provides airport information.

Precautions

Although the pump has multiple safety alarms, it cannot notify you if the set is leaking or the insulin has lost its potency. It is essential, therefore, that you test your blood glucose levels at least four times per day. If your blood glucose is out of range, check the pump and the infusion set to ensure that the necessary amount of insulin is being delivered.

Avoid extreme temperatures

- 1 Avoid exposure of your pump and remote control to temperatures above 42°C (108°F) or below 1°C (34°F).
- 2 Insulin solutions freeze near 0°C (32°F) and degrade at high temperatures. If you are outside in cold weather, wear your pump close to your body and cover it with warm clothing. If you are in a warm environment, take measures to keep your pump and insulin cool.
- 3 Do not steam, sterilize or autoclave your pump or remote control.

Infusion sets and sites

Avoid using an infusion set insertion site that will be irritated by clothing and accessories, or by rigorous stretching and exercise.

Sensor

Prior to exercising, make sure the sensor is firmly attached.

Adverse reactions

Operation of the sensor feature requires the insertion of a glucose sensor into the skin. Bleeding, swelling, bruising, or infection at the sensor insertion site are possible risks of sensor use. The sensor should be removed if redness, pain, tenderness or swelling develop at the insertion site. The optional occlusive dressing should be removed if irritation or a reaction to this develops. Contact your doctor and your local representative in the event of any adverse reaction.

Notice

CAUTION: Any changes or modifications to the devices not expressly approved by Medtronic Diabetes could void your ability to operate the equipment.

Insulin pump and RF accessories

The pump, meter, transmitter and remote control comply with the United States Federal Communications Commission and international standards for electromagnetic compatibility.

Do not use the RF meter to send your blood glucose reading to the pump while on an aircraft. Manually enter your blood glucose.

These devices comply with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesirable operation. It does not interfere with any RF signals transmitted from outside sources.

These standards are designed to provide reasonable protection against excessive radio frequency interference and prevent undesirable operation of the device from unwanted electromagnetic interference. Operation is subject to the following two conditions:

- 1 This device has been tested and found to comply with the regulations governing such devices in your area. For the specific regulation and test results for your area, please contact your local representative.
- 2 This device generates, uses, and can radiate radio frequency energy and, if installed and used in accordance with the instruction, may cause harmful interference to radio communications. If the device does cause interference to radio or television reception, you are encouraged to try to correct the interference by one or more of the following measures:
 - Reorient or relocate the insulin pump/remote control/transmitter/meter
 - Increase the separation between the insulin pump/remote control/transmitter/meter and the device that is receiving/emitting interference

The meter and the transmitter send information to the pump using radio frequency. If other devices that use radio frequency are in use, such as cell phones, cordless phones and wireless networks, they may prevent communication between the pump and the meter. This interference will not cause any incorrect data to be sent and will not cause any harm to your pump or meter. Moving away from or turning off these other devices may allow communication. Refer to the *Troubleshooting, alerts and alarms* chapter to correct interference problems you may have.

Wireless transmission between the pump and transmitter within the 1.8 meters (6 feet) operating range may be interrupted due to the transmitter orientation and the pump position on the body. Move the pump closer to the transmitter or to another position. If a Lost Sensor alarm has occurred retry:

Main > Sensor > Link to Sensor > Find Lost Sensor

If you have questions, please contact your local Medtronic Diabetes representative.

RF interference from other devices

Common consumer electronic devices that transmit in the same frequency band used by the MMT-7703 MiniLink transmitter may prevent the pump from receiving the glucose information sent by the transmitter. Most cellular (mobile) phones and 900 MHz cordless phones, when transmitting or receiving, may cause significant interruption of transmitter-receiver communication. It is likely that other devices operating in similar frequency ranges will have a similar effect. This interference, however, will not cause any incorrect data to be sent and will not cause any harm to your transmitter.

The pump includes a programmable WEAK SIGNAL alert that notifies you when one or more expected transmissions were not received as expected by the pump. (The pump will also issue a LOST SENSOR alert if communication is interrupted for an extended period of time.)

Communication problems can usually be resolved by:

- ensuring that there is less than 1.8 meters (6 feet) between the receiving device and transmitter, to include RF glucose meters;
- turning off or moving away from other RF transmitting devices; or
- reorienting or relocating the transmitter or receiving device, or both.

Communications problems may also occur due to cellular telephone interference. Testing conducted with several different cellular telephones indicates that, while a glucose value is being transmitted, using a cell phone within 31 cm (12 inches) of receiving devices, transmitters or RF glucose meters can interfere with reception of the transmitted values. In the event of such interference, normal communication can be re-established by:

- turning the cell phone off; or
- keeping the cell phone at least 31 cm (12 inches) away from the receiving device, transmitter or glucose meter when a glucose measurement is being transmitted.

If you have questions, please contact your local Medtronic Diabetes representative.

Introduction to pump therapy

Theory of insulin pump therapy

Insulin pumps deliver insulin closer to the way the human pancreas delivers insulin than any other method of treating diabetes.

All people, with or without diabetes, need background insulin for normal functions of the body without food. They also need a dose of insulin on demand - when food is eaten. People without diabetes can trust that their pancreas will produce this insulin for them. People with diabetes need to take insulin in a way that is most like the way a pancreas produces it.

Most people with diabetes who take more than one daily injection, use long-acting insulin. This takes care of their background insulin needs. They take fast-acting insulin for food. When using an insulin pump, only fast-acting insulin is used. You control when and at what rate your insulin is delivered.

Insulin pump therapy allows you to set a **basal rate**, or background insulin. This is delivered all throughout the day and night for normal body function without food. When you exercise, you can reduce the basal rate so that your blood glucose (BG) does not drop too low. When you are sick or have an infection, you can increase the basal rate. This will keep your BG from going up too high.

Insulin pump therapy allows you to give a **bolus**, or dose of insulin on demand - when you eat. You can also increase or decrease your meal bolus based on the foods you choose to eat. A bolus may also be used to lower an elevated BG. This is called a correction bolus.

The Paradigm pump also contains an optional feature called the Bolus Wizard. This does the math for the required bolus amount based on your personal settings. The Bolus Wizard will use your BG reading, carbohydrate intake, and active insulin when coming up with your bolus amount.

Remember, the pump uses only fast-acting insulin. This means you will not need to follow as rigid a schedule as you did before. There is no long-acting insulin telling you when you should eat or when you will need more insulin.

Sound exciting? It is exciting because not only can you manage your blood glucose more easily, your lifestyle can be easier too.

As long as you are ready to learn how to use the pump, the pump can work for you.

Understanding pump therapy

If you are new to pump therapy, you will need the following data from your healthcare professional. Gather this before going to your pump start appointment. If you are unsure, contact your healthcare professional or pump trainer for guidance.

Basal rate

Basal insulin is needed to maintain your target glucose values when you are not eating. Your pump allows up to 48 basal rates to be programmed. Your healthcare professional may only have you program one or two basal rates when you start pump therapy.

Meal bolus

Deliver a meal bolus before eating carbohydrates. Your insulin to carbohydrate ratio is the amount of insulin required to cover a given number of carbohydrates or exchanges. This ratio gives you an idea of how much of an insulin bolus you should take when eating carbohydrates.

Gram counting

1.0 unit of insulin for _____ grams of carbohydrate

Exchange counting

_____ units of insulin for each carbohydrate exchange

BG targets

The Bolus Wizard will use your BG targets to calculate a correction dose. Keeping your BG within target range is key to living well with diabetes. Your healthcare professional should help you decide these targets.

Insulin sensitivity

Insulin sensitivity is used to decide a correction bolus for a high BG. This is how much your blood glucose will be dropped with 1 unit of insulin.

1 unit of insulin will lower BG _____ mmol/L (mg/dL).

Active Insulin

Active insulin is the bolus insulin that has already been delivered to your body, but has not yet been used. The pump considers your active insulin time setting in determining any active insulin still in your body from prior boluses. This may help prevent hypoglycemia caused by over-correcting for high blood glucose.

Blood glucose and A1C testing

When you check your blood glucose with a meter, you measure blood glucose at the moment you perform the test. This number gives you key feedback for making present and day-to-day adjustments in your diabetes management. The A1C test, done at your doctor's, shows you your average blood glucose over the last 60 to 90 days. Both blood glucose and A1C testing are necessary for good diabetes management.

BG testing

With any insulin therapy, you must monitor your blood glucose four to six times a day. With insulin pump therapy, blood glucose testing gives you correct feedback. It also allows you to make prompt changes based on the results. It alerts you to high blood glucose readings that need changing. It allows you to adjust your insulin to carbohydrate ratio for certain foods. Your pump only uses fast-acting insulin; as a result, you have no long-acting insulin as a back-up. Therefore, if your insulin delivery is disrupted on pump therapy, your blood glucose can go dangerously high fairly fast. This can happen much faster than it could when you were using daily injections with long-acting insulin. Blood glucose testing is needed to alert you to high blood glucose so you can prevent diabetic ketoacidosis (DKA).

A1C

A1C testing also plays a key part in diabetes management. Sustained high BG levels can cause serious long-term problems. These problems may be prevented and/or delayed if you maintain your blood glucose levels close to normal. The best measure of your overall blood glucose is the A1C test. It has been proven that an A1C level of 7.0 or less greatly reduces the risk of problems from diabetes. But, any reduction in A1C is a plus. ADA Clinical Practice Recommendations state that your A1C should be tested at least every three months by your healthcare professional.

Your healthcare professional will rely on the results of your blood glucose tests to make changes in your pump settings. The results of your A1C tests are used to assess your overall blood glucose control.

If you are on pump therapy, you must test your BG at least:

- Upon waking up
- Before each meal

- At bedtime
- If you feel nauseous
- One hour after any BG above 13.9 mmol/L (250 mg/dL)

Below are target guidelines from the American Diabetes Association for adults with diabetes. Consult your healthcare professional to find out what your targets should be.

Summary of recommendations for adults with diabetes mellitus	
Glycemic control A1C	<7.0 percent
Preprandial plasma glucose (before meals)	5.0-7.2 mmol/L (90-130 mg/dL)
Peak postprandial plasma glucose (1-2 hours after meals)	<10.0 mmol/L (<180 mg/dL)

Correlation between A1C level and blood glucose levels		
	Mean plasma glucose	
A1C (percent)	mg/dL	mmol/L
6	135	7.5
7	170	9.5
8	205	11.5
9	240	13.5
10	275	15.5
11	310	17.5
12	345	19.5
American Diabetes Association. Diabetes Care. “Tests of Glycemia in Diabetes.” & “Standards of Medical Care.” Vol. 31, Supplement 1, January 2008. S18.		

Low blood glucose (hypoglycemia)

Low blood glucose can occur while using the pump for the same reasons it can while using daily injections:

- **not enough food**
- **too much insulin**
- **more exercise than usual**
- **drinking alcoholic beverages**

As you know, you cannot always avoid low blood glucose. It is important that you create a routine to follow when your BG is low. If you have a routine, you will have something ready to treat a low BG. This makes it less likely that you will over-treat your low BG, and cause your BG to go up too much.

It is best to use something that you can always carry with you. Choose a food that is all carbohydrate because it will be fast acting. Avoid using high fat foods such as chocolate. They will not work fast enough, and you may also end up over-treating your lows. When treating hypoglycemia, start with 15 grams of fast acting carbohydrates. Do not keep eating until you feel better. Eat the required amount and stop. The feelings will pass and you will be grateful that you did not overdo it.

Some people with diabetes know when their blood glucose is low and some do not. If you are not aware of hypoglycemia, it is important to test your blood glucose more often. Everyone with diabetes should test before driving a car to assure safety on the road. Blood glucose should be above 5.6 mmol/L (100 mg/dL) before driving and at bedtime.

Tip: We suggest that you use glucose tablets to treat a low. These tablets contain a known, pre-measured amount of carbohydrate, and are in handy packages that you can carry in your pocket, purse, or car.

Hypoglycemia protocol: the rule of 15

Check with your healthcare professional for guidelines on treating your low BG. The following instructions are commonly used to treat a BG that is 3.9 mmol/L (70 mg/dL) or lower:

- 1 Eat 15 grams of fast acting carbohydrate.
- 2 Check BG again in 15 minutes; if not above 3.9 mmol/L (70 mg/dL), repeat treatment.
- 3 Check BG every 15 minutes, and repeat treatment until your BG is above 3.9 mmol/L (70 mg/dL).
- 4 Contact your healthcare professional if your BG does not go above 3.9 mmol/L (70 mg/dL) after repeated treatment.

These items have 15 grams of fast acting carbohydrates:

- Glucose tablets (three, 5-gram tablets or four, 4-gram tablets)
- 120 mL of juice or soda (not diet)
- 6-7 Life Savers® (hard candies)
- 15 mL of table sugar or honey

High blood glucose (hyperglycemia)

High blood glucose can occur while using the pump for the same reasons it can while using daily injections:

- too much food
- not enough insulin
- loss of insulin strength
- disruption of insulin delivery from the pump

The goal of treating hyperglycemia is to prevent Diabetic Ketoacidosis (DKA) and delay or prevent diabetes problems due to high blood glucose over a lengthy period of time.

If for any reason you are not getting the proper amount of insulin, your blood glucose rises quickly. This can occur with insulin pump therapy from the disruption of insulin delivery. This happens from the infusion set coming out, clogs, or leaks, or insulin not being absorbed right.

Since the pump only delivers fast-acting insulin, hyperglycemia can occur rapidly. Your healthcare professional will give you data needed for you to determine your correction bolus. This correction dose is based on your insulin sensitivity.

It is vital that you know these guidelines in the hyperglycemia protocol.

Hyperglycemia protocol

If one blood glucose reading is above 13.9 mmol/L (250 mg/dL):

- Take a correction bolus right away.
- Test BG in one hour.

If the second blood glucose is above 13.9 mmol/L (250 mg/dL):

- Take an insulin injection by syringe (not through the pump). The amount should be the same as if you were taking a correction bolus.
- Change the entire infusion set system (new reservoir, infusion set, and cannula). Consider changing the insulin vial if you believe that the insulin is no longer stable.
- Check urine for ketones and call your healthcare professional if ketones are present.
- Drink liquid with no calories every 30 minutes (for example: 240 mL diet ginger ale, broth, water).
- Test blood glucose every two hours and keep taking correction insulin until blood glucose reaches target.
- Call your healthcare professional if your blood glucose and urine ketones remain high or you are unable to drink.

Treating hyperglycemia from illness or infection through insulin pump therapy will be discussed on the following pages.

Diabetes ketoacidosis (DKA)

DKA results from untreated hyperglycemia. DKA is a serious medical problem that needs treatment right away. The insulin pump therapy uses only fast-acting insulin, so DKA can occur rapidly if insulin delivery is disrupted. It is important that you understand these guidelines.

Preventing DKA

If you have nausea or vomiting, immediately check your blood glucose and ketones. If your blood glucose is above 13.9 mmol/L (250 mg/dL) and/or ketones are present:

- Take an injection of fast-acting or regular-acting insulin with a syringe (not through the pump). The amount should be the same as if you were taking a correction bolus. If ketones are present, more insulin may be needed. Contact your healthcare professional for guidelines for insulin doses especially when ketones are present.
- Change the entire infusion set system (new reservoir, infusion set, and cannula). Consider changing the insulin vial if you believe that the insulin is no longer stable.
- Call your healthcare professional.
- Troubleshoot the pump. If you need help, please refer to the Medtronic Diabetes International Contacts list at the beginning of this user guide for your local representative.
- Drink liquids with no calories every 30 minutes (for example: 240 mL diet ginger ale, broth, water).
- Check your blood glucose and ketones in one hour.
- Continue to take insulin as discussed with your healthcare professional.
- Call your healthcare professional right away if your blood glucose and ketones are not decreasing, or you cannot drink fluids.
- Continue to test blood glucose and ketones. Continue to take correction insulin injections as told until your blood glucose reaches your target level.
- If your blood glucose is less than 11.1 mmol/L (200 mg/dL) and ketones are present, drink liquids with calories. Some examples of caloric fluids are juice and non-diet soda. Also, extra insulin may be needed. Contact your healthcare professional for guidelines for insulin doses when ketones are present.

Tip: You should have a meter, glucose strips, and ketones strips handy at all times. This way you are always ready to test your blood glucose and ketones.

Sick day management

Managing diabetes during an illness or infection requires frequent blood glucose and urine ketone testing. Illness and infection put extra stress on the body and often raise blood glucose. The insulin pump allows you to make changes. It lets you quickly and easily respond to blood glucose changes due to illness and infection.

Even if you are not able to eat, you need insulin. Depending on blood glucose testing, you may or may not need to adjust your basal insulin to cover your needs. You may also need to take frequent correction boluses during illness.



Sick day protocol

- Test your blood glucose every two hours, 24 hours a day.
- Check urine ketones every time you urinate.
- If you are throwing up and/or have ketones, you must call your healthcare professional. They will give you the proper treatment with fluids and insulin to avoid DKA.
- Consult your healthcare professional about taking more insulin if your blood glucose is 13.9 mmol/L (250 mg/dL) or higher. Also take more insulin when ketones are fair or large.
- Keep exact records of your blood glucose values, ketones, medication, temperature and all other signs.
- Keep in mind that you need more insulin and fluids when there are urine ketones. **This is true even if your blood glucose is within your target range.**

Sick day supplies

You should have the supplies listed below in the house at all times, and carry them with you when you travel:

- Fluids that have sugar (regular soda, juice, jello) to replace solid food.
- Sugar-free liquids (diet drinks, bouillon, water) for replacing lost fluids.
- Thermometer.
- Sugar free medicine for fever, cough, congestion, nausea and vomiting.
- Extra blood glucose and ketone strips.
- Glucagon emergency kit in case of severe hypoglycemia. This can be used if you cannot eat or are unconscious.

Caring for your diabetes when you are sick is worthwhile and can be a challenge. If you are too sick to watch your diabetes carefully, ask a friend or family member to help. If there is no one to help you, ask your healthcare professional for help.

Nutrition

Good nutrition and making healthy choices are important for everyone. Additionally, for those with diabetes, it is more important to know the nutrient content in food, and how the nutrients affect blood glucose.

Carbohydrates have the most effect on blood glucose, notably within a few hours of being eaten. Counting carbohydrates allows you to match your insulin dose to the food you are eating. This is your insulin to carbohydrate ratio. Although fat and protein can affect your blood glucose when eaten in large amounts, carbohydrates affect blood glucose the most.

Consult with your healthcare professional for your personal nutritional guidelines.

Carbohydrate counting

With insulin pump therapy, you must take a bolus for the carbohydrates you eat. You need to count the carbohydrates you plan to eat to see how much insulin to take.

There are two basic methods to count carbohydrates, and in truth, they are very much alike. Many people use parts of both of the two methods.

Carbohydrate gram counting

With this method, you add up the exact number of grams of carbohydrate in each meal. Food labels, food list and meal planning books are great tools for carbohydrate gram counting.

Carbohydrate exchange system

This method uses food groups called exchange groups. One exchange with carbohydrate has about 15 grams of carbohydrate.

1 bread exchange	1 fruit exchange	1 milk exchange
15 grams carbohydrate	15 grams carbohydrate	12 to 15 grams carbohydrate

Foods with carbohydrates include:

- starches and starchy vegetables
- fruits and fruit juices
- milk and milk products
- sugars and foods with sugar

Facts on the nutrition labels on packaged food will help you to figure how much carbohydrate is in a certain food. There are a number of paperback books that list the grams of carbohydrates in many foods. These books are very helpful when you are eating foods without labels or are eating out. Also, if you know the exchange system, you can count carbohydrate exchanges.

Your healthcare professional will give you an insulin to carbohydrate ratio at the start. This will be a place to begin. It may need to be changed based on the results of your blood glucose readings after you eat.

Determining your meal bolus based on carbohydrate counting is not an exact science. Sometimes it takes trial and error to come up with the correct meal bolus for certain foods. If you are not sure how some foods affect your blood glucose, test your blood glucose two hours after eating. If your reading is too high or too low, adjust your meal bolus the next time you make that food choice.

As you start to learn carbohydrate counting, try to assess the carbohydrates in the food you are eating. Then look at the food labels to see how close you are.

Reading a food label

Use the Serving Size and Total Carbohydrates information on the food label to determine how many total grams of carbohydrate you will be eating. Then, use the insulin to carbohydrate ratio prescribed by your healthcare professional. This allows you to figure out how much insulin to take for the foods you eat.

The following food label shows 12 grams of carbohydrates for each 1 piece serving. If your carbohydrate ratio is 1 unit of insulin for every 12 grams of carbohydrate, you would need 1 unit of insulin for each 1 piece serving of this food.

Nutrition Facts Serving Size 1 piece (20g) Serving per Container 22
Amount Per Serving:
Calories: 50
Fat Calories: 0
Total Fat 0g
Sodium 0g
Total Carbohydrates 12g
Sugars 11g
Protein 0g

Fat and blood glucose

Although fat does not have carbohydrate, it may affect your blood glucose level. Your blood glucose can be normal two hours after a meal with high fat foods. However, four to six hours later your blood glucose will increase. That is because the fat slows down the rate the carbohydrates are absorbed in that meal.

Foods with fat include:

- margarine
- salad dressing
- butter
- cream cheese
- oils
- nuts
- sour cream
- bacon
- mayonnaise
- fried foods
- whole milk
- high fat meats

Doctors and health experts suggest that all people, with or without diabetes, eat less total fat—namely less saturated fat. If you have high cholesterol or are trying to lose weight, it is key to limit the amount of fat you eat per day.

There are bolus type options available on your pump that can be helpful in handling the rare high-fat meal. Before using bolus type options, you should be at ease using basic pump features and consult with your healthcare professional.

Protein and blood glucose

Protein has little effect on blood glucose when eaten in moderate amounts, and normally does not require any insulin.

Sources of protein include:

- meat
- eggs
- tofu
- fish
- cheese
- dried beans
- poultry
- peanut butter
- dried peas

Protein needs are the same for people with and without diabetes, unless kidney disease is present. Most people eat more protein than they need. Protein does not need to be part of every meal. In general, 140-200 grams of protein a day is enough. If you have high cholesterol or are trying to lose weight, eat leaner proteins such as chicken and fish.

Other things that affect blood glucose

Fiber may slow down the rise in blood glucose after a meal. Fiber grams are not absorbed. They can be subtracted from the Total Carbohydrate amount for foods or meals with more than 5 grams of fiber. Fiber is an important part of a healthy diet. Suggestions for fiber intake are the same for people with and without diabetes.

Sodium (salt) has no effect on blood glucose. Moderate sodium intake is suggested for people with or without diabetes. If you have high blood pressure or are affected by sodium, limit the amount of sodium in your diet.

Alcohol can cause low blood glucose by slowing down the release of glucose from the liver. This glucose release keeps blood glucose in target between meals. Therefore, food should be eaten while having a drink with alcohol. The carbohydrate content of different types of alcohol can be found in carbohydrate counting books. Consult your healthcare professional for guidelines on bolusing for alcohol.

Vitamins and minerals have no effect on blood glucose.

Healthy eating habits

Insulin pump therapy gives you the choice of when, what and how much to eat without impacting blood glucose control. This gives you the flexibility to sleep-in, postpone or skip meals, and eat snacks—the same as people who do not have diabetes. Be careful to not abuse this freedom and compromise your nutritional health. All people need to make healthy food choices to be mentally and physically healthy.

It is not an exact science to figure your meal bolus based on carbohydrate counting. Sometimes it takes trial and error to figure the correct meal bolus for certain foods. If you are unsure how certain foods will affect your blood glucose, test your blood glucose two hours after you eat. Keep records of foods that you eat on regular basis. Note the bolus amounts and post-prandial readings for each entry. If your reading is too high or too low, adjust your meal bolus the next time you make that food choice.

Exercise

Your body needs insulin and carbohydrate to provide cells with energy during exercise. How much insulin and carbohydrates you need varies with the types of exercise and with how often you exercise. Using an insulin pump during exercise allows you to reduce your basal rate instead of eating extra carbohydrate. Before, you needed to eat extra carbohydrate to make up for the glucose lowering effect exercise often has on blood glucose. It takes trial and error to get the balance right.

The body needs insulin during exercise. Therefore, it is not suggested that you stop the pump unless the exercise is for a short amount of time (one hour or less). If you need to remove the pump for more than one hour, refer to the *If you remove your pump* section in *The basics* chapter.

If you do an exercise that on average lowers your blood glucose, you can use the temporary basal rate feature. You can use this feature to lower the amount of insulin delivered and reduce the risk of hypoglycemia. You may also choose to change your basal rate pattern for the days you exercise.


It takes a period of time to determine how much to adjust your basal rate during exercise. Record your blood glucose before, during and after exercise to figure the proper reduction in your basal rate. Record any carbohydrates that you have eaten during the exercise period. There is no magic formula that will tell you just what to do. You have to test often to figure out the basal rate that you need during exercise.

The basics

Your pump

CAUTION: Never use sharp objects to press the buttons on your Paradigm pump as this can damage the buttons or compromise the seal of the pump. Some examples of sharp objects that may damage your keypad are fingernail files, pens or pencils, paper clips, knives, scissors, and keys.

Take a look at your pump. The reservoir window allows you to view the insulin in the reservoir. The reservoir, with the tubing connector attached, is inserted into the reservoir compartment of the pump.

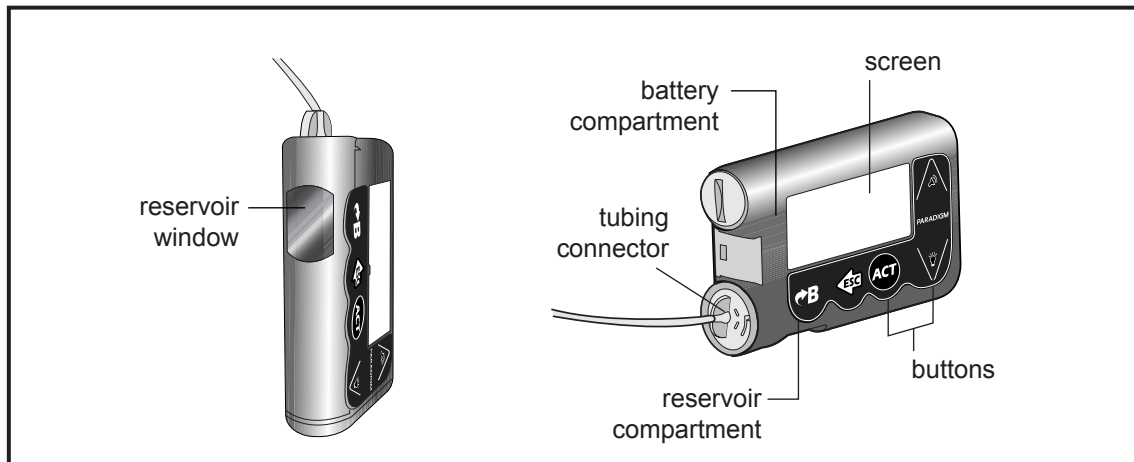


Install battery

CAUTION: Do not use a rechargeable or carbon zinc battery in your pump. For best results use a new Energizer AAA alkaline battery, size E92, type LR03.

Medtronic Diabetes designed the pump to only accept a new battery. As a safety measure, if you install a battery that does not have full power, the WEAK BATTERY or FAILED BATT TEST alarm may sound. If you receive a WEAK BATTERY alarm, respond to the alarm and continue. The pump will still operate normally, but with a decreased battery life. The pump uses one AAA alkaline battery.


- 1 Make sure all the following apply:
 - Clear (ESC, ACT) any alarms and/or alerts before removing and replacing the battery.
 - Make sure the pump is at the HOME (idle) screen when you remove the battery.
 - Do NOT remove the battery during a bolus or Fill Cannula delivery.
- 2 Use the edge of a coin to remove the battery cap. Turn the cap in a counter-clockwise direction.

- 3 Remove the old battery and dispose of it per the disposable requirements of your state or country. Put the new battery in the pump with the negative end [(-) symbol] going in first. Check the label on the back of the pump to make sure the battery is inserted correctly.

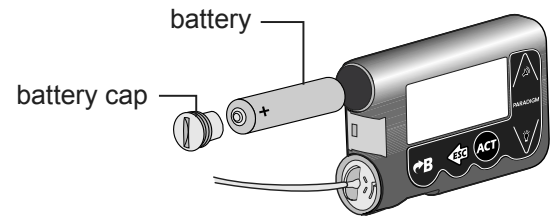

NOTE: Do not use batteries that have been in cold storage, such as in the refrigerator or in your car during winter in cold climates.

- 4 Place the battery cap in the pump and tighten so the slot is aligned horizontally with the pump as shown here:

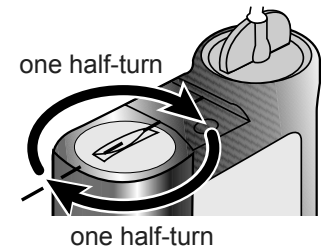

CAUTION: Do NOT over-tighten the battery cap. You should not turn the cap more than four half turns. If you over-tighten the cap you may not be able to remove it, and you can damage your pump.

- 5 While the pump turns on, it will show one or more screens until the HOME screen appears.

If the HOME screen does not appear, do these steps:


- a. Check that the battery is inserted correctly. If the battery has been installed backwards, remove the battery and install it properly.
 - b. If your pump still does not turn on or you get a FAILED BATT TEST alarm, remove and replace the battery with a new one.
 - c. If the pump is still not on, call your local help line or representative.
- 6 Check to make sure the time and date are correct. If more than five minutes have passed since you removed the battery, you will be prompted to check the time and date. Refer to the *Setting the time and date* section in the *Basic programming* chapter for programming instructions.
- 7 Press **ESC** to view the STATUS screen, making sure no alarms are active. If an alarm is active, follow the instructions on the screen.








HOME screen

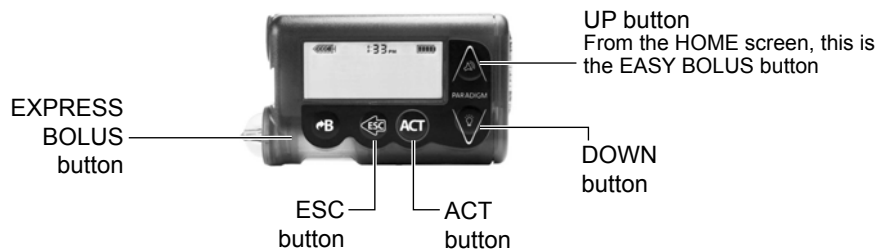Pump buttons






The buttons on the pump are used to navigate through the menus and screens, and to program the features of the pump.



The following table describes how to use the buttons on your pump from the **HOME** screen:

Button	Description
	EASY BOLUS™ button - Shortcut to set and deliver an Easy Bolus.
	Turns the backlight on or off.
	Opens the MAIN MENU.
	<ul style="list-style-type: none"> • Opens the pump STATUS screen if the sensor feature is turned off. • Opens the following screens and graphs from the HOME screen if the sensor feature is turned on: <ul style="list-style-type: none"> • 1 press opens the last sensor graph selected. • 2 presses open the pump STATUS screen. • 3 presses open the SENSOR STATUS screen.
	EXPRESS BOLUS button - Shortcut to set any bolus.

The following table describes how to use the buttons on your pump from the **menus and programming screens**:

Button	Description
	Increases or decreases the value of a flashing item. Scrolls up or down the items in a list.
	Accepts a selected menu item or activates a selected setting.
	<ul style="list-style-type: none"> • Returns to previous screen or exits the menu. • Backs out of unintentional menu selections if the ACT button has not been pressed yet.
	<ul style="list-style-type: none"> • Press simultaneously with  to turn on backlight when in the menus. • Use as a Shift button by pressing it in combination with another button to access certain features.

The pump screen


The screen shows five lines of text at one time. The first is the operating mode. The second is the current open menu or function. The last three lines show either information or text that you can select for the current function.

NOTE: The screen text in the examples used in this guide might not exactly match the text on your pump screen. Please follow your pump screen instructions. If you have any questions, contact your local representative.

HOME screen

The HOME screen serves as the starting point to access the programming screens. When no buttons are pressed for about 30 seconds, the pump returns to this screen.

When you press **ACT** from the HOME screen, the MAIN MENU will appear.

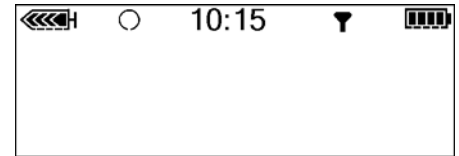

When the pump is on, the following icons always appear across the top of the screen: reservoir volume icon, the time (12- or 24-hr), and the battery icon. If these do not appear, the pump is not operating.

Screen icons

There are various icons that appear at the top of your pump screen, like the time, battery and reservoir icons mentioned above. The next sections describe what the icons mean.

Battery

The battery icon tells you how much usable life is left in your battery. There are four segments in the icon. Each segment represents approximately 25 percent of the usable battery life you have left until you reach Low Battery point. So if you only have one segment left, make sure you have a new battery available.




At least 75% left



Low battery


Time display

The current time of day is displayed across the top of the pump screen in the format you select—12-hour or 24-hour. The AM or PM is only displayed for the 12-hour format. For instructions on setting the time on your pump, see the *Setting the time and date* section in the *Basic programming* chapter.

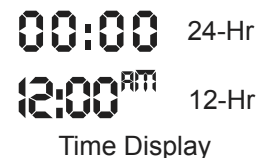

Reservoir volume


The reservoir volume icon tells you how much insulin is in your pump. This icon is also divided into four segments. Each segment represents approximately 25 percent of the reservoir volume you have left. This allows you to be aware of how much insulin you have left in your pump. Refer to your STATUS screen to view the number of units left in the reservoir.

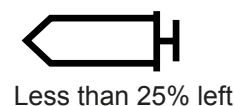

NOTE: For the Paradigm 754 pump, your reservoir icon will only appear full if using a filled 300-unit Paradigm reservoir. The icon will not appear full if you use the 176-unit Paradigm reservoir.


Alert and alarm icons

An open circle (alert) or a solid circle (alarm) are displayed in the upper part of your pump screen only when there is an alert or alarm condition on your pump. For alarm and alert information, see the *Troubleshooting, alerts and alarms* chapter.

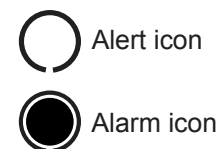

Sensor icons

A dark sensor icon means that the sensor is on and the pump is receiving sensor data. A light-colored sensor icon with a dark circle around it means that the sensor is on, but the pump is not receiving sensor data. For information about the sensor, see the *Sensor features* chapter.

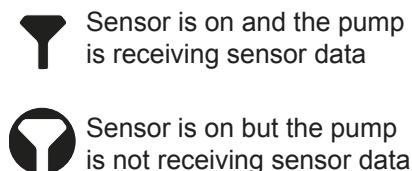




Scroll bar

If there is more text than the screen can show, a scroll bar appears on the right side of the screen. Press  to view any additional text.



Screen backlight

When you press  from the HOME screen, the light on the screen turns on or off. During programming, the backlight can be turned on by pressing the  and the  together. The light will stay on while you are pressing any of the pump buttons. It will stay on as long as the current screen is active.

To conserve your battery, the backlight will turn off automatically while the pump is vibrating. After the vibration is finished, the light will turn back on. The backlight cannot be turned on in a low or dead battery condition.

Beep/vibrate


Your pump will beep or vibrate to indicate activity. Refer to the *Setting your alert type* section in the *Utilities* chapter for setup instructions.

Operating modes

The screen lets you know when a special feature is active or if there is a condition that needs your attention. The active features and pump status will determine the operating mode. The screens for the three modes are shown in the following sections.


Normal mode

Mode for standard pump operations for normal basal and bolus delivery. No special features are active (such as basal patterns, temp basal, and others). No alarms and no alerts exist.

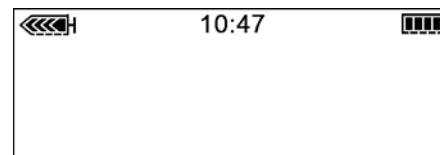

Special mode

Indicates a special feature is active or an alert condition(s) exists. Special mode does not restrict any of the pump functions. When the pump is in Special mode, an open circle appears at the top of the screen and it beeps or vibrates periodically to remind you of the condition. The conditions and features that put the pump in Special mode are:

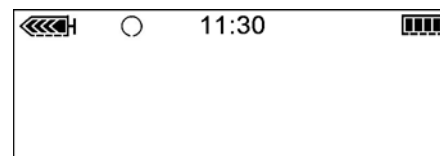


- Low Reservoir condition
- Low Battery condition
- Block feature is on
- Dual or Square Wave bolus delivery
- Basal pattern A or B is active
- Sensor alerts
- Temporary basal delivery
- Alert Silence is active

Attention mode

WARNING: When the pump is in Attention mode, showing a solid circle icon, all insulin deliveries are stopped.

Indicates insulin delivery has stopped. This can mean that the pump is in Suspend mode. It can also mean an alarm is active or an alarm condition exists that needs immediate attention for insulin delivery to resume. A solid circle appears at the top of the screen and the pump will beep periodically until either the pump is taken out of Suspend mode or the condition is cleared. The screen will show text describing the condition that put the pump in Attention mode. For example, if the reservoir is empty, EMPTY RESERVOIR will appear on the screen.

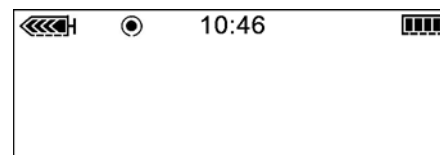

When the pump is in Attention mode, it will beep or vibrate periodically to remind you of the condition. The beep or vibrate frequency varies depending on the condition that put the pump in Attention mode. Refer to the *Alarm conditions* section in the *Troubleshooting, alerts and alarms* chapter for alarm conditions that will trigger the Attention mode. See the *Stopping your pump* section in the *Basic programming* chapter to learn about Suspend mode.

Menus


The MAIN MENU is the highest level menu. There are submenus, functions, status and programming screens in the lower menu levels. The menus are described in the following paragraphs.

Tip: If a screen item is flashing (blinking), during programming, press  or  to change the value.

MAIN MENU

Highest menu level in the menu system. When you press **ACT** from the HOME screen, the MAIN MENU will appear.

BOLUS MENU

Contains the settings and functions for bolus deliveries. The  button allows direct access to MANUAL BOLUS or to the BOLUS WIZARD feature without having to navigate through the menus. Refer to the *Basic programming* chapter for the manual bolus information or to the *Using the Bolus Wizard feature* chapter for bolusing using the Bolus Wizard feature.

SUSPEND

Stops all current insulin deliveries. Refer to the *Stopping your pump* section in the *Basic programming* chapter for more information.

SENSOR

Contains the functions to setup the interface between the sensor and the pump and access to the pump's sensor features.

CAPTURE EVENT

Contains the options you use to enter information about different events into the system. Examples of events include: the amount of insulin used for injections, and the carbohydrates consumed during meals or snacks. The Capture Event option appears in the MAIN MENU after this feature has been turned on in the UTILITIES MENU menu. See the *Capturing Events* section in the *Basic programming* chapter for more details.

BASAL MENU

Contains the functions to setup and deliver your basal. Refer to the *Basal* section in the *Basic programming* chapter for more information.

RESERVOIR + SET

Contains the functions required to change your reservoir and fill the infusion set with insulin. Refer to *Changing your infusion set* section in the *Starting on insulin* chapter for more information.

UTILITIES MENU



Contains features for your safety and convenience. Refer to the *Utilities* chapter for more information.

STATUS screen

The STATUS screen shows information about what your pump is doing. Only check your pump status (press **ESC**) when you are not programming your pump. If you press **ESC** during programming, you will cancel the settings you are trying to enter.



The information that shows on the STATUS screen depends on the current activities and conditions of your pump.

- To open the STATUS screen, press **ESC** until the STATUS screen appears.
- To view more text on STATUS screen, press  or  to scroll and view all of the information.
- To exit the STATUS screen, press **ESC** until the STATUS screen disappears.

The screen includes information about:

- Last bolus insulin delivery
- Current basal insulin delivery
- The amount of active insulin, if any

- Special features that are turned on
- Last alarm/alert that was cleared within 24 hours
- Most recent BG meter reading received within 24 hours
- Reservoir status
- Battery status
- Time and date

Refer to the *Pump specifications* chapter for a complete list of the information that is available in the STATUS screen.

If you remove your pump

You may have an occasion when you need or want to remove your pump. If you have to remove and store your pump, it is recommended that you do the following:

- Store your pump with the battery in place
- Keep a record of your current basal rates and use the *Save Settings* feature (see the *Saving the settings* section in the *Utilities* chapter).
- To preserve battery life, reset the basal rates to 0 (zero), turn off the RF options (meter, remote), and set Auto-off to dashes or zeroes.

Remember, your body still needs insulin while your pump is removed.

It is important that you consult with your healthcare professional to determine an alternate method of receiving your insulin. You can remove your pump for up to one hour without taking insulin. If you remove your pump for more than one hour, you will have to use another way to take your insulin, such as injections of fast-acting insulin, or reconnecting your pump to take boluses. Take the injection or bolus approximately every four hours. Calculate the amount of insulin to take based on the total of your basal insulin in four hours. Include the amount you need for meal and correction boluses. If you will have your pump off for several days, you will need to return to a multiple injection regimen.

Basic programming

Setting the time and date


Setting the correct time and date in your pump is necessary for accurate basal insulin delivery and allows you to keep an accurate record of your insulin delivery and other pump functions. You can select a 12-hour or 24-hour clock. You must reset the time and date if you receive a CHECK SETTINGS alarm or you clear your settings (Clear Settings function).

- 1 Go to the TIME/DATE SETUP screen.

Main > Utilities > Time/Date



- 2 Select 12-Hour Setup or 24-Hour Setup and press ACT.

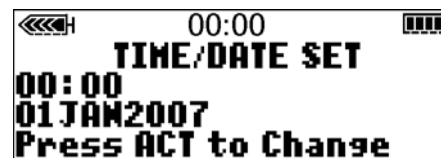



- 3 Press ACT again to change the settings.

4 Change each of the settings as follows:

a. Hour


Change the hour. Press **ACT**.

For 12-hour setups, press  or  until the correct A (am) or P (pm) appears.



b. Minutes

Change the minutes. Press **ACT**.

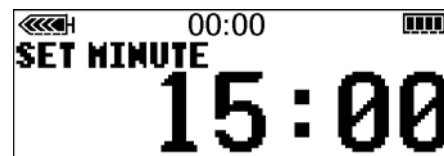


c. Year

Change the year. Press **ACT**.




d. Month

Change the month. Press **ACT**.

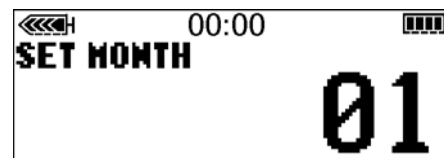

e. Day

Change the day. Press **ACT**.

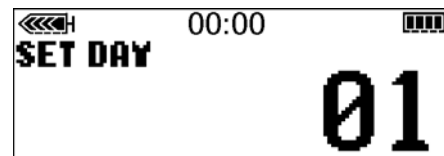

- 5 The TIME SET AT screen will show the settings that you programmed. Press **ACT** and exit the menus.
Your time/date settings are complete.



Selecting the language

The language shown on the pump screens can be changed. Some languages may not be available on all pumps. Before you can select another language, you need to set the time using the English screens. Refer to the previous section.

To change the language for your pump:

- 1 Go to the LANGUAGE MENU screen.
Main > Utilities > Language



- 2 Select your language, then press **ACT**.



- 3 The language setting is now changed. It may take several seconds for the new language to display on the screen. Exit the menus.

Bolus

There are three bolus types: Normal, Square Wave[®], and Dual Wave[®]. This section gives instructions for a Normal bolus using the **EXPRESS BOLUS** button and navigating through the menus. (For information about Square Wave and Dual Wave boluses, refer to the *Optimizing pump therapy* chapter.)

The Normal bolus delivers an immediate food or correction bolus. It can be delivered at any time except during another Normal bolus. During a Normal bolus, most pump features are disabled until after all the bolus has been delivered. The suspend function and the STATUS screen, however, are always available.

Setting the Normal bolus

Normal bolus can be used to cover the carbohydrate in a meal or snack and/or to correct a blood glucose that is higher than your blood glucose target.

The following instructions are for a Normal bolus when the Bolus Wizard feature is turned off.

- 1 Go to the BOLUS MENU.

Main > Bolus

Select **Set Bolus** and press **ACT**. Go to step 2.

You can also use the  **EXPRESS BOLUS** button to get to step 2. Press  from your HOME screen.

- 2
 - a. If the **SET BOLUS** screen appears: (Dual/Square option is off) Go to step 3.
 - b. If the **BOLUS TYPE** screen appears: (Dual/Square Wave is on) Select **Normal Bolus** and press **ACT**. Go to step 3.

- 3 The **SET BOLUS** screen (or **SET NORMAL BOLUS** screen if the Dual/Square option is on) appears with 0.0 unit bolus amount flashing. Enter your bolus amount and press **ACT**.

If some of the bolus insulin delivered by the pump has not been used and is active in your body, the screen will show this amount as Act. Ins. (Active Insulin). This amount is based on the active insulin setting selected in the Bolus Wizard feature. The Active Insulin amount shows even if the Bolus Wizard feature is turned off. Take this amount into consideration to determine your bolus. See the *About active insulin* section in the *Using the Bolus Wizard feature* chapter for more information.

NOTE: *If you have BG Reminder turned On, the BG REMINDER DURATION screen displays. It allows you to set the duration before you are reminded to check your blood glucose after a bolus. See the BG Reminder section in this chapter for information about this feature.*

- 4 The **BOLUS DELIVERY** screen appears and the Normal bolus starts. The pump will beep/vibrate at the start of the bolus. As the bolus delivers, the amount shown on the screen will increase until the entire bolus has been delivered. When the bolus is finished, the pump will beep/vibrate again and the HOME screen will appear.

The following practice lessons will help you understand this pump feature.

Normal meal bolus using the exchange system

Normal bolus can be used to cover the carbohydrate in a meal or snack and to correct a blood glucose that is higher than the target that was chosen for you.

Fred has been taught that he needs to take 1 unit of insulin for every carbohydrate exchange that he eats (every milk, every starch or every fruit). For lunch today he will eat:

Turkey sandwich with two slices of bread **2** starches

1 small apple **1** fruit

240 mL of non-fat milk **1** milk

Total carbohydrate exchanges = 4

Fred's lunch has a total of 4 carbohydrate exchanges so he will take a meal bolus of 4 units for his lunch.

Bolus practice:

NOTE: Make sure that you are not connected to your pump while practicing.

Going through the menus, program a 2.0 unit Normal bolus now.

Check here if you were able to program it. ☐

Using the **EXPRESS BOLUS** button **⬇B**, program a 2.0 unit Normal bolus now.

Check here if you were able to program it. ☐

Normal meal bolus practice using exchanges

Choose a meal you might eat and fill in the blanks.

Food: _____

exchange: _____

exchange: _____

exchange: _____

total exchanges: _____

You will take _____ units of insulin for each exchange. Your total bolus is _____ for this meal.

Normal meal bolus using carbohydrate counting

Lydia has been taught that she needs to take 1 unit of insulin for every 10 grams of carbohydrate. This is her insulin to carbohydrate ratio. For dinner she will have:

112 grams broiled chicken 0 grams

160 ml of rice 30 grams

120 ml cooked broccoli 5 grams

28 grams dinner roll 15 grams

5 ml margarine 0 grams

total grams of carbohydrates = 50 grams

Lydia's dinner totals 50 grams of carbohydrate. Her insulin to carbohydrate ratio is 1 unit to 10 grams. She will take a meal bolus of 5 units for her dinner. She determined this by dividing 50 (total grams of carbohydrate) by 10 (insulin to carbohydrate ratio).

Normal meal bolus practice using carbohydrates

Choose a meal you might eat and fill in the blanks.

Food: _____

grams of carbohydrate: _____

grams of carbohydrate: _____

grams of carbohydrate: _____

total grams of carbohydrate: _____

Your insulin to carbohydrate ratio: 1 unit of insulin for _____ grams carbohydrate.

Divide your total carbohydrates by your insulin to carbohydrate ratio and take _____ units of insulin for your meal.

Meal bolus, correction bolus and insulin sensitivity

Jason is ready to eat his breakfast. He has calculated that he will need 4.0 units for his food.

He tests his blood glucose and finds that it is 11.1 mmol/L (200 mg/dL). Jason knows that his blood glucose level is above his blood glucose target and will need additional insulin before he eats.

Jason's healthcare professional has determined the following for him:

BG Target: 6.1 mmol/L (110 mg/dL)

Insulin sensitivity¹: 2.0 mmol/L/u (36 mg/dL/u)

Jason determines that he will need a correction bolus of 2.5 units insulin to lower his elevated blood glucose. The 2.5 correction bolus will lower his current blood glucose of 11.1 mmol/L (200 mg/dL) to his target of 6.1 mmol/L (110 mg/dL).

- elevated blood glucose: $11.1 - 6.1 = 5$ mmol/L ($200 - 110 = 90$ mg/dL)
- correction bolus: $5.0 / 2.0$ mmol/L/u (insulin sensitivity) ($90 / 36$ mg/dL/u) = 2.5 units

(He will add this 2.5 correction bolus to the 4.0 units of insulin that he will need for his meal bolus. Jason will take a total bolus of 6.5 units.)

1. *Insulin sensitivity is the amount (in mmol/L or mg/dL) by which blood glucose will be lowered after taking 1 unit of insulin.*

Consult with your healthcare professional to determine your insulin sensitivity.

Practice: Meal bolus

You have determined your meal bolus as: _____ units.

Your target blood glucose range is: _____ to _____ (average is _____).

Your current blood glucose level is: _____.

Your correction factor is: 1 unit of insulin will drop your blood glucose _____.

You will take _____ unit(s) of insulin to correct your high blood glucose level.

Your total bolus (meal bolus plus correction bolus) is _____.

Review your bolus deliveries

You can view a list of your bolus deliveries in the BOLUS HISTORY screen. This screen shows a list of the dates, times, units, and types for your last 24 boluses. This feature is helpful for record keeping or to check if you bolused for your last meal.


If a bolus was stopped before delivery was complete, the BOLUS HISTORY screen will show only the amount actually delivered. Refer to the *Bolus details* section in this chapter for instructions about viewing bolus details.

Do the following steps to view the BOLUS HISTORY screen:

- 1 Go to the BOLUS HISTORY screen and scroll through the bolus deliveries.

Main > Bolus > Bolus History

If you used the Bolus Wizard feature to deliver any of these boluses, the BOLUS HISTORY screen shows the carbohydrate/food (CH) and blood glucose values that the Bolus Wizard feature used to calculate the boluses.

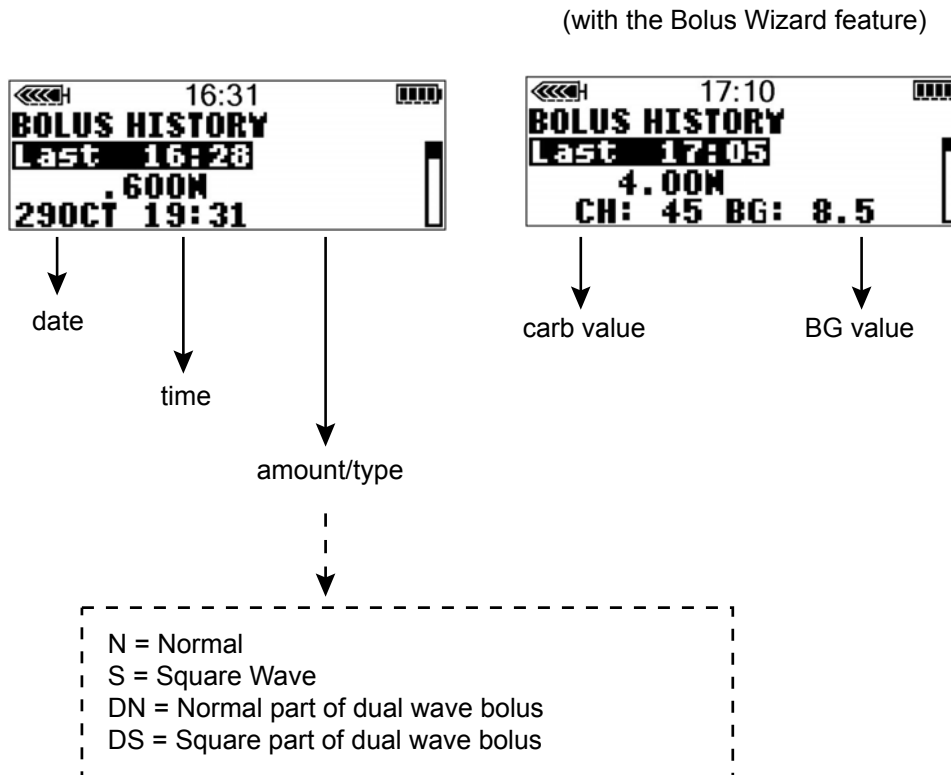

2 Refer to the instructions in the *Bolus details* section to see the details for any of these boluses.

Bolus details

You can view the details of any of the deliveries in the BOLUS DETAIL screen. The details include:

- bolus types: normal, square, and dual
- programmed bolus amount

- delivered bolus amount
- Bolus Wizard feature information (if used)

To see the details of any bolus, do these steps:

- 1 In the BOLUS HISTORY screen, select the bolus that you want to review and press **ACT**.
- 2 The details for that bolus will appear on the screen. Scroll through the details.

Without the Bolus Wizard feature

delivery is active

```

<<<H 12:48
BOLUS DETAIL
Type:      Normal
User entry: 4.00U
Delivered: 1.40U
  
```

delivery stopped

```

<<<H 12:48
BOLUS DETAIL
Type:      Normal
User entry: 4.00U
Stopped at: .650U
  
```

delivery completed

```

<<<H 12:48
BOLUS DETAIL
Type:      Normal
Delivered: 4.00U
  
```

With the Bolus Wizard feature

If the Bolus Wizard feature calculated your bolus, more information will appear in the BOLUS DETAIL screen.

```

<<<H 17:10
BOLUS DETAIL
Type:      Normal
Delivered: 4.00U
Bolus wizard:
  
```

Press the down key to see the details

- 3 Exit when you are done.

Maximum bolus limit

The Maximum bolus (Max bolus) is a safety feature that limits the amount of insulin that can be delivered in a single bolus. The factory setting is 10.0 units. You can specify the limit from 0.0 to 75.0 units. It is important to discuss this feature with your healthcare professional to determine your Max bolus amount.

To set the Maximum bolus limit, do these steps:

- 1 Go to the MAX BOLUS SETUP screen.
Main > Bolus > Bolus Setup > Max Bolus
- 2 Set your Max bolus limit and press **ACT**. If you set your Max bolus limit to more than 25.0 units, the pump displays a **WARNING** screen because the set limit is above the typical limit. To accept the set Max bolus limit, press **ACT**. To change this limit, press **ESC**, and then set a lower Max bolus limit.
- 3 Your Max bolus is set. Exit the menus.

Example 1: Max bolus

Shelby takes very small doses of insulin for her meal boluses. As a safety limit, she and her healthcare professional reset her pump with a Maximum bolus of 5.0 units.

Example 2: Max bolus

David is a growing teenager. He loves to eat big meals and requires very large doses of insulin for his food. He reset his pump with a Maximum bolus of 35.0 units so he can take more insulin when he needs to.

Scroll rate

The scroll rate feature allows you to program bolus insulin delivery in preset increments. These increments are effective on all bolus screens, except for the SET EASY BOLUS screen. See the *Easy bolus setup* section in the *Optimizing pump therapy* chapter for more information on how to set up the steps for your Easy Bolus. The factory setting for the scroll rate feature is 0.10 units, and the available options are 0.025, 0.05, or 0.10 unit increments.

When you set the scroll rate to 0.025, you can use the following increments to program a bolus:

- 0.025 unit increments for a bolus of 0.975 units or less
- 0.05 unit increments for a bolus between 1 unit and 9.95 units
- 0.10 unit increments for a bolus of 10 units or more

When you set the scroll rate to 0.05, you can use the following increments to program a bolus:

- 0.05 unit increments for a bolus of 9.95 units or less
- 0.10 unit increments for a bolus of 10 units or more

When you set the scroll rate to 0.10, the increments will increase or decrease by 0.10 units during bolus programming regardless of the bolus value.

To program your scroll rate:

- 1 Go to the SET SCROLL RATE screen.
Main > Bolus > Bolus Setup > Scroll Rate
- 2 The scroll rate value appears flashing. Change the value and press **ACT**.
- 3 The pump returns to the BOLUS SETUP screen. Your scroll rate is now programmed and ready to use. Exit the menus.

Example: Scroll rate

You have set your scroll rate to 0.025 units. You need to set your pump to deliver a 2.5 unit bolus. The bolus values on the SET BOLUS screen will change in 0.025 unit increments until you reach 1.0 unit. After this the value will change in 0.05 unit increments until you reach 2.5 units in the SET BOLUS screen.

BG Reminder

After you deliver a bolus you may want to check your blood glucose. The BG Reminder is an optional feature that makes the pump beep or vibrate to remind you to check your blood glucose after a bolus. However, this reminder is not available after an Easy Bolus. Your pump is set at the factory with the BG Reminder feature turned off.

When setting a bolus, if the BG Reminder is on, your pump will ask you to set up the BG Reminder Duration. This sets how long after bolus delivery you will be reminded to check your BG. This time can be from 30 minutes to 5 hours, or NONE. The BG Reminder is not available after an Easy Bolus.

- 1 Go to the BG REMINDER SETUP screen.
Main > Bolus > Bolus Setup > BG Reminder
- 2 Select **On** and press **ACT**. The BG Reminder is now enabled. Exit the menus.
Now, the next time you program a bolus, your pump will ask you for the amount of time after your bolus before you want to be reminded to check your blood glucose.

When the BG Reminder goes off, your pump will beep or vibrate and CHECK BG will appear on the screen. Your pump will beep or vibrate periodically until it is cleared (ESC, ACT).

When you set a BG Reminder after a bolus, the STATUS screen will show the amount of time remaining before the reminder goes off. Here, the STATUS screen indicates that the BG Reminder will go off in 18 minutes.



Capturing Events

The Capture Event feature electronically saves certain types of information. Some examples include:

- Blood glucose (BG) measurements
- The amount of insulin you use
- The amount of carbohydrates you eat or drink
- The exercise you do
- Other (this can be used for anything else as needed)

Make sure that you enter events when they happen since the system records the time of the entry. You cannot change entries after you have put the information into your device. The entered information can be sent to Medtronic CareLink® Personal software. Here it can be used to generate treatment reports that you can share with your healthcare professional. You can view up to 10 of the latest saved events in the CAPTURE HISTORY screen.

Turning on the Capture Event feature

The Capture Events option appears in the MAIN MENU after this feature has been turned on in the UTILITIES MENU.

To turn on the Capture Event feature, do these steps:

- 1 Go to the CAPTURE EVENT ON/OFF screen:

Main > Utilities > Capture Option

- 2 Select **On**, then press **ACT**.

You will now have a Capture Event option in the MAIN MENU.

Entering BG measurements

You can save non-calibrating BG readings in device history. Non-calibration meter BG measurements are those that are taken at times when you will not use them for calibration, such as after eating or when your blood glucose is rising or falling rapidly. You also have the option of using the entered reading for calibration purposes. Do **not** calibrate if your sensor is not communicating with your device.

Before you can enter BG readings, you must select the BG units for the Capture Event features. You can set the BG units to mg/dL or mmol/L.

To set the BG units:

- 1 Go to the BLOOD GLUCOSE UNITS screen.
Main > Capture Event > BG Units
- 2 Select either **mg/dL** or **mmol/L**, then press **ACT**. The CAPTURE EVENT menu shows.
- 3 You can now enter your BG reading.

To enter the BG measurement:

- 1 Make a note of the BG measurement.
- 2 Go to the ENTER BG screen.

Main > Capture Event > Enter BG

The ENTER BG screen flashes with dashes or the most recent meter BG measurement you entered.

- 3 Enter your new measurement, then press **ACT**. The BG measurement must be between 1.1-33.3 mmol/L (20-600 mg/dL). A message asks if you want to save the new measurement that is shown on the screen.
 - If the information is correct, select **Save BG**, then press **ACT**. Your BG is saved.
 - If the Bolus Wizard feature is off, the CAPTURE EVENT screen appears.
 - If the Bolus Wizard is on and the saved BG is below your BG Target range, the BG Saved screen appears. This screen indicates that your saved BG is below your BG Target set in the Bolus Wizard feature, and the pump beeps three times. See *Setting the BG Targets, on page 83* for more information.
 - If the Bolus Wizard is on, the saved BG is above your BG Target range, and the calculated correction estimate is more than the preset scroll rate increment, the BG Saved screen appears. This screen indicates that your saved BG is above your BG Target set in the Bolus Wizard feature, and the pump beeps three times. Go to the ENTER BG screen of the Bolus Wizard feature to enter a correction bolus. See *Setting the BG Targets, on page 83* for more information.

- If you want to use the new BG measurement for calibration, select **Save and Calibrate**. This option is available when the Auto Calibration feature is **Off** and your sensor is communicating with the pump. The BG measurements must be between 2 mmol/L (40 mg/dL) and 22.2 mmol/L (400 mg/dL).
 - If the Bolus Wizard is on and the saved BG is below your BG Target range, the BG Saved screen appears. This screen indicates that your saved BG is below your BG Target set in the Bolus Wizard feature, and the pump beeps three times. See *Setting the BG Targets, on page 83* for more information.
 - If the Bolus Wizard is on, the saved BG is above your BG Target range, and the calculated correction estimate is more than the preset scroll rate increment, the BG Saved screen appears. This screen indicates that your saved BG is above your BG Target set in the Bolus Wizard feature, and the pump beeps three times. Go to the ENTER BG screen of the Bolus Wizard feature to enter a correction bolus. See *Setting the BG Targets, on page 83* for more information.
- If the information is not correct, select **Cancel**, then press **ACT**. The CAPTURE EVENT screen reappears with **Enter BG** selected. Repeat the procedure and enter the correct information.

Entering insulin injection information

Follow these steps to enter the amount of insulin you use.

- 1 Go to the ENTER INSULIN screen.

Main > Capture Event > Insulin Marker

The ENTER INSULIN screen flashes with dashes or the insulin amount that you most recently entered.

- 2 Enter the amount of insulin you used, then press **ACT**. A message appears asking you if you want to save this information. The **Yes** option is selected.
- 3 Make sure that the amount of insulin shown on the screen is correct.
 - If the information is correct, press **ACT**. The information you entered is saved to the system and can now be used in reports.
 - If the information is **not** correct, select **No**, then press **ACT**. The CAPTURE EVENT menu shows. Repeat the procedure to enter the correct information.

Entering carbohydrate information

This section shows you how to enter information about the carbohydrates you eat or drink.

Before you begin

Before you can enter carbohydrate information, you must select the carbohydrate units (Carb Units) for the Capture Event feature. You can set the carbohydrate units to grams or exchanges.

To set the carbohydrate units:

- 1 Go to the CARB UNITS screen.
Main > Capture Event > Carb Units
- 2 Select either **Grams** or **Exchanges**, then press **ACT**. The CAPTURE EVENT menu shows. You are now ready to enter your carbohydrate information.

To enter carbohydrate information:

- 1 Determine the total units of carbohydrates in the meal or snack that you plan to eat.
- 2 Go to the ENTER FOOD screen.
Main > Capture Event > Meal Marker
- 3 The ENTER FOOD screen flashes with dashes or with the number of carbohydrate grams or exchanges you entered last time.
- 4 Enter the carbohydrate grams or exchanges, then press **ACT**. A message asks if you want to save the information that is displayed on the screen. The **Yes** option is selected.
- 5 Make sure the number shown on the screen is correct.
 - If the information is correct, press **ACT**. The information you entered is saved to the system and can now be used in reports.
 - If the information is **not** correct, select **No**, then press **ACT**. The CAPTURE EVENT menu shows. Repeat the steps above to enter the correct information.

Entering exercise information

Follow these steps to enter exercise information. Be consistent and enter the marker either before or after each time you exercise.

- 1 Follow the path below to save your exercise information:
Main > Capture Event > Exercise Marker
- 2 A message asks if you want to save this exercise.
- 3 Make your selection, then press **ACT**.

Entering Other markers

This section shows you how to enter markers other than Enter BG, Insulin Marker, Meal Marker, or Exercise Marker. Examples of Other markers include: when you take medications, when you feel ill, when you are under stress, and so on. The Other marker appears in CareLink Personal reports, and shows at what time you have entered it in the Capture Event feature.

To enter Other markers:

- 1 Follow the path below to save your Other marker:
Main > Capture Event > Other
A message asks if you want to enter this marker.
- 2 Select **Yes**, then press **ACT**. Your pump saves your Other marker and returns to the CAPTURE EVENT screen.

Viewing Capture Event history

You can view up to 10 captured events saved into the device. The most recent event will be displayed at the top of the CAPTURE HISTORY screen.

To view the Capture Event history:

- 1 Go to the CAPTURE HISTORY screen:
Main > Capture Event > History
- 2 The CAPTURE HISTORY screen displays the saved events. Each event has the date and time when it was saved, the name, and any detail associated with that event.


Basal

Basal insulin is required to maintain your target glucose values when you are not eating. Your healthcare professional will calculate this rate for you. Your basal insulin accounts for approximately one half of the body's total daily insulin requirements. Your pump mimics your pancreas by delivering insulin continuously over 24 hours.

You can set your insulin pump to change rates during the day to match your needs. Your needs depend on your lifestyle and insulin requirements. Some people only use one rate throughout the day, while others find they need more. Your basal rates are made up of insulin deliveries that have start and stop times. Once set, these rates make up your 24-hour basal pattern and are repeated daily.

Start and stop times

When you set your basal rate(s) in the BASAL MENU, your pump prompts you to set the start time for each basal delivery. The start time of one basal rate is the stop time of the previous rate. This gives you continuous basal insulin through a 24-hour period. For basal rate 1, the start time is midnight (00:00 or 12:00 a.m.) and cannot be changed. See the following sample graph of basal rates.

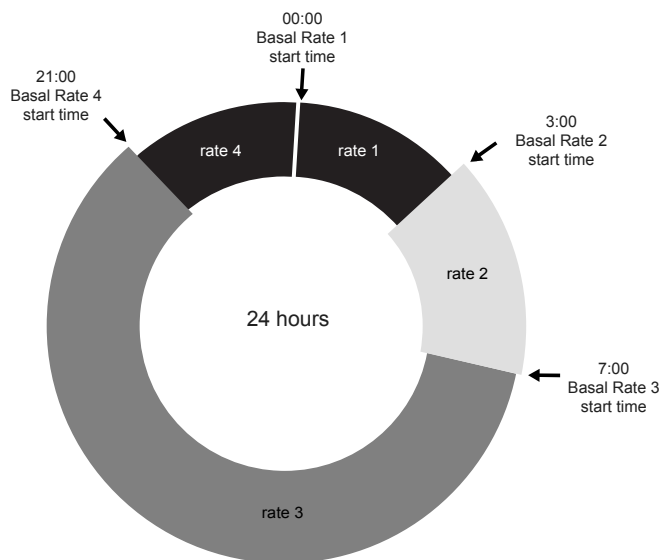

The start and stop times of basal rates for this example are:

- 00:00 (12:00 a.m.) to 2:59 (2:59 a.m.) for basal rate 1
- 3:00 (3:00 a.m.) to 6:59 (6:59 a.m.) for basal rate 2
- 7:00 (7:00 a.m.) to 20:59 (8:59 p.m.) for basal rate 3
- 21:00 (9:00 p.m.) to 23:59 (11:59 p.m.) for basal rate 4

If you need one basal insulin that starts before midnight and stops after midnight, you need to set up two basal rates. These basal rates have the same rate but different start and stop times. For example, in this graph a basal rate is needed from 21:00 to 2:59 (9:00 p.m. to 2:59 a.m.) at one rate. Basal rate 1 always starts at midnight and should be set to end at 2:59 (2:59 a.m.). After setting basal rates 2 and 3, rate 4 is set from 21:00 (9:00 p.m.) to 23:59 (11:59 p.m.) at the same rate as basal rate 1. This gives continuous basal insulin from 21:00 to 2:59 (9:00 p.m. to 2:59 a.m.) at one rate.

You cannot set a start time for one basal rate to overlap the next basal rate. The addition of a new basal rate will erase any basal rates that follow.

For best results, setting or changing your basal rate(s) should be discussed with your healthcare professional.

Your basal settings

You must program your basal settings before you can deliver basal insulin. Keep a written record of your basal settings.

It is recommended that you set your basal rates with the assistance of your healthcare professional.

If you plan to take off your pump for an extended period of time, such as more than a day, set the basal rate to 0.000 U/H. This will ensure that the insulin delivery records in your pump are accurate. Refer to the *If you remove your pump* section in *The basics* chapter for more information.

Basal programming and delivery

NOTE: *You cannot make changes to your basal rate settings while a percent temp basal is active.*

To set your basal rates:

- 1 Go to the SET BASAL RATE 1 screen.
Main > Basal > Basal Setup > Set/Edit Basal
- 2 The SET BASAL RATE 1 screen flashes the basal rate in U/H.
- 3 Enter your first basal rate amount and press **ACT**.
 The start time for your first basal rate is midnight and cannot be changed.
- 4 The SET START TIME 2 screen appears. The dashes under the screen name flash. The first basal rate is now set.
 If you do not need a second basal rate for the day, press **ESC**. If you need to set up a second basal rate for the day, follow steps 5 and 6.
- 5 In the SET START TIME 2 screen, enter the start time for the next rate.
- 6 Press **ACT**. The SET BASAL RATE 2 screen appears. Enter the rate.
- 7 Press **ACT**. The SET START TIME 3 screen appears. The second basal rate is now set. If you do not need to set up any more basal rates for the day, press **ESC** and skip to step 9. If you need to set up more basal rates, follow steps 5 and 6 for each basal rate.
- 8 After you program your last basal rate, press **ESC**.

- 9 The BASAL RATE screen appears. Your basal rate(s) will now deliver as programmed. Exit the menus.
- NOTE:** The pump delivers basal amounts in 0.025 U/H increments for basal rates of less than 1.00 U/H, and in 0.05 U/H increments for basal rates of 1.00 U/H or more.
- 10 When you have finished programming your pump, you can save your settings. Refer to the *User settings* section in the *Utilities* chapter for instructions.

Current basal delivery

The STATUS screen shows your current basal information.




Daily basal rate(s)

The BASAL REVIEW screen shows your daily basal rates programmed for delivery from midnight to midnight. Compare your daily insulin deliveries to your blood glucose records to help you and your healthcare professional identify your optimal daily basal insulin rate(s).

To review your basal rates, do these steps:

- 1 Go to the BASAL MENU screen.
Main > Basal
- 2 Select **Basal Review** and press **ACT**.


- 3 If you do not use patterns, the STANDARD screen appears.
The delivery details for your standard basal will appear.

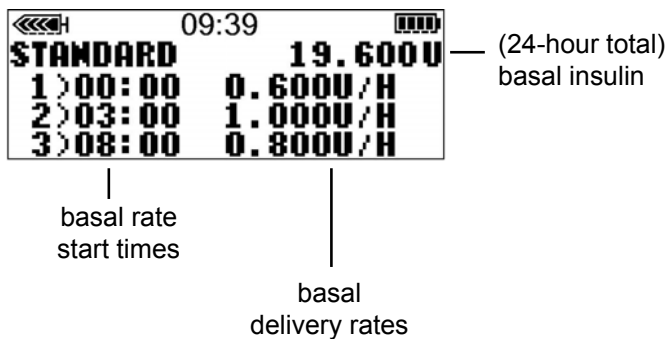

If you use patterns, the BASAL REVIEW screen appears.

The current basal pattern will be highlighted. Select the pattern you want to view. Press **ACT**.

The delivery details for that pattern will appear.

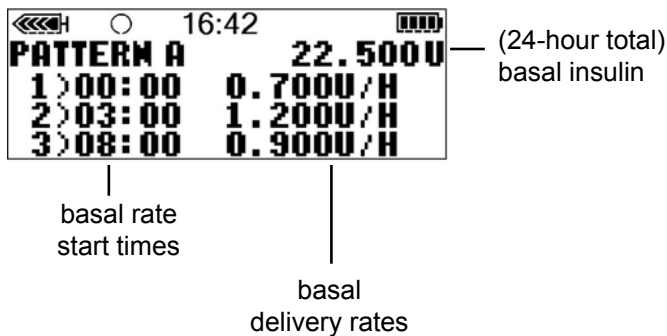

- 4 Exit the menus when you are done.

Setting the Max basal rate

Maximum basal (Max basal) rate is a safety limit for the amount of basal insulin that is able to be delivered per hour. Your pump is sent from the factory with the Max basal set to 2.00 units per hour. It is important to determine your Max basal rates with your healthcare professional. This safety feature will not allow you to program any basal rates, including patterns and temporary basal rates, that are greater than the Max basal rate.

If you are setting your Max basal rate after your basal rates have been set, you **cannot** set a Max basal that is less than any of the programmed basal rates.

To set your Max basal rate, do these steps:

- 1 Go to the MAX BASAL RATE screen. The Max basal rate will be flashing.
Main > Basal > Basal Setup > Max Basal Rate
- 2 Change the rate and press **ACT**.
- 3 Your Max basal rate is now set. Exit the menus.

Example 1: Max basal

Helen has a very low insulin requirement. Her highest basal rate is only 0.400 units per hour. As a safety measure, Helen's healthcare professional set her pump with a Maximum basal rate of 1.00 units per hour.

Example 2: Max basal

Rusty needs large amounts of insulin to control his blood glucose levels. His new pump was delivered from the factory with a Maximum basal rate of 2.00 units per hour, but he needs 2.80 units per hour in the early morning. Rusty will reprogram his Maximum basal to 3.00 units per hour to accommodate his needs.

Stopping your pump

Suspend stops all insulin delivery, including the current basal and any bolus or Fill Cannula deliveries that are in progress. While suspended, your pump will not deliver insulin until you resume your pump. When you resume your pump, the basal delivery will continue.

The pump will beep or vibrate about every 15 minutes on the hour to remind you that it is not delivering insulin. Example: You suspend your pump at 11:20. The pump will beep/vibrate at 11:30, 11:45, 12:00, and so on until you resume your pump (basal resumes).

When suspended, your pump is in Attention mode with a solid circle showing on the pump. When in **Suspend**, you can only resume your basal or view the STATUS screen. If the Sensor feature is turned on, you can also view the sensor graphs and the SENSOR STATUS screen. No other functions are available.

To suspend the pump:

- 1 Select **Suspend** from the MAIN MENU, and press **ACT**.
Main > Suspend
- 2 SUSPEND will flash on your screen. Press **ACT** to stop your pump.



- 3 The screen will show that the pump is suspended and the time that it stopped. After less than one minute, the pump will return to the HOME screen with a solid circle.

You can verify on your pump STATUS screen that your pump is suspended.



Resume pump delivery

When the pump is suspended, it defaults to the HOME screen with a solid circle.

To resume your pump and basal delivery, do these steps:

- 1 From any screen, press **ACT** until the RESUME screen appears. Press **ACT** again.



- 2 Your pump will beep once, then the HOME screen will appear without the solid circle.

NOTE: A bolus or fill cannula that was stopped by suspend will **not** restart when you resume your pump. You must reprogram and activate it to finish delivery.

Example:
Suspend function

Helen is ready to eat her lunch. She has just programmed her pump to deliver a meal bolus when the phone rings. Helen wants to talk on the phone and not eat her lunch right away. She knows that if she lets the bolus continue and she does not eat her lunch soon, she may be at risk for low blood glucose. Helen suspends delivery of her pump to stop the bolus, but then resumes delivery to restart her basal insulin. When she is off the phone and ready to eat, she checks her STATUS screen to see how much insulin she received from the partially delivered bolus before she suspended her pump. She will reprogram a new bolus for the remainder of her bolus amount.

Practice:
Suspend function

- 1 **Make sure you are NOT connected to your pump while practicing.**
Program your pump to deliver a Normal bolus of 3.0 units. Once the bolus begins, stop the bolus by suspending your pump.
Remember, when you stop the bolus delivery with SUSPEND, *all* insulin delivery will stop.
- 2 Now, RESUME delivery, so that your basal insulin will continue.
- 3 Check the STATUS screen.
- 4 How much insulin did the bolus deliver before you Suspended delivery? _____.
- 5 If you wanted to take the rest of the bolus later, how much would you take to equal 3.0 units? _____.

Practice:
Resume basal delivery after a suspend

Make sure you are NOT connected to your pump while practicing.

- 1 Give a 3.0 unit bolus now. While it is delivering, suspend the bolus.
- 2 Check here if you were able to suspend the bolus. ☐
- 3 Now restart the pump.
- 4 Check here if you were able to restart the pump. ☐

Starting on insulin

The Paradigm pump is intended for use with U100 insulin.

Prepare your pump for use

Before continuing with the steps in this chapter, we recommend that you watch your pump training CD-ROM if available, and complete your pump start training.

When you are done practicing and ready to use your pump with insulin, you must make sure the time and date are correct on your pump. You must also program your settings as instructed by your healthcare professional.

You will need these items:

- Pump
- Insulin (U100)
- Paradigm reservoir and user guide
- Paradigm compatible infusion set and user guide



Filling the reservoir

WARNING: Using cold insulin can cause air bubbles in the reservoir and tubing. If your insulin is stored in the refrigerator, allow it to reach room temperature before filling your reservoir. When filling the reservoir, take care to remove air bubbles.




CAUTION: Insulin comes in different vial types, depending on your country. If your vial is different than the one described in the User Guide, consult with your healthcare professional on how to fill your reservoir.

Filling the reservoir

- 1 Remove reservoir from package. Make sure plunger rod is fully extended.
- 2 Swab vial with alcohol.
- 3 Making sure you do not push down on the plunger, press the transfer guard onto the vial.

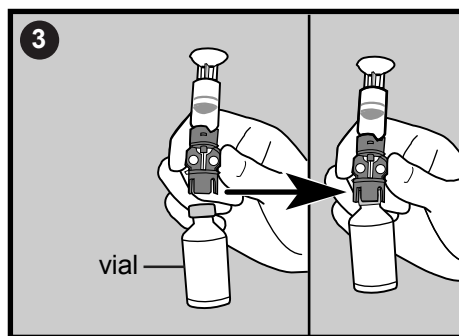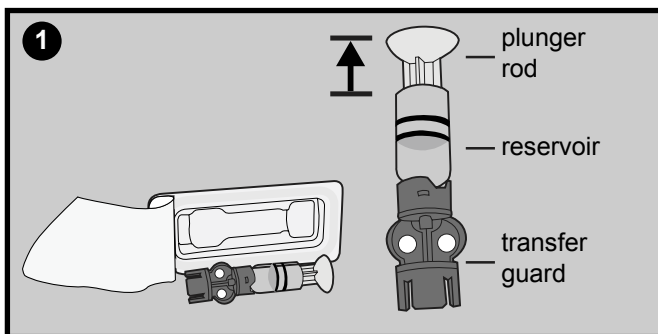


- 4 Push down on the plunger to pressurize the vial.
- 5 While still holding the plunger rod, flip the vial over so the vial is up, slowly pull down on the plunger to fill the reservoir.
- 6 Gently tap the side of the reservoir to make any air bubbles rise to the top of the reservoir.

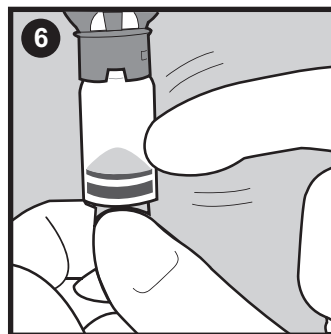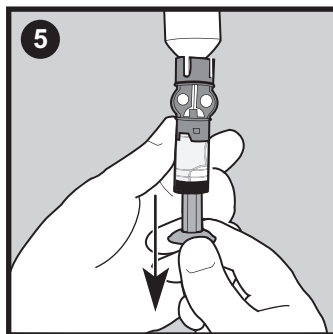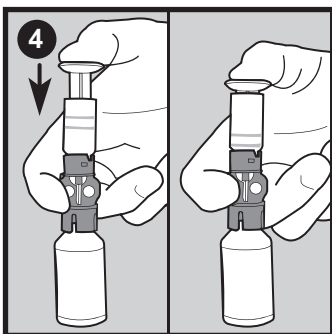


- 7 Slowly push up on the plunger just enough to remove any air bubbles from the reservoir.
- 8 Slowly pull down on the plunger to fill the reservoir to the number of units desired.
- 9 With the vial down, hold the transfer guard. Turn the reservoir counter-clockwise, then pull straight up to remove the vial from the transfer guard.

- 10** Place the tubing connector onto the reservoir. Turn the connector clockwise, pressing gently against the reservoir until you feel it slide in. Push in and continue turning until the reservoir and the connector lock with a click.

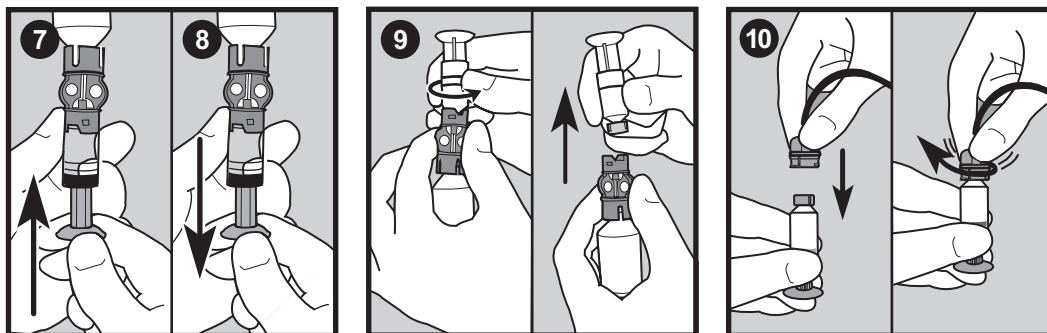

- 11** Tap the side of the reservoir to remove any air bubbles.
- 12** To purge air bubbles that have risen to the top of the reservoir, push up on the plunger until you see insulin in the tubing.
- 13** Without pulling, turn the plunger counter-clockwise to remove it from the reservoir.

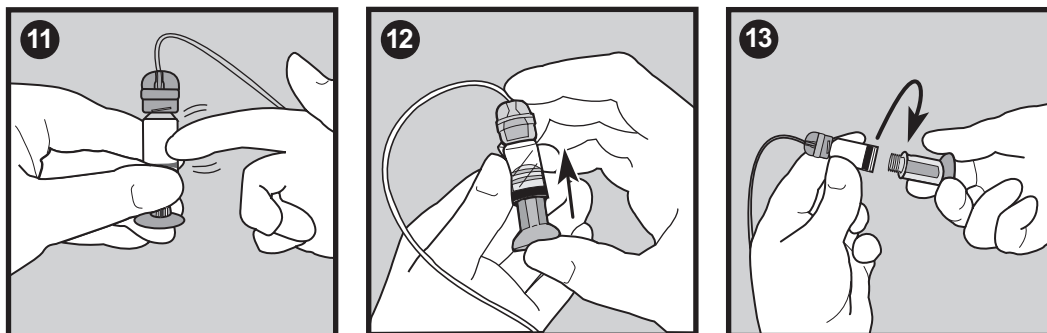


- 14** Use your reservoir immediately after you fill it. Do not store it filled.

Changing your infusion set

Removing the reservoir

Each time you remove and replace a reservoir in your pump, you have to rewind your pump, and fill the infusion set with insulin.

- 1 Remove the entire infusion set from your body.
- 2 If attached, remove the activity guard.
- 3 Turn the tubing connector half-turn counter clockwise, then pull the reservoir and connector out from the pump.

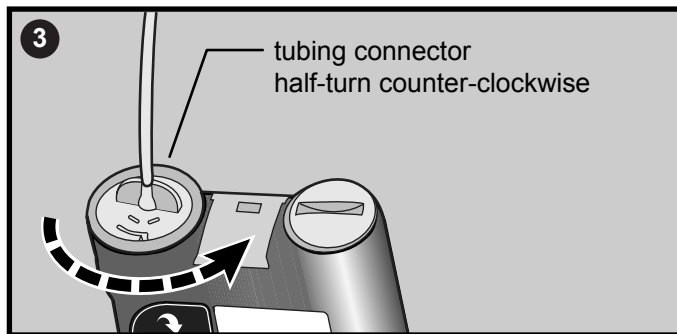

- 4 Safely dispose of the used reservoir and infusion set item in a sharps container.
- 5 You must now rewind your pump as described in the next section.

Rewinding your pump

Before you continue, make sure the infusion set is NOT connected to your body and the reservoir is NOT in the pump.

WARNING: Make sure the infusion set is disconnected from your body before you rewind your pump or fill the infusion set tubing. Never insert the reservoir into the pump while the tubing is connected to your body. Doing so could result in an accidental infusion of insulin.

- 1 If you removed your reservoir and are replacing it, go to the REWIND screen.

Main > Reservoir + Set > Reservoir Setup

- 2 In the REWIND screen, press **ACT** to start the rewind process. The REWINDING screen will appear while the pump rewinds.
- 3 After the pump rewinds the REWIND COMPLETE screen will appear.

If you are practicing:

- a. Do **NOT** insert the reservoir in your pump. Make sure the red shipping cap is installed in the reservoir compartment as shipped.
- b. Press **ACT**, then continue with the instructions in the *Filling the tubing* section in this chapter.

If you are **not practicing**, continue to the next section to insert the reservoir in your pump.

Inserting the reservoir in your pump

If your reservoir is already inserted in your pump, continue to the next section.

You must do these steps in the order described. If you are practicing, do **NOT** insert the reservoir in your pump.

CAUTION: You must rewind your pump before installing a new reservoir. As part of the pump's function, it calculates the reservoir volume. To ensure correct volume calculations, your pump has been designed to require a rewind before you insert your reservoir.


- 1 If you are using the pump for the first time, remove the red shipping cap from the reservoir compartment.

WARNING: Do not insert the reservoir in the pump if you did not rewind. Doing so could result in inaccurate insulin delivery.

Never insert the reservoir into the pump while the tubing is connected to your body. Doing so could result in an accidental infusion of insulin.

- 2 Insert the reservoir into the top of the pump case.

- 3 Turn the tubing connector approximately 1/2 -turn clockwise until the connector is seated. The tubing connector should be aligned horizontally with the pump case as shown here.

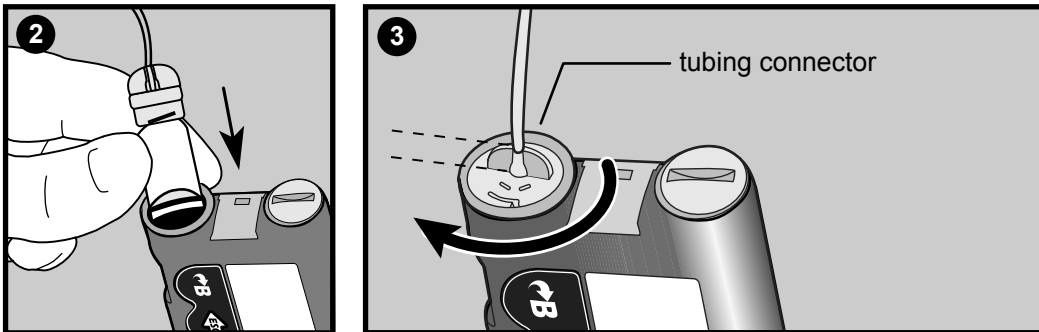

- 4 Attach the activity guard, if desired.
- 5 If the pump has returned to the HOME screen, press **ACT** to show the REWIND COMPLETE screen. Press **ACT** again to go to the FILL TUBING screen.
- 6 You must now fill the infusion set tubing as described in the next section.

Filling the tubing

You need to fill the infusion set tubing with insulin before you insert the set into the body.

WARNING: Make sure the infusion set is disconnected from your body before you fill the infusion set tubing. Never insert the reservoir into the pump while the tubing is connected to your body. Doing so could result in an accidental infusion of insulin.

- 1 After you rewind your pump, the FILL TUBING screen will appear.
 - a. If the infusion set is **NOT** disconnected from your body, press **ACT** to select **No** on the screen. After you disconnect the infusion set from your body, press **ACT** on the next screen to continue.
 - b. If the infusion set is disconnected from your body, select **Yes** on the FILL TUBING screen. Press **ACT** to continue to the next FILL TUBING screen.
- 2 Press and hold **ACT**. The pump will beep six times indicating that the pump is positioning the reservoir.
- 3 After the reservoir is positioned, you need to fill the infusion set tube with insulin. Press and hold **ACT** until insulin droplets form on the tip of the infusion set needle, then release. Follow the instructions on the screen. You will hear beeps as the pump fills the tubing with insulin. Be sure no bubbles are in the tubing.

WARNING: If the FILL TUBING screen does not display, do NOT continue. Do NOT insert the infusion set into your body. Please contact your local representative for assistance.

Take care to remove air bubbles when filling the infusion set tubing.

If you use more than 30 units of insulin to fill the infusion set tubing, the MAX FILL REACHED alarm screen will appear. If you get this alarm, do these steps:

- a. **Make sure that you are not connected to the pump.**
 - b. Read the message on the screen, then press **ESC**, **ACT** to clear.
 - c. To continue filling the infusion set tubing, select **Yes, Continue** and press **ACT**. See the next step to continue. If you have made a mistake, select **No, Rewind** and press **ACT**. The REWIND screen will appear. Refer to the pump rewinding and insulin loading instructions in this chapter to continue.
- 4 A message will appear: DO YOU SEE DROPS AT END OF TUBING?
Do not press **ESC** because this will trigger a FINISH LOADING alarm in 10 minutes. Make sure to complete this Fill Tubing step. Select **Yes** or **No**.
- a. If you do NOT see drops at the tip of the needle, select **No** and press **ACT**. Make sure your infusion set is NOT connected to your body. Follow the instructions on the FILL TUBING screens to continue filling the infusion set tubing with insulin.
 - b. If you see drops formed at the tip of the needle, select **Yes** and press **ACT**. The FILL CANNULA screen will appear.
- 5 You can now insert the infusion set into your body as described in the next section.

Inserting the infusion set

WARNING: While the infusion set is connected to your body, do not unscrew and retighten the tubing connector on the reservoir.


After you complete all of the following, you will be ready to insert the infusion set into your body:

- fill your reservoir
- rewind your pump
- insert the reservoir into pump


- and fill the infusion set with insulin

Shown here are the best body areas (shaded) for infusion set insertion. Avoid the 5.0 cm (2-inch) area around the navel.

It is important that you change your infusion set every two to three days. Be sure to rotate the infusion set insertion sites so that they do not become overused. The abdominal area is the most common site for set insertions because absorption is very consistent. To keep abdominal sites healthy, some people find it helpful to use a visual scheme to help them rotate their insertion sites in an organized way. Here are two commonly used methods. For maximum effectiveness, use both methods, alternating between them:

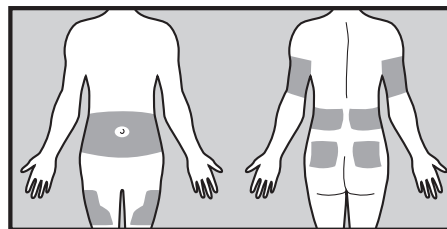

- Visualize an imaginary clock drawn on your abdomen surrounding your belly button. Rotate infusion set insertion sites by starting at 12 o'clock and then rotate the site clockwise to 3 o'clock, 6 o'clock, and so on.

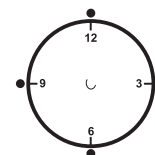

- Imagine a letter M or a letter W on either side of your belly button. Start at the end of one letter and proceed through the letter, rotating to each intersection in turn.




Medtronic Diabetes offers different infusion sets for your pump. Instructions for the Quick-set® begin on the next page as an example. Always refer to the instructions that shipped with your infusion set.

After your infusion set is inserted, see the *Filling the cannula* section in this chapter to fill the infusion set cannula.

Quick-set infusion set (with Quick-serter®)

There are different infusion sets that you can use with your pump. As an example, the following procedure shows how to insert the Quick-set infusion set. Always refer to the instructions that shipped with your infusion set.

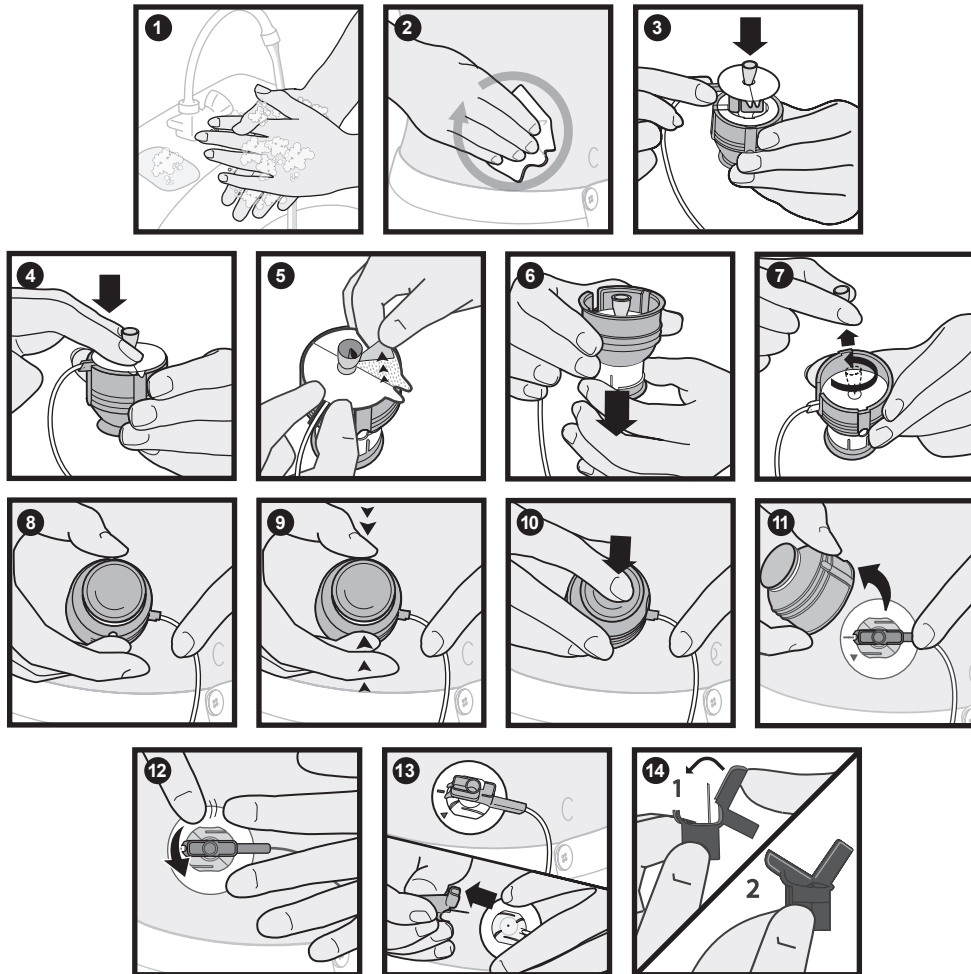

Filling the cannula

Filling the soft cannula with insulin is required after the infusion set is inserted into your body and the introducer needle is pulled out. The insulin amounts required to fill the cannula depend on the type of infusion set you are using. Refer to your infusion set instructions for this information. If you are using an infusion set with a needle, press **ESC** to skip this step.

If you do not complete the Fill Cannula steps, your pump triggers a FINISH LOADING alarm. If you need to skip the Fill Cannula step, press **ESC** to avoid getting this alarm.

- 1 After you have filled the infusion set tubing, the FILL CANNULA screen appears. It displays a message to connect the infusion set to your body. Press **ACT** to fill the cannula. To skip this step, press **ESC** until you see the MAIN MENU screen. If you do not complete this step, the RESERVOIR + SET MENU screen appears after 15 minutes. After 10 minutes, your pump triggers a FINISH LOADING alarm.
- 2 Enter the amount for your type of infusion set, then press **ACT**. To skip this step, press **ESC** until you see the MAIN MENU screen. If you do not complete this step, the RESERVOIR + SET MENU screen appears after one minute. After 10 minutes, your pump triggers a FINISH LOADING alarm.
- 3 As the cannula starts filling, the FILLING CANNULA screen will count up the units as they are delivered. The pump will beep or vibrate when the cannula is full.


To see how much insulin was delivered to fill the infusion set:

- 1 Go to the HISTORY screen.
Main > Reservoir + Set > History
- 2 Scroll through the list of insulin deliveries. The letter T at the end of the text line indicates that insulin was delivered to fill the tubing. The letter C at the end of the text line indicates insulin was delivered to fill the cannula. Exit the menus.

Disconnecting Quick-set


The Quick-set allows you the freedom to temporarily disconnect from your pump without removing the infusion set from your body.

- 1 Hold the side grips of the connector part with your fingers.
- 2 Twist the connector counter-clockwise.
- 3 Remove the connector from the site.

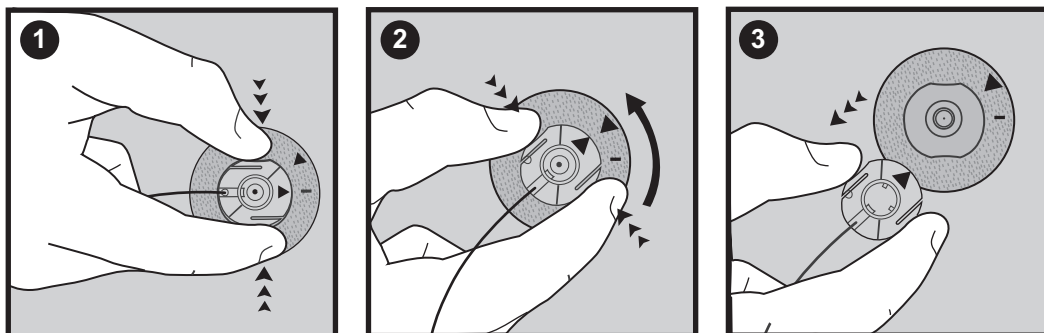

Reconnecting Quick-set

Place the connector part (flat side facing down) on the infusion site until it is fully seated. Do not squeeze the connector part by the flat side grips.

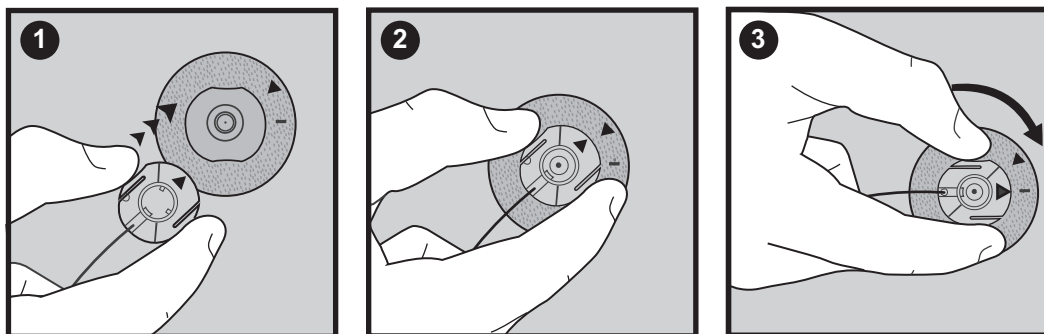

Record keeping for diabetes management

Now that you are using the pump, we will be asking you to test your blood glucose regularly. It is important to test often and write down your blood glucose readings, the food you eat, any exercise you perform and any other notes to explain your blood glucose results.

You must test at the recommended times and any other time that you feel your blood glucose is high or low. Be sure to include your meal boluses, correction boluses, the amount of carbohydrate you eat, basal rate and any other information that will be helpful in assisting your healthcare professional in adjusting your pump settings.

It is very important to look at your blood glucose readings as feedback regarding your diabetes management, not as statements about you or your self-worth. Try not to have an emotional reaction to the numbers and do not judge them too harshly. You will soon learn how to modify the numbers easily and precisely through insulin pump therapy.

Test at least four to six times per day. These are the recommended times to test to determine control:

- Overnight (occasionally, at approximately 2:00 - 3:00)
- Pre-breakfast (fasting)
- Post-breakfast (approximately two hours after eating)
- Pre-lunch
- Post-lunch (approximately two hours after eating)
- Pre-dinner
- Post-dinner (approximately two hours after eating)
- Bedtime
- Before driving

Determining your pump settings

Your healthcare professional will use your daily blood glucose journal records to program your pump. It is very important to keep good records during the first weeks after you start on pump therapy. Not only must you record your blood glucose readings, but it will be important to eat regularly scheduled meals and to keep your activity as consistent as possible.

Until you and your healthcare professional determine the pump settings that will work best for you, it is important to eat meals for which it is easy to count the carbohydrates. After your correct basal rate is determined, you will be able to experiment with varied food choices and amounts.

After you and your healthcare professional are satisfied with your initial pump settings, you may begin to experiment with different food choices, meal times and exercise schedules.

Using the Bolus Wizard feature

What is it?

Bolus Wizard is a feature that calculates an estimated bolus to support your food intake or to correct a high blood glucose. To use this feature, you need the following information:

Carbohydrate counting

You need to know which foods contain carbohydrates and how to count these.

Your blood glucose reading

You need to know your blood glucose (BG) reading. When using the Bolus Wizard feature, the pump can work with blood glucose meter powered by MWT1 technology to automatically receive your blood glucose reading. MWT1 is the wireless Radio Frequency (RF) technology that is used to transmit information from the meter to the pump. You can program your pump to automatically receive your BG reading from this meter. All meters referenced in this user guide are blood glucose meters supported by MWT1 technology. The *Meter option* section in this chapter has more information. If you are not using this meter, you will manually enter your blood glucose.



Your personal Bolus Wizard feature settings

In addition to your blood glucose reading and/or your food entry, the Bolus Wizard feature uses personal settings that you program into the pump. (For instructions see the *How to program the Bolus Wizard feature* section in this chapter.)

- carb units (grams or exchanges)
- carb ratios (in carbohydrate grams/unit of insulin or insulin units/carb exchanges)
- BG units (mmol/L or mg/dL)
- insulin sensitivity

- target blood glucose range
- active insulin time (hours)


Get this information from your healthcare professional, and, for best results, talk to your healthcare professional before making any changes. Keep a record of your settings in the following *Bolus Wizard feature settings* table:

Bolus Wizard feature settings		
Information	Setting	
Carb units:	_____ grams or _____ exchanges	
Carb ratios: Bolus Wizard feature uses this for your food bolus calculations. If you count carbs: this ratio is the amount of carbohydrate grams covered by one unit of insulin. range: 1-200 grams/unit If you count exchanges: this ratio is the amount of insulin you need to cover one (carb) exchange. range: 0.075-15.0 units/exchange NOTE: <i>Your carb ratios may vary throughout the day. Your pump allows you to program up to eight different carb ratios.</i>	#1: _____ #2: _____ #3: _____ (additional settings, if needed) #4: _____ #5: _____ #6: _____ #7: _____ #8: _____	start time (midnight)
BG units: (how you measure your blood glucose)	_____ mmol/L or _____ mg/dL	

Bolus Wizard feature settings		
Information	Setting	
<p>Insulin sensitivity: This ratio is used for your correction bolus calculations. This ratio is the BG units reduced by 1.0 unit of insulin. range: 0.5-22.2 mmol/L/u or 10-400 mg/dL/u</p> <p>NOTE: <i>Your insulin sensitivity may vary throughout the day. Your pump allows you to program up to eight different insulin sensitivities.</i></p>	<p>BG units reduced / 1 unit of insulin</p> <p>#1: _____</p> <p>#2: _____</p> <p>#3: _____ (additional settings, if needed)</p> <p>#4: _____</p> <p>#5: _____</p> <p>#6: _____</p> <p>#7: _____</p> <p>#8: _____</p>	<p>start time (midnight)</p>

Bolus Wizard feature settings		
Information	Setting	
<p>BG Target range: If your current blood glucose is above the BG Target range, the Bolus Wizard feature will calculate a correction dose. If your current blood glucose is below the BG Target range, the Bolus Wizard feature will calculate a negative correction and subtract it from your food bolus.</p> <p>range: 3.3-13.9 mmol/L or 60-250 mg/dL</p> <p>NOTE: Your pump will allow you to program up to eight different BG Target ranges.</p>	#1: _____ #2: _____ #3: _____ (additional settings, if needed) #4: _____ #5: _____ #6: _____ #7: _____ #8: _____	start time (midnight)
<p>Active insulin time: The Bolus Wizard feature uses this time to calculate the active insulin in your system (see the <i>About active insulin</i> section in this chapter). Use your healthcare professional's recommendation for the active insulin time that best represents the insulin type you use and your physiological insulin absorption rate.</p> <p>range: 2-8 hours</p>	Number of Hours: _____	

How the Bolus Wizard feature works

- If you want your current blood glucose to be factored in, enter your BG Reading.
 - automatically from the meter (refer to the *Meter option* section in this chapter), or
 - manually by selecting the  button.
- If you are going to eat, enter your food amount in grams or exchanges.
- The Bolus Wizard feature will calculate a bolus for you. An ESTIMATE DETAILS screen will appear with your estimated total bolus amount.

Bolus Wizard warnings

When using the Bolus Wizard feature, the pump may display a HIGH BG, LOW BG, and MAX BOLUS EXCEEDED warnings.

HIGH BG

If your blood glucose in the ENTER BG screen is above 13.9 mmol/L (250 mg/dL), the Bolus Wizard feature displays a HIGH BG warning. Read the instructions, then press **ACT** or **ESC** to clear the message. You can then continue programming to deliver your bolus.

LOW BG

If your blood glucose in the ENTER BG screen is below 3.9 mmol/L (70 mg/dL), the Bolus Wizard feature displays a LOW BG warning. Read the instructions, then press **ACT** or **ESC** to clear the message. You can then continue programming to deliver your bolus.

MAX BOLUS EXCEEDED

Bolus Wizard feature will not deliver more than the limit set for your maximum bolus. If the Bolus Wizard feature calculates a bolus amount that is larger than your max bolus limit setting, the message, MAX BOLUS EXCEEDED will appear. If this happens, do these steps:

- 1 In the MAX BOLUS EXCEEDED screen, press **ACT** to continue your bolus programming. The EST : MAX screen appears with the estimated and maximum bolus amounts. Continue to the next step.
If you do not want to continue, press **ESC** to cancel and the screen will return to the ENTER BG screen.
- 2 In the EST : MAX screen, press **ACT** again to continue your bolus programming.
- 3 The SET BOLUS screen appears with the maximum bolus amount flashing. Enter the bolus amount. This amount cannot be more than the max bolus amount. Press **ACT**.
- 4 The BOLUS DELIVERY screen appears showing the insulin units being delivered.
- 5 The pump will beep/vibrate after it has completed the insulin delivery.

How to program the Bolus Wizard feature

You need your personal settings from the Bolus Wizard feature settings table to setup the Bolus Wizard feature. Your Bolus Wizard feature settings are programmed in the EDIT SETTINGS screen.

Main > Bolus > Bolus Setup > Bolus Wizard Setup > Edit Settings

Once the settings are programmed, you do not have to program them again unless the values change. After you program one setting, the screen will automatically go to the next required setting. After you program all your settings, review them as described in this section to make sure they are set correctly.

Instructions for programming the Bolus Wizard feature settings are in the next paragraphs. Program your settings in the order described to make sure you program all the settings. If you have not completed setting up all of the required settings, MISSING INFO screen appears. It lists the required settings for this feature. You must program the listed settings before you can use the Bolus Wizard feature.

Turning on the Bolus Wizard feature

- 1 Go to the EDIT SETTINGS screen.
Main > Bolus > Bolus Setup > Bolus Wizard Setup > Edit Settings
- 2 The EDIT SETTINGS screen appears with **Wizard: Off** selected. Press **ACT**.
- 3 The WIZARD ON/OFF screen appears. Select **On**, then press **ACT**.
- 4 The EDIT SETTINGS screen appears, showing that the wizard is now turned on. You are now ready to select your carb units.

Selecting the Carb units

The carb unit setting lets the pump know which way to count your carbohydrates (grams or exchanges).

- 1 Make sure the EDIT SETTINGS screen is open.
Main > Bolus > Bolus Setup > Bolus Wizard Setup > Edit Settings
- 2 Select **Carb Units**, then press **ACT**.
- 3 The CARB UNITS screen appears. Select **Grams** or **Exchanges**, then press **ACT**.
- 4 The EDIT SETTINGS screen shows the carb units you selected. You are now ready to set your carb/exchange ratios.

Setting the Carb/Exch ratios

Your pump allows you to set up to eight Carb/Exch ratios because this ratio may vary throughout the day. Your healthcare professional may only have you program one or two carb ratios when you first start using the Bolus Wizard feature.

To set the Carb/Exch ratios:

- 1 Make sure the EDIT SETTINGS screen is open.
Main > Bolus > Bolus Setup > Bolus Wizard Setup > Edit Settings
- 2 Select **Carb Ratios**, then press **ACT**.
 - **If you use grams as your carb units:** Carb ratio is the number of carb grams that are covered by one (1.0) unit of insulin.
 - **If you use exchanges as your carb units:** Carb ratio is the number of insulin units that are needed to cover one (1.0) carb exchange.
- 3 The SET CARB RATIO 1 (if you are using grams) or SET EXCH RATIO 1 (if using exchanges) screen appears. The default ratio flashes on the screen.
- 4 Set your first ratio, then press **ACT**. Carb ratio values are normally between 5-50 grams/u or 0.3-3.0 u/exch. If your ratio value is outside the range, a warning message will appear on the screen. This message warns that the entered carb ratio is valid but outside the usual range. Press **ESC** to correct or **ACT** to continue.
 The start time for your first ratio is midnight and cannot be changed.
- 5 The SET START TIME 2 screen appears. The dashes under the screen name flash. The first Carb Ratio or Exchange Ratio is now set.
 If you do not need a second ratio, press **ESC** and skip to the next section. If you need to set up another ratio, follow steps 6 through 9.
- 6 In the SET START TIME 2 screen, enter the time of the day you want this ratio to become active.
- 7 Press **ACT**. The SET CARB RATIO 2 screen (if using grams) or SET EXCH RATIO 2 screen (if using exchanges) appears.
- 8 The default ratio flashes. Select your ratio.
- 9 Press **ACT**. The SET START TIME 3 screen appears. The second carb or exchange ratio is now set.
- 10 If you do not need to set up any more ratios, press **ESC**. If you need to set up more ratios, repeat steps 6 through 9 above for each ratio.
 You are now ready to set up the BG units.

Setting the BG units

You can select **mmol/L** or **mg/dL** as your Blood Glucose Unit (measurement type). You can also set these units in the Sensor and Capture Event menus.

- 1 Make sure the EDIT SETTINGS screen is open.
Main > Bolus > Bolus Setup > Bolus Wizard Setup > Edit Settings
- 2 Select **BG Units**, then press **ACT**.
- 3 The BLOOD GLUCOSE UNITS screen displays. Select **mmol/L** or **mg/dL**, then press **ACT**.
- 4 The EDIT SETTINGS screen shows the BG units you selected. You are now ready to set your insulin sensitivity.

Insulin sensitivity

Your insulin sensitivity is the amount your blood glucose (BG) level is reduced by one unit of insulin. This value is used to calculate a suggested insulin dose to correct a high BG. Because this sensitivity may vary throughout the day, your pump lets you set up to eight sensitivity settings. Your healthcare professional may only have you program one or two insulin sensitivities when you first start using the Bolus Wizard feature. Record your settings in the *Bolus Wizard feature settings* table in this chapter.

Insulin sensitivity values are normally between 1.1-5.6 mmol/L/u (20-100 mg/dL/u). If your value is outside this range, a warning message will appear on the screen.

- 1 Make sure the EDIT SETTINGS screen is open.
Main > Bolus > Bolus Setup > Bolus Wizard Setup > Edit Settings
- 2 Select **Sensitivity**, then press **ACT**.
- 3 The INS SENSITIVITY 1 screen appears. The default sensitivity value flashes on the screen.
- 4 Set the value for your first insulin sensitivity setting, then press **ACT**.
The start time for your first insulin sensitivity is midnight and cannot be changed.
- 5 The SET START TIME 2 screen appears. The dashes under the screen name flash. The first insulin sensitivity is now set.
If you do not need a second insulin sensitivity, press **ESC** and skip to the next section. If you need to set up another insulin sensitivity, follow steps 6 through 9.
- 6 In the SET START TIME 2 screen enter the time of the day you want this insulin sensitivity to become active.
- 7 Press **ACT**. The INS SENSITIVITY 2 screen appears.
- 8 The default sensitivity value flashes. Select the value for this insulin sensitivity.


- 9 Press **ACT**. The SET START TIME 3 screen appears. The second insulin sensitivity is now set.
 - 10 If you do not need to set up any more insulin sensitivities, press **ESC**. If you need to set up more insulin sensitivities, repeat steps 6 through 9 above for each insulin sensitivity.
- You are now ready to set up the BG Target.

Setting the BG Targets

The BG Target setting allows you to set glucose targets. The Bolus Wizard will use these targets to calculate a correction dose. Because the targets may vary throughout the day, your pump allows you to set up to eight BG Targets each day. If you want to set just one target value, instead of a range, set both the low and high values to the same number.

If your current BG is above the BG Target range, the Bolus Wizard feature may calculate a correction dose. The correction dose will deliver enough insulin to bring your BG down to your current high end of the BG Target range. If your current BG is below the BG Target range, the Bolus Wizard may calculate a negative correction and subtract it from your food bolus. This will bring your BG to the low end of the BG Target range.

Example:

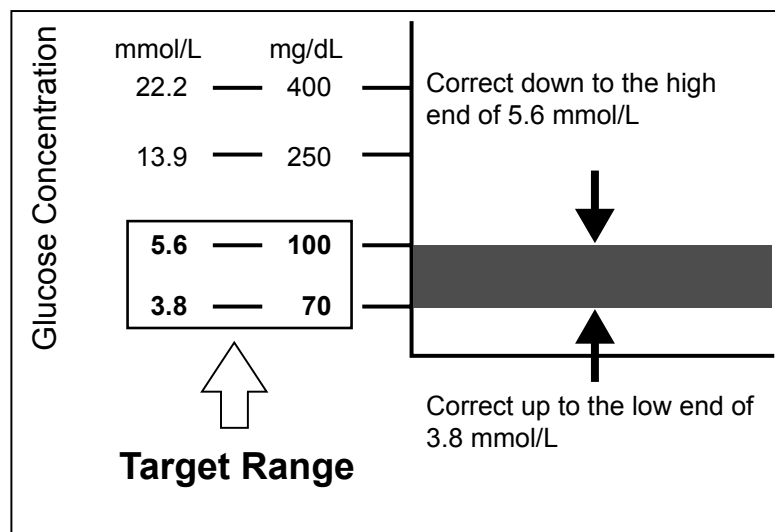

Pumps are sent from the factory with default BG Target range of 5.6-5.6 mmol/L (100-100 mg/dL).

- 1 Make sure the EDIT SETTINGS screen is open.

Main > Bolus > Bolus Setup > Bolus Wizard Setup > Edit Settings

- 2 Select **BG Target**, then press **ACT**.
- 3 The TARGET RANGE 1 screen appears. The low end of your BG Target range flashes on the screen.
- 4 Set the BG Target, then press **ACT**.
- 5 The high end of your BG Target range flashes on the screen. Enter the BG Target, then press **ACT**.
The start time for your first BG Target is midnight and cannot be changed.
- 6
 - a. If your BG Target is outside of 5.0-7.8 mmol/L (90-140 mg/dL), the pump screen displays a warning that the values are acceptable but outside normal range. Press **ESC** to change your BG Target or press **ACT** to set this range.
 - b. If your BG Target is within 5.0-7.8 mmol/L (90-140 mg/dL), the SET START TIME 2 screen appears. The dashes under the screen name flash. The first BG Target range is now set.
If you do not need a second BG Target range, press **ESC** and skip to the next section. If you need to set up another BG Target range, follow steps 7 through 11.
- 7 In the SET START TIME 2 screen enter the time of the day you want this BG Target range to become active.
- 8 Press **ACT**. The TARGET RANGE 2 screen appears.
- 9 The low end of your BG Target range flashes on the screen. Set the BG Target, then press **ACT**.
- 10 The high end of your BG Target range flashes on the screen. Enter the BG Target, then press **ACT**.
- 11 The SET START TIME 3 screen appears. The second BG Target range is now set.
- 12 If you do not need to set up any more BG Target ranges, press **ESC**. If you need to set up more BG Target ranges, repeat steps 7 through 11 above for each one.

You are now ready to set up the Active insulin time.

About active insulin

Active insulin is the bolus insulin that has already been delivered to your body, but has not yet been used. The pump considers your active insulin time setting in determining any active insulin still in your body from prior boluses. This may help prevent hypoglycemia caused by over-correcting for high blood glucose.

The pump shows the active insulin amount in the ESTIMATE DETAILS screen during the Bolus Wizard programming steps, the STATUS and the SET BOLUS screens. However, the active insulin amount is calculated differently in the ESTIMATE DETAILS screen and appears with an asterisk (*Active Insulin). The active insulin amount calculated in the Bolus wizard steps includes the insulin that has already been delivered and insulin that is going to be delivered by the active Square Bolus.

Your Paradigm pump is shipped from the factory with an active insulin time setting of six hours, which most closely matches the published scientific data. If your healthcare professional prescribes a different time for you, the active insulin time setting can be adjusted in the Bolus Wizard menu in one-hour increments from two to eight hours.

For more details about active insulin, see the *Bolus Wizard feature specifications* section in the *Pump specifications* chapter.

CAUTION: If you give yourself insulin by using a syringe, the Bolus Wizard feature will not be able to correctly determine the active insulin in your system. Consult with your healthcare professional on how long you need to wait after a manual injection before you can rely on the active insulin calculation of your Bolus Wizard feature.

Active insulin time

The active insulin time setting lets the pump know which active insulin time to use in calculating the amount of active insulin to subtract before estimating a bolus. Your healthcare professional should determine the active insulin time that is best for you.

To set the Active insulin time, do the following steps:

- 1 Make sure the EDIT SETTINGS screen is open.
Main > Bolus > Bolus Setup > Bolus Wizard Setup > Edit Settings
- 2 Select **Active Ins Time**, then press **ACT**.
- 3 The **ACTIVE INS TIME** screen appears. The default time of 6 hours flashes on the screen.
- 4 Set the number of hours for the active insulin time, then press **ACT**.
- 5 The EDIT SETTINGS screen shows the new Active insulin time setting. You have now completed the Bolus Wizard setup. Press **ESC** or wait to see the message: **Bolus Wizard setup is complete**.

Review your Bolus Wizard feature settings

Check your Bolus Wizard feature settings in the REVIEW SETTINGS screen. If necessary, compare this information with your information in the Bolus Wizard feature settings table.

- 1 Go to the REVIEW SETTINGS screen.
Main > Bolus > Bolus Setup > Bolus Wizard Setup > Review Settings
- 2 Scroll through the text to view your Bolus Wizard settings.
- 3 Exit the menus when you are done.

Meter option

You can set up your pump to automatically receive your blood glucose reading from the linked meter. This meter may not be available in all countries. Check with your local Medtronic Diabetes representative. Your pump is set at the factory with the meter option turned off. See the *Add, delete, review meter IDs* section in this chapter to turn the meter option on, and to enter your meter ID. Programming your meter ID links your pump to the meter. Each meter has its own unique ID. You can link up to three meters to your pump. If you do not link the meter to your pump, you must enter your blood glucose readings manually.



When the pump is idle (at the HOME screen), it will beep or vibrate when it receives a blood glucose reading from the meter. The reading will appear on the pump screen. If your BG reading is not within your set BG Target range, the pump beeps three times. For BG readings that are above the set BG Target range, use the Bolus Wizard feature to enter a correction bolus.

NOTE: *The use of RF (radio frequency) devices with the pump reduces pump battery life.*

You have to turn on the meter option to add, delete or review the meter ID(s) programmed in your pump. The meter ID is printed on the back of the meter. See the user guide that came with your meter for detailed information on how to use it.

Meter rules

If you want your pump to communicate with the meter, the following conditions must apply:

- 1 The meter option must be turned on and programmed. Refer to the instructions in this section.
- 2 Your pump must be within 1.2 meters (4 feet) of your meter to receive the blood glucose reading.
- 3 The pump cannot have a LOW BATTERY alert condition.
- 4 When programming a bolus, the blood glucose measurement from the meter will appear as the default blood glucose value on the ENTER BG screen. The pump will not display a reading that is older than 12 minutes on the ENTER BG screen.
- 5 Do not use the RF meter to send your blood glucose readings to the pump while onboard aircraft. Manually enter your blood glucose.

CAUTION: The pump will not receive signals from the meter while it has a LOW BATTERY condition. To ensure the meter communicates with the pump, make sure the pump does not have a low battery. (Replacing the low battery with a new battery will restore meter-pump communication.)

Add, delete, review meter IDs

The meter programming screens are very similar to those for the remote control. Make sure to select **Meters** in the DEVICE OPTIONS screen when programming your meter.

You have to turn on the meter option to add, delete, or review the meter ID(s) programmed in your pump.

To turn on the meter option:

- 1 Go to the METER OPTION screen.
Main > Utilities > Connect Devices > Meters
- 2 Select **On**, then press **ACT**. The METER ID MENU appears.
- 3 Add, delete or review your meter ID(s) as desired.
 Adding meter IDs
 - a. Select **Add ID**, then press **ACT**.
 - b. Use the up and down arrow buttons to enter each of the six ID characters. Press **ACT** after each entry.
 - c. After you set the last character of the ID, the screen will return to the METER ID MENU.
 Deleting meter IDs
 - a. Select **Delete ID**, then press **ACT**.
 - b. Select the meter ID that you want to delete and press **ACT**.
 - c. The selected ID is now deleted.
 Reviewing meter IDs

If you are not sure that your meter ID is entered in your pump, check the REVIEW METER ID screen.

 - a. Select **Review ID**, then press **ACT**.
 - b. The programmed IDs will show in the REVIEW METER ID screen.
- 4 Exit the menus when you are done.

Normal bolus using Bolus Wizard feature

After the Bolus Wizard feature is turned on and programmed, this feature can calculate an estimate of insulin you need for your correction bolus and/or your food bolus. You have the option of using the estimate or changing it as necessary. Additionally, your pump can receive your blood glucose reading from the meter, if they are linked.

Use the **⬅B** button to deliver a Normal bolus at any time except during another Normal bolus. A Normal bolus will temporarily interrupt a Square Wave or Dual Wave bolus that is delivering. After the Normal bolus is finished, the Square Wave or Dual Wave bolus delivery will resume.

NOTE: *If you want to use the pump-to-meter link, make sure the meter option is on. Refer to the Meter option section in this chapter for instructions.*

- 1 If you want a correction bolus, check your blood glucose with your blood glucose meter and go to step 2. If you want to bolus for food, go to step 2.
- 2 Press **⬅B** on your pump, or go to the BOLUS MENU, select **Use Bolus Wizard**, and press **ACT**.
- 3 The ENTER BG screen will appear.
 - a. If you are **not** using the meter:
Enter your blood glucose value. Press **ACT** and continue to step 4. If you are not entering a blood glucose and want to bolus for food, select the dashes in the ENTER BG screen. The Bolus Wizard feature will calculate the insulin needed for your food entry without considering your blood glucose level. Press **ACT** and continue to step 4.
 - b. If you are using the meter, you must program your bolus within 12 minutes of the pump receiving the reading from the meter. If more than 12 minutes have passed, the reading will no longer be available from the screen and you must enter your blood glucose manually.
Your pump checks if the entered BG is within your target range. Press **ACT** to accept the blood glucose value. You can also change this blood glucose value, if necessary, then press **ACT**.
- 4 The ENTER FOOD screen appears.
 - a. If this is a food bolus, enter the food value you will eat, then press **ACT**.
 - b. If this is a correction bolus, select 0 (zero) as the value, then press **ACT**.
- 5 The ESTIMATE DETAILS screen appears. Review the information on this screen. If you need to make any changes, press **ESC** to return to the ENTER BG screen (step 3) and make changes as necessary.
- 6 Press **ACT** in the ESTIMATE DETAILS screen. The SET BOLUS screen appears with the estimated bolus amount flashing. Change the amount if desired. Press **ACT** to accept and start delivering your bolus.

NOTE: *If you have BG Reminder turned On, the BG REMINDER DURATION screen displays. It allows you to set the duration before you are reminded to check your blood glucose after a bolus. See the BG Reminder section in the Basic programming chapter for information about this feature.*

- 7 The BOLUS DELIVERY screen appears. The pump will beep or vibrate at the start and end of the bolus. As the bolus is delivered, the screen shows the bolus type and amount until the total units have been delivered. The screen then returns to the HOME screen.

Bolus Wizard feature examples

For the scenarios that follow, Michael has his Bolus Wizard feature turned on with the following settings:

- Carb ratio: 15 grams per unit of insulin
- Insulin Sensitivity: 2.2 mmol/L (40 mg/dL) per unit of insulin
- BG Target: 5.0-6.6 mmol/L (90-120 mg/dL)
- Active Insulin Time: 6 hours

NOTE: If you want to see details of the formulas the Bolus Wizard feature uses to calculate estimate boluses like the ones in the following examples, see the Bolus Wizard feature specifications section in the Pump specifications chapter.

Example 1: Blood glucose on target (normal blood glucose) and no active insulin

Michael awakens in the morning before school and his mother has breakfast waiting for him. Before he begins eating, he tests his blood glucose with his meter and his blood glucose result of 6.6 mmol/L (120 mg/dL) is automatically sent to his pump.

He estimates that his meal consists of 60 grams of carbohydrates. When prompted by the Bolus Wizard feature, he enters this amount in the ENTER FOOD screen. Based on his Bolus Wizard feature settings, the pump will suggest that he take 4.0 units of insulin.

Food estimate: $\frac{60 \text{ g}}{15 \text{ g/u}} = 4 \text{ units}$

Correction estimate: Correction is 0 because the current blood glucose reading is within the BG Target range.

Bolus estimate = Food estimate + Correction estimate

Bolus estimate = 4 units + 0 = 4 units

Example 2: Blood glucose above target (high blood glucose) and no active insulin

The next day, Michael wakes up before school. Before eating the same breakfast, he tests his blood glucose with his meter and finds it to be 11.1 mmol/L (200 mg/dL), which is above his target of 6.6 mmol/L (120 mg/dL). His blood glucose reading is automatically sent to his pump.

When prompted by the Bolus Wizard feature, he enters his carbohydrate amount of 60 grams in the ENTER FOOD screen. Based on his settings, the pump will suggest that he take 6.0 units of insulin.

$$\text{Food estimate:} \quad \frac{60 \text{ g}}{15 \text{ g/u}} = 4 \text{ units}$$

$$\text{Correction estimate:} \quad \frac{11.1 \text{ mmol/L} - 6.6 \text{ mmol/L}}{2.2 \text{ mmol/L/u}} = 2 \text{ units}$$

or

$$\frac{200 \text{ mg/dL} - 120 \text{ mg/dL}}{40 \text{ mg/dL/u}} = 2 \text{ units}$$

Bolus estimate = Food estimate + Correction estimate

Bolus estimate = 4 units + 2 units = 6 units

Example 3: Blood glucose below target (low blood glucose) and no active insulin

On another morning, Michael sits down before eating the same breakfast. He tests his blood glucose with his meter and finds it at 3.9 mmol/L (70 mg/dL), which is below his Low BG Target of 5.0 mmol/L (90 mg/dL). His reading is automatically sent to his pump.

When prompted by the Bolus Wizard feature, he enters his carbohydrate amount of 60 grams in the ENTER FOOD screen. Based on his settings, the pump will suggest that he only take 3.5 unit of insulin.

$$\begin{array}{lcl} \text{Food estimate:} & \frac{60 \text{ g}}{15 \text{ g/u}} & = 4 \text{ units} \\ & & \\ & \frac{3.9 \text{ mmol/L} - 5.0 \text{ mmol/L}}{2.2 \text{ mmol/L/u}} & = -0.5 \text{ units} \\ \text{Correction estimate:} & \text{or} & \\ & \frac{70 \text{ mg/dL} - 90 \text{ mg/dL}}{40 \text{ mg/dL/u}} & = -0.5 \text{ units} \\ & & \\ \text{Bolus estimate} & = \text{Food estimate} + \text{Correction estimate} & \\ \text{Bolus estimate} & = 4 \text{ units} + (-0.5) & = 3.5 \text{ units} \end{array}$$

Example 4: Blood glucose above target (high blood glucose) with active insulin

Michael is at school and wants to eat a snack in the late morning. He tests his blood glucose with his meter and finds it to be at 11.1 mmol/L (200 mg/dL), which is above his target of 6.6 mmol/L (120 mg/dL). He estimates that his snack contains 60 grams of carbohydrate, so he enters 60 into the pump when prompted by the Bolus Wizard feature. Based on his settings, and as a result of 1.5 units of active insulin, his pump will suggest that he take 4.5 units.

$$\text{Food estimate:} \quad \frac{60 \text{ g}}{15 \text{ g/u}} = 4 \text{ units}$$

$$\text{Active insulin:} \quad 1.5 \text{ units}$$

$$\text{Correction estimate:} \quad \frac{11.1 \text{ mmol/L} - 6.6 \text{ mmol/L}}{2.2 \text{ mmol/L/u}} - \text{Active insulin} = 0.5 \text{ units}$$

or

$$\frac{200 \text{ mg/dL} - 120 \text{ mg/dL}}{40 \text{ mg/dL/u}} - \text{Active insulin} = 0.5 \text{ units}$$

$$\text{Bolus estimate} = \text{Food estimate} + \text{Correction estimate}$$

$$\text{Bolus estimate} = 4 \text{ units} + 0.5 \text{ units} = 4.5 \text{ units}$$

Example 5: Blood glucose below target (low blood glucose) with active insulin

Another day at school, Michael is getting ready to eat lunch. He tests his blood glucose with his meter and finds it at 3.9 mmol/L (70 mg/dL), which is below his Low BG Target of 5.0 mmol/L (90 mg/dL). His reading is automatically sent to his pump.

When prompted by the Bolus Wizard feature, he enters his carbohydrate amount of 60 grams in the ENTER FOOD screen. Based on his settings, and despite 1.5 units of active insulin, his pump will suggest that he take 3.5 units of insulin.

$$\text{Food estimate:} \quad \frac{60 \text{ g}}{15 \text{ g/u}} = 4 \text{ units}$$

Active insulin: 0* units

$$\text{Correction estimate:} \quad \frac{3.9 \text{ mmol/L} - 5.0 \text{ mmol/L}}{2.2 \text{ mmol/L/u}} = -0.5 \text{ units}$$

or

$$\frac{70 \text{ mg/dL} - 90 \text{ mg/dL}}{40 \text{ mg/dL/u}} = -0.5 \text{ units}$$

Bolus estimate = Food estimate + Correction estimate

Bolus estimate = 4 units + (-0.5) units = 3.5 units

*When the current blood glucose is below the Low BG Target, an active insulin amount is not considered in the Bolus Wizard feature calculations.

Optimizing pump therapy


Square Wave and Dual Wave bolus

Square Wave bolus delivers a bolus evenly over a period of time (30 minutes to 8 hours). This bolus can be used for insulin delivery when you have eaten a long meal with extended snacking. It can also be useful if you have delayed food digestion due to gastroparesis or meals high in fat. A Square Wave bolus can be useful if a Normal bolus drops your blood glucose too rapidly. Since the Square Wave portion extends over a period of time, the insulin is more likely to be available to match your individual needs.

NOTE: *During delivery of a Square Wave bolus, you will not be able to do the following pump functions: change the max bolus amount, change the scroll rate, disable or deliver Dual and Square Wave boluses, rewind or fill the cannula, change the active insulin time, run a self-test, or access the User Settings menu. All other pump functions are still available during the Square Wave bolus.*

Dual Wave bolus delivers a combination of an immediate Normal bolus followed by a Square Wave bolus. The Square Wave portion is delivered evenly over a period of time. A Dual Wave bolus is useful for meals with both rapidly and slowly absorbed carbohydrates. For example, a Dual Wave bolus would be appropriate for fruit and crackers followed by pasta. The Dual Wave option meets both immediate and extended insulin needs. A Dual Wave bolus is also useful for correcting elevated blood glucose before a meal.

See the following graphic for a description of the different bolus types:

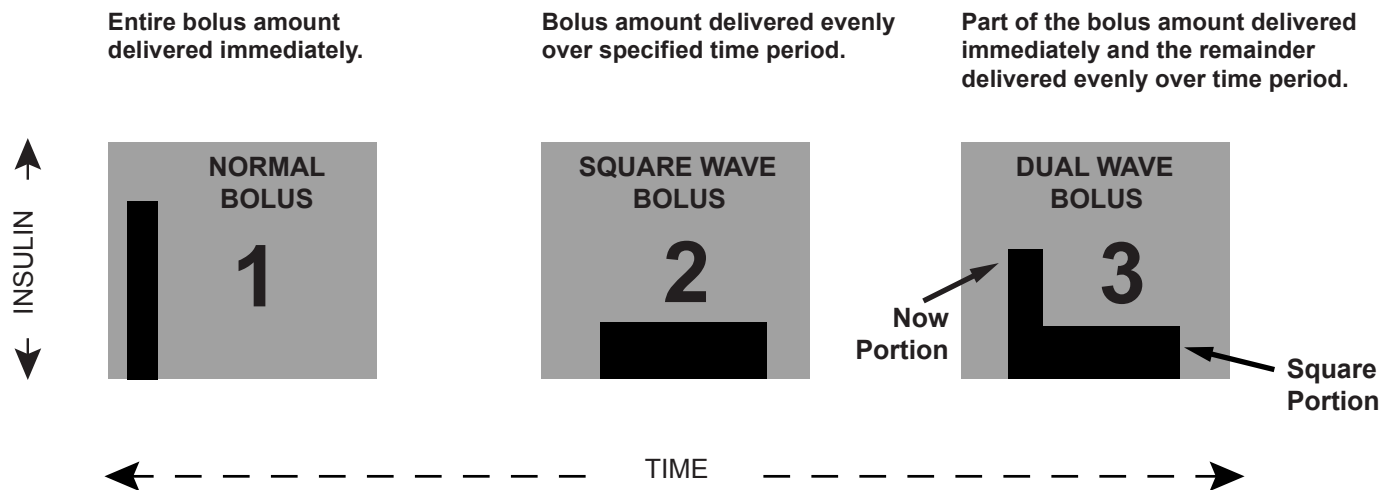

Turning on the Dual Wave/Square Wave option

It is important that you consult with your healthcare professional before using a Square Wave or Dual Wave bolus. You should be familiar with the basic functions of your pump before exploring these options.

To set up a Dual Wave or Square Wave bolus, you must first turn on the dual/square bolus option.

- 1 Go to the DUAL/SQUARE OPTION screen.
Main > Bolus > Bolus Setup > Dual/Square Bolus
- 2 Select **On**, then press **ACT**. The feature is now on. Exit the menus.

Square Wave or Dual Wave bolus without Bolus Wizard feature

- 1 Make sure the dual/square option is on.
- 2 Calculate your food and/or correction bolus amount.
- 3 Go to the BOLUS TYPE screen.



Press **⬅B** on your pump, or follow this path:

Main > Bolus > Set Bolus

4 For a Square Wave bolus, do these steps:

- a. Select **Square Wave Bolus**, then press **ACT**. The SET SQUARE BOLUS screen appears.
- b. Enter the desired amount for the Square Wave bolus units, then press **ACT**.
- c. Continue to step 5.

For a Dual Wave bolus, do these steps:

- a. Select **Dual Wave Bolus**, then press **ACT**. The SET DUAL BOLUS TOTAL screen appears.
 - b. Enter the desired amount for the total dual bolus units. This amount is the total of Normal and Square Wave bolus units. Press **ACT**.
 - c. The next screen flashes the amounts of Now (Normal) and Square Wave portions of the Dual Wave bolus. The screen also shows the percentage amount of each portion. Press  or  to change the percentage/number of units. Continue to step 5.
- 5** The SQUARE DURATION screen appears. Enter the amount of time you want the Square Wave bolus to last, then press **ACT**.

NOTE: *If you have BG Reminder turned On, the BG REMINDER DURATION screen displays. It allows you to set the duration before you are reminded to check your blood glucose after a bolus. See the BG Reminder section in the Basic programming chapter for information about this feature.*

- 6** The BOLUS DELIVERY screen appears with an open circle indicating that your pump is in Special mode. The pump beeps/vibrates at the start of the bolus. During bolus delivery, the pump will return to the HOME screen. The pump beeps/vibrates at the end of the bolus and the open circle disappears.

Square Wave bolus practice

Your target pre-meal blood glucose range is _____ to _____.

Check your pre-meal blood glucose. Are you within your target? _____. If yes, continue. If no, wait to try the following test until your pre-meal blood glucose is within your target range:

TEST: Choose a meal that is high in fat (hot dogs, pizza, cheese enchiladas). Determine your meal bolus amount. Set the Square Wave bolus to deliver the determined amount of insulin over two hours. (This duration time is an example. As always, consult with your healthcare professional for guidance.)

Check your blood glucose and record:

Pre-meal _____

1 hour post meal _____

2 hours post meal _____

3 hours post meal _____

4 hours post meal _____

Did your blood glucose return to your pre-meal target within 4 hours post meal? _____

If yes, then repeat this test with the same meal on another day to verify your results.

If no, discuss this with your healthcare professional for guidance.

Dual Wave bolus practice

Can you think of any meals where this feature would help you with blood glucose control?

Your target pre-meal blood glucose range is _____ to _____

Check your pre-meal blood glucose. Are you within your target? _____ If yes, continue. If no, try this test when your pre-meal blood glucose is within your target range:

TEST: Choose a meal that has a combination of both rapidly absorbed and slowly absorbed carbohydrates. Determine your meal bolus amount. Set the Dual Wave bolus to deliver the determined amount of insulin. Program your pump to deliver one-half over 2-hours*, and the other half immediately.

(* This duration of time and ratio is an example. As always, consult with your healthcare professional for guidance.)

Check blood glucose and record:

Pre-meal_____

1 hour post meal_____

2 hours post meal_____

3 hours post meal_____

4 hours post meal_____

Did your blood glucose return to your pre-meal target within 4 hours post meal? _____*

If yes, repeat this test with the same meal on another day to verify results.

If no, discuss this with your healthcare professional for guidance.

Using the Bolus Wizard feature for a Square Wave or Dual Wave bolus

If you are using the Bolus Wizard feature to calculate your Square Wave or Dual Wave bolus amounts, you will be prompted to enter your blood glucose reading and/or the (gram or exchange) units you will eat. The Bolus Wizard feature will use this input to calculate your suggested correction and/or food bolus amount. If you do not want to use the Bolus Wizard feature estimate, you can change it.

The Bolus Wizard feature must be turned on and the settings must be programmed (see the *How to program the Bolus Wizard feature* section in the *Using the Bolus Wizard* chapter). Also, make sure the dual/square option is turned on (see the *Turning on the Dual Wave/Square Wave option* section in this chapter).

If you want to use the pump-to-meter link, make sure the meter option is on. Refer to the *Meter option* section in the *Using the Bolus Wizard feature* chapter for instructions.

- 1 Go to the ENTER BG screen.

Press  on your pump, or use the following path:

Main > Bolus > Use Bolus Wizard



- 2 Enter your blood glucose value, then press **ACT**.
- 3 The ENTER FOOD screen appears. Enter your food, then press **ACT**.
- 4 The ESTIMATE DETAILS screen appears. Scroll down to review the information there. Press **ACT** to continue to step 5.

If you need to make any changes, press **ESC** to return to the ENTER BG screen. Make changes as necessary.

- 5 The BOLUS EST screen appears with Normal Bolus, Square Wave Bolus, and Dual Wave Bolus options. If your Bolus Wizard feature calculates that your bolus includes a portion to correct your high blood glucose, the Square Wave Bolus option will not be available. This helps you to select a bolus type (Normal or Dual Wave) that has an immediate delivery option to cover your high blood glucose.
- 6 **To set a Square Wave bolus, do these steps:**

- a. In the BOLUS EST screen, select **Square Wave Bolus**, then press **ACT**.
- b. The SET SQUARE BOLUS screen flashes the estimated bolus amount. Change the amount if needed. Press **ACT** to accept the bolus amount.

To set a Dual Wave bolus, do these steps:

- a. The SET DUAL BOLUS TOTAL screen flashes the estimated bolus amount. This amount is the total of both the Normal and Square Wave bolus units. Change the amount if needed. Press **ACT** to accept the bolus amount.
- b. The next screen flashes the amounts of Now (Normal) and Square portions of the Dual Wave bolus. The screen also shows the percentage amount of each portion. Press **ACT** to accept the Bolus Wizard's suggested portions. You can also press  or  to change these portions, then press **ACT**.


The Bolus Wizard feature recommends splitting the food portion of your bolus 50/50 between the Square and Now portions. The entire correction amount is always recommended to the Now portion. In this example the NOW portion consists of half of the food insulin plus the correction amount less the active insulin (1.5U + 2.5U - 1.5U). This gives 2.5U or 62% of total insulin of 4.0U. The Square portion consists of the other half of the food insulin (1.5U), which is 38% of total insulin of 4.0U.


- 7 The SQUARE DURATION screen appears. Enter the amount of time you want the Square Wave bolus to last, then press **ACT**.


NOTE: *If you have BG Reminder turned On, the BG REMINDER DURATION screen displays. It allows you to set the duration before you are reminded to check your blood glucose after a bolus. See the BG Reminder section in the Basic programming chapter for information about this feature.*

- 8 Press **ACT** to accept and deliver the bolus. The BOLUS DELIVERY screen appears with an open circle indicating that your pump is in Special Mode. The pump beeps or vibrates at the start of the bolus. During bolus delivery, the pump returns to the HOME screen. If you want to see the progress of the delivery, press **ESC** to see the STATUS screen. The pump beeps or vibrates at the end of the bolus, and the open circle disappears.

Easy bolus

The **EASY BOLUS** button  allows a quick way to deliver a Normal bolus. You will preset the settings for this feature in the EASY BOLUS OPTION screen in the BOLUS MENU. Your pump is sent from the factory with the Easy Bolus feature set to off. If you want to use Easy Bolus, turn it on.

After you set up Easy Bolus, with each  press, you can increase the Normal bolus amount by a fixed amount, called a step. Before you can deliver an Easy Bolus, you must set the amount in the EASY BOLUS ENTRY screen. This amount equals the number of units of insulin for each step. The maximum number of steps can equal up to your maximum bolus limit. When using vibrate mode, EASY BOLUS is limited to 20 steps or maximum bolus, whichever comes first.

After you set your step amount, you can program your Easy Bolus. When you are in the HOME screen, each time you press  the Easy Bolus amount increases by one step. You will hear a beep or feel a vibration for each step increase. Each beep is a different tone. This makes it easy for you to count the beeps while you are programming your Easy Bolus.

Easy bolus setup


- 1 Go to the EASY BOLUS OPTION screen.
Main > Bolus > Bolus Setup > Easy Bolus

Select **On/Set**, then press **ACT**.

NOTE: *If you are using the remote control, the Easy Bolus must be set to on.*

Step value setup

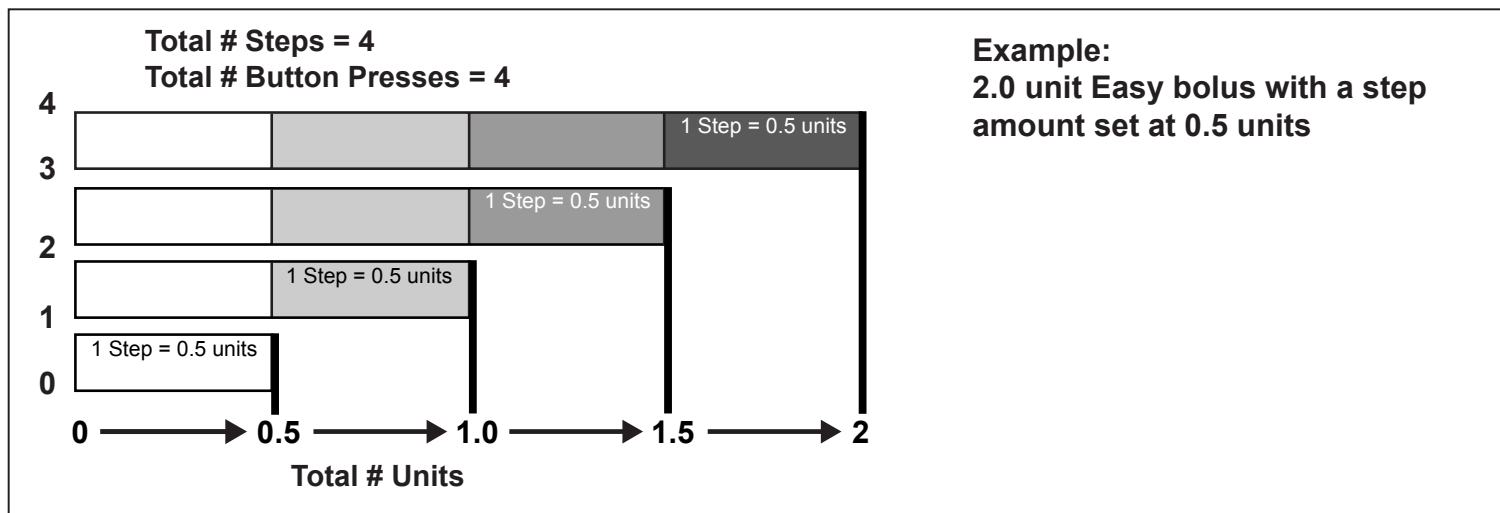


You can set the step value from 0.1 to 2.0 units, or to the set maximum bolus limit of less than 2.0 units (factory setting: 0.1). Set the step to a number that is convenient to use and easy to multiply.





- 1 The EASY BOLUS ENTRY screen flashes the step value. Change the value, then press **ACT**. The step value is the increment you will use for your Easy bolus.
- 2 The screen will return to the BOLUS MENU. Your step amount is now programmed and Easy Bolus is ready to use. Exit the menus.

Delivering Easy Bolus


Practice using the Easy Bolus feature while looking at the pump screen as you count the beeps. This feature only works from the HOME screen. After you are familiar with Easy Bolus, you can use the audible tones for bolus delivery without having to look at the screen.

- 1 From the HOME screen, press . The SET EASY BOLUS screen appears with the single step value flashing.
- 2 Press  the number of times needed for your bolus amount. Watch the amount change on the screen with each press. The pump will vibrate or sound a different tone for each  press.


- 3 When your total bolus amount appears on the SET EASY BOLUS screen, press **ACT**. Listen/feel to count the steps without looking at the screen.


For example: You need to deliver a 2.0 unit bolus, and your step size is set to 0.5. Each time you press , the units increase by the step amount of 0.5 units. To deliver 2.0 units, you need to press the  button four times. ($4 \times 0.5 = 2.0$) The screen will show 2.0 units.



NOTE: For your safety, you cannot use the  to select the Easy Bolus values. Pressing  or **ESC** will cancel the Easy bolus.

- 4 If this amount is correct, press **ACT** to start the Easy Bolus delivery. The BOLUS DELIVERY screen will show the units being delivered. When the total bolus is finished, the pump will beep or vibrate. If this amount is wrong, press **ESC** or  to start over. The pump will return to the HOME screen.

Example 1: Easy bolus

Alexander is a busy executive with an accounting firm. He wears his Paradigm pump on his belt and does not want to take it off to give himself a bolus. Alexander can easily reach down and feel for the Easy Bolus button  to give a bolus.

He previously programmed his pump to deliver an Easy Bolus in steps of 0.5 unit increments. From the HOME screen, with each press of the , the pump will sound a different tone so he can keep track of the number of button presses.

He wants to give himself 2.0 units for a snack, so he will press  4 times ($4 \text{ presses} \times 0.5 \text{ units/press} = 2.0 \text{ units}$) and then press the **ACT** button. The pump counts back 4 beeps because he pressed  4 times. He simply presses **ACT** to confirm the amount, and his pump delivers the 2.0 units.

When Alexander wants to be more discreet, or does not want his pump to beep in an important meeting, he can set the pump to Vibrate mode and feel for vibrations rather than listening for the tones. (See the *Setting your alert type* section in the *Utilities* chapter for details.)

Your turn: Easy bolus practice

The factory default setting for the Easy Bolus feature is 0.1 unit steps. You can change the step level as necessary to a value that is more convenient for you to use and easier to multiply.

Give your next bolus by using the Easy Bolus feature on your pump.

Your step level is _____.

How many units did you give? _____

How many tones did you count? _____

It might be a good idea to look at your pump's screen to see the bolus amount as well as counting the steps the first few times you try this until you become familiar and comfortable with the feature.

Missed Bolus Reminder

Missed Bolus Reminder feature helps you remember to give yourself a meal bolus. It triggers a MISSED BOLUS alert if you do not deliver a bolus during the time periods that you set. Set these time periods for the meals for which you want to receive a reminder. The pump is set at the factory with the Missed Bolus Reminder feature turned Off.

Add, delete, and review bolus reminders

You have to turn on the Missed Bolus Reminder option to add, delete and review the programmed reminders.

- 1 Go to the MISSED BOLUS REMINDER screen.
Main > Bolus > Bolus Setup > Missed Bolus Reminder
- 2 Select **On/Set**, then press **ACT**.

Add

You can set up to four Missed Bolus Reminders.

- 1 Select **Add Reminder** on the MISSED BOLUS REMINDER screen, then press **ACT**.
- 2 The Start Time flashes. Times can be entered in 30 minute increments. Set the Start Time, then press **ACT**.
- 3 After this time is set, the End Time flashes. Set the End Time, then press **ACT**.
If you enter the same Start Time and End Time, you will get one Missed Bolus Reminder in a 24 hour period. If you add a Missed Bolus Reminder with a start time earlier than the current pump time, you will not get this Missed Bolus Reminder until the next day.

Delete

- 1 Select **Delete Reminder** on the MISSED BOLUS REMINDER screen and press **ACT**.
- 2 In the DELETE REMINDER screen, highlight the bolus reminder that you want to delete and press **ACT**.

Review

- 1 Select **Review Reminders** on the MISSED BOLUS REMINDER screen and press **ACT**.
- 2 Review your Missed Bolus reminders.

Basal patterns

The Basal Patterns feature is optional for pump users. You can set your pump to deliver a standard basal and two additional basal patterns to meet your individual daily, weekly, or monthly needs. Keep a paper copy of your programmed patterns with you at all times in case you need to reprogram your pump. To select and use pattern A or pattern B, the patterns option must be turned on and programmed.

NOTE: *You may want to explore this option after you become familiar with the basic pump functions. It is important that you consult your healthcare professional before using a pattern other than your standard pattern.*

- Standard pattern: Your normal basal that supports your usual day-to-day activity. When the Patterns feature is off, the pump uses your standard basal pattern.
- Pattern A/B: Basal pattern that supports activity levels that are not a part of your day-to-day routine, but are normal in your lifestyle. Such activities could be a sport that you do once a week or a change in your sleep pattern over the weekend, extended periods of higher or lower activity, or menses.

Turning on the patterns

Your pump is set at the factory with the basal patterns feature turned off. After you turn on patterns, you still have to program and select a pattern (A or B), as described in the next sections, before the patterns feature is active. If you turn off the patterns feature, your pump will automatically select your standard basal pattern.

- 1 Go to the PATTERNS OPTION screen.

Main > Basal > Basal Setup > Patterns

- 2 Select **On**, then press **ACT**. The patterns feature is now on. Exit the menus.

Programming a pattern

Your pump will keep your pattern settings even when the Patterns option is turned off. However, the patterns feature must be on to program a basal pattern.

To program your patterns, do these steps:

- 1 Go to the EDIT BASAL screen.

Main > Basal > Basal Setup > Set/Edit Basal

- 2 Select the basal pattern you want to program, then press **ACT**.
- 3 The SET BASAL RATE 1 screen appears. The basal rate flashes, indicating that it can be changed. Set your first rate, then press **ACT**.
The start time for your first basal rate is midnight and cannot be changed.
- 4 The SET START TIME 2 screen appears. The dashes under the screen name flash. The first basal rate is now set.
If you want to use the same rate for the whole day, press **ESC** and skip to step 10. If you want to program another basal rate, follow steps 5 through 8.
- 5 In the SET START TIME 2 screen enter the time of the day you want this basal rate to become active.
- 6 Press **ACT**. The SET BASAL RATE 2 screen appears.
- 7 The previously set basal rate or dashes flash on the screen. Select the value for this basal rate.
- 8 Press **ACT**. The SET START TIME 3 screen appears. The second basal rate is now set.

- 9 If you do not need to set up any more basal rates, press **ESC**. If you need to set up more basal rates, repeat steps 5 through 8 above for each rate.
- 10 After you press **ESC**, the BASAL RATE screen appears. The screen will show:
 - the current basal pattern and basal rate,
 - time it started, and
 - the 24-hour basal total.

After you make changes to a pattern, the pump will use that pattern as the current basal. Make sure the basal you want is selected in the SELECT PATTERN screen.

Select a pattern

Before you select a pattern to be active, make sure the Patterns feature is turned on. After your standard pattern and/or pattern A or B are set, do these steps to select a pattern to be the active one:

- 1 Go to the SELECT PATTERN screen.
Main > Basal > Select Patterns
- 2 Select the desired pattern, then press **ACT**.
- 3 The screen will return to the BASAL MENU. Your basal pattern is now active. Exit the menus.

NOTE: If pattern A or B is active, the pump is in Special mode. An open circle appears at the top of the screen.

Example 1: Basal patterns

Ken has had his insulin pump for about a month. He tests his blood glucose 4 - 6 times a day and records his results in his logbook. He is happy with his glucose control during the week but on the weekends, he noticed that he has to eat more food to prevent his blood glucose from running too low.

Ken has realized that during the week while he is at work, he is very inactive and sits at a desk most of the time. On the weekends, though, he is busy with yard work, running errands and playing with his kids. He determines that he needs to have lower basal settings to receive less insulin during active times, such as his weekend.

He can use the Basal Patterns feature to support his weekend change in activity. During the week, he can set his pump to deliver in the standard setting, and on Saturday morning, he can switch over to Pattern A, which he can set with lower basal rates for the weekend. On Monday morning, he can return his pump to the Standard setting for his weekday insulin needs.

Example 2: Basal patterns

Cynthia has had diabetes for about 12 years and has been on her Paradigm pump for several weeks. Every Monday, Wednesday and Friday, Cynthia goes on a 3 kilometer walk in the morning. To prevent hypoglycemia on these days, she uses the patterns feature. For those days, she simply switches over to Pattern A, which she has programmed with a lower set of basal rates. Before she learned to use the patterns feature, she would have to eat more food throughout the day to keep her blood glucose at a safe level. Cynthia has also noticed that a few days prior to menstruation, her blood glucose levels seem to rise, requiring more insulin. She has programmed Pattern B on her Paradigm pump with higher basal rates for this time. For her usual schedule, she uses the standard basal pattern.

Your turn:

Can you think of situations where you might require different basal rate settings on different days?

Temp basal rates

Temp basal should be used based on the guidance of your healthcare professional. This feature is useful to manage blood glucose levels during **unusual** short-term activities or conditions. These conditions could be an illness or physical activity.

A temporary basal rate allows an immediate short-term change to your basal insulin for a specified period of time (30 minutes to 24 hours). This rate can be up to your Maximum basal rate setting. It offers an easy way to immediately meet short-term insulin needs for temporary activities or situations. When your blood glucose is temporarily high or low, a temp basal rate allows you to set a temporarily higher or lower basal to accommodate your blood glucose. For ongoing periods of increased or decreased activity, the patterns feature may be more suitable.

How does temp basal work?

During a temp basal delivery, all other basal programming is temporarily overridden. After the temp basal delivery is completed, your pump will return to the programmed basal. A temp basal is delivered only once and does not repeat. If you want another temp basal, you must program the temp basal again. This feature may be useful to temporarily increase or decrease basal insulin during illness, exercise or similar unusual situations.


Temp basal types

Based on your preference, you can select either Percent of Basal or Insulin Rate.

Insulin rate (U/H)

Insulin rate is a fixed basal in units per hour (U/H). This temp basal type is independent of your current basal. When you select **Insulin rate (U/H)** for your temp basal type, your pump will deliver the fixed amount you have set for the duration as set. The amount of your temp basal insulin rate can be set up to your Maximum basal rate setting.

If you make changes to your normal basal rate, your U/H temp basal is not affected and will continue to deliver as programmed.

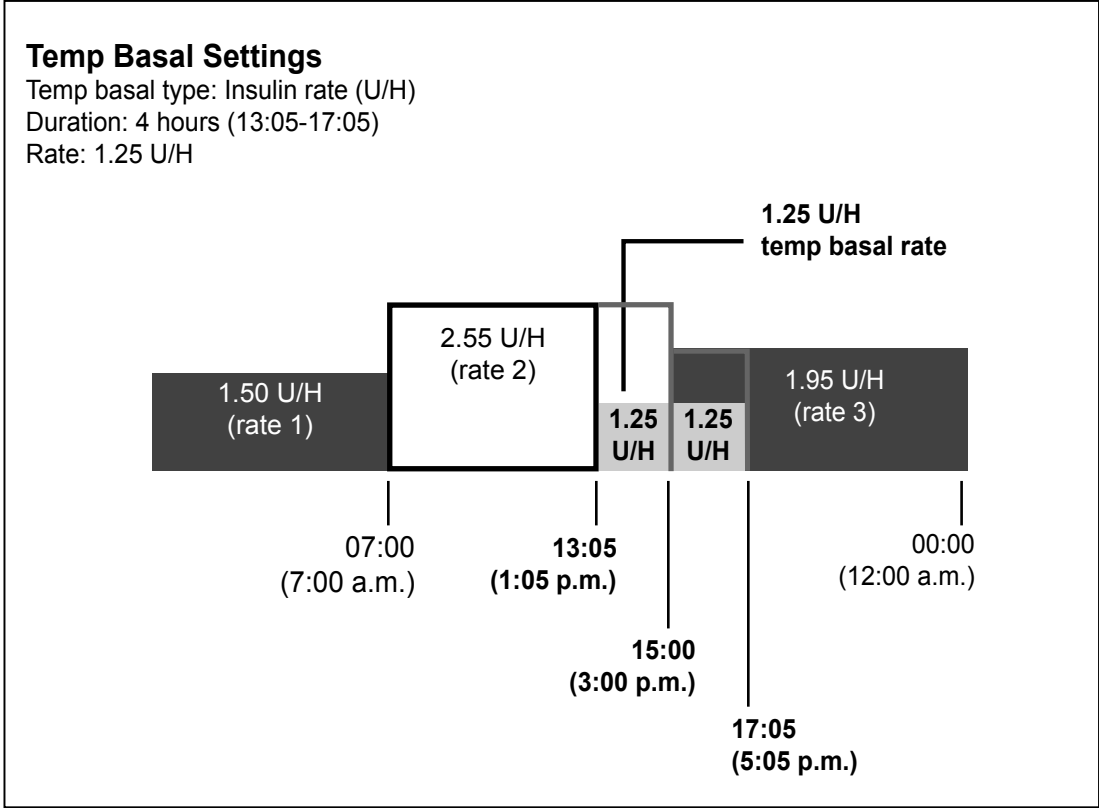

Percent of basal

Unlike insulin rate, the percent of basal is dependent on your current basal rates. Percent temp basal is a percentage of your current basal (0-200 percent limited to your Maximum basal rate setting). The temp basal amount is rounded down to the next 0.025 or 0.05 U/H increment based on the basal rate.

The maximum percent limit is based on the largest basal rate segment with your programmed temp basal rate time.

For example: It is 06:00 and your current basal rate is 1.50 U/H. You want to set a temp basal rate of 130 percent for seven hours. The maximum percent temp basal rate you can set is 125 percent. Anything larger would make #2 segment exceed your Maximum basal setting of 2.0 U/H.

Your current basal rates:	Your Maximum basal rate setting: 2.0 U/H
Segment #1: 00:00	1.50 U/H
Segment #2: 11:00	1.60 U/H (largest)
Segment #3: 16:00	1.30 U/H

If your current basal changes (for example, from rate 1 to rate 2), your percent temp basal amount will also change. The pump will deliver the percentage for the duration that you have set.


You cannot make changes to your normal basal rate while a percent temp basal is active. You must either wait until the temp basal is finished or cancel the temp basal in order to reprogram your normal basal rate setting(s).

Temp Basal Settings

Temp basal type: Percent of basal

Duration: 4 hours (13:05-17:05)

Rate: 50%

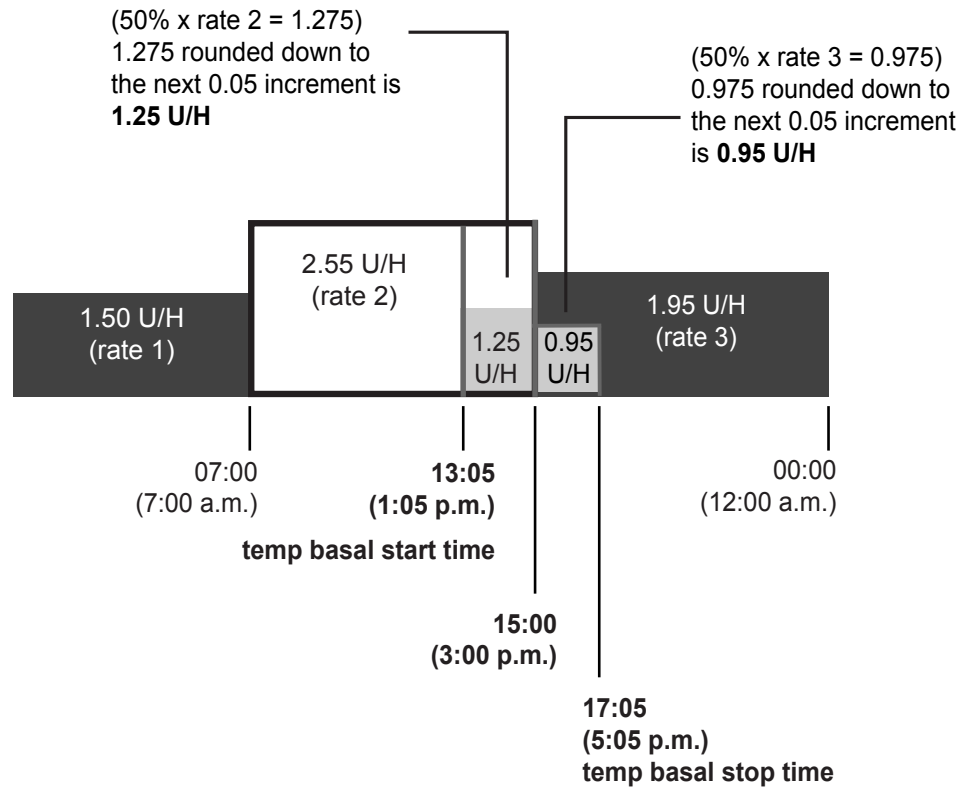

Selecting temp basal type

Your pump will remember the temp basal type setting. Once the type is set, you do not have to set it again. To select a temp basal type, do these steps:

- 1 Go to the SET TEMP BASAL AS screen.
Main > Basal > Basal Setup > Temp Basal Type
- 2 The SET TEMP BASAL AS screen appears. Select **Insulin Rate (U/H)** or **Percent of Basal**, then press **ACT**.
- 3 The screen will return to the BASAL SETUP screen. The temp basal type is now set. Exit the menus.
If you select **Percent of Basal** as your temp basal type, changes to your basal rate are not allowed until after temp basal is completed or cancelled.

Delivering a temp basal


A temp basal cannot exceed your programmed Max basal rate.

- 1 Go to the BASAL MENU.
Main > Basal
- 2 Select **Set/Edit Temp Basal**, then press **ACT**.
- 3 The SET DURATION screen appears. The duration will flash. Duration is the amount of time it will take for the pump to deliver the temporary basal. Enter the desired minutes or hours (30 minutes to 24 hours), then press **ACT**.
- 4 If you have selected insulin rate as your temp basal type, the SET TEMP BASAL U/H screen appears. If you have selected percent of basal as your temp basal type, the SET TEMP BASAL % screen appears. The temporary basal rate will flash. Enter your temp basal rate, then press **ACT**.
- 5 The BASAL MENU screen appears with an open circle at the top of the screen. Your pump is in Special mode; temp basal is now set and delivering. Exit the menus.

Verifying temp basal delivery

Temporary basal information is available in the STATUS screen only.

During a temporary (temp) basal, the pump is in Special mode (an open circle appears). This open circle will remind you that a temp basal is active. Additionally, your pump will beep/vibrate three times every hour during delivery. During delivery, go to the STATUS screen to see the current temp basal information.

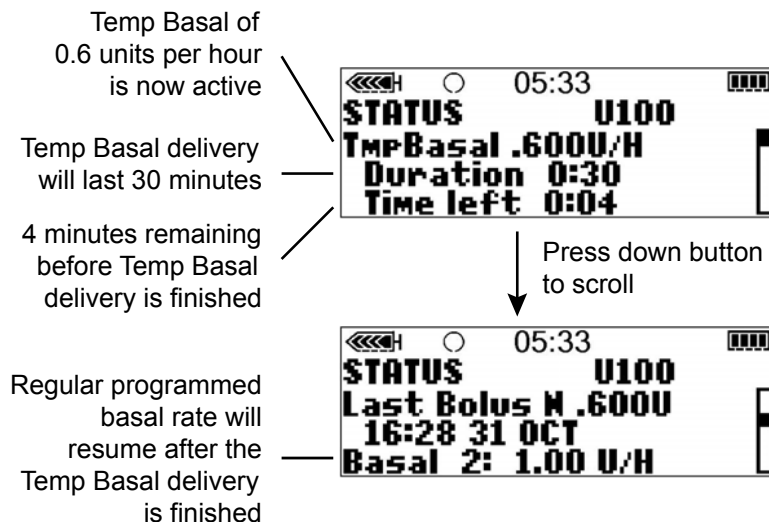

Canceling a temp basal

Use the Cancel Temp Basal function in the BASAL MENU to cancel a temporary basal. This function immediately stops the temp basal and resumes the regular programmed basal delivery.

To cancel a temp basal, do these steps:

- 1 Go to the BASAL MENU.
Main > Basal
- 2 Select **Cancel Temp Basal**, then press **ACT**.

The screen will return to the BASAL MENU. Your temp basal is cancelled and open circle has disappeared from the top of the screen. The programmed basal is now active again. Exit the menus.

Example 1:**Temp Basal for a decreased temporary basal rate**

Ramon and his friends got together for an unplanned game of soccer. Before using the pump, he was taking shots to manage his diabetes. Ramon experienced frequent low blood glucose reactions sometimes during, and very often after, he played games with his friends. Now that he is using his Paradigm pump, he can use the Temporary Basal Rate feature to help prevent low blood glucose. He simply programs his pump to temporarily deliver less basal insulin during the time that he is playing, and often for several hours after play, as well.

Ramon was able to determine how to set his Temporary Basal rates by frequent blood glucose testing, both during and after activity, and recording his results. The first time he tried using the pump, his healthcare professional advised him to program his pump to deliver one-half of his usual basal rate for the amount of time that he was playing and for an hour after he was done. He made small adjustments of the temporary basal rate and the duration of time, each time he tried to use the feature. After several different attempts with similar activity for the same amount of time, (such as his soccer game that lasted two hours), he found a temporary basal rate that worked well for him.

Example 2:**Temp Basal for an increased temporary basal rate**

Gail has had a cold with a cough for a couple of days. Because she is not feeling well, she tests her blood glucose more frequently. She finds that her blood glucose levels are running above target range before meals and she has needed several correction boluses to keep her blood glucose levels within her normal limits. Gail decides to use the Temporary Basal Rate to increase her basal rate during the day today. As advised by her healthcare professional, she will continue to check her blood glucose more frequently until she is feeling well.

Your turn:

Think of an activity where you might need to use a Temporary Basal Rate.

At what rate is your current basal rate running? _____

What Temporary Basal Rate would you try using at this time? _____

How long will you be active? _____

What duration will you set for the Temporary Basal Rate? _____

Test your blood glucose before and during activity and several times after as well. What are your blood glucose results?

Pre-activity _____

During activity _____

1 hour after activity _____



Several hours after activity _____

What Temporary Basal Rate changes will you make for the next time you try this?

Sensor features

The optional sensor and transmitter can provide continuous glucose measurements to help you control your glucose levels better. The sensor measures the glucose levels in the fluid under your skin. The transmitter receives this measurement from the sensor and sends it wirelessly to the pump. To take advantage of your pump sensor features, contact your local representative to purchase your sensor and transmitter.

Entering your sensor settings

Enter your pump sensor settings in the order they appear below, as some of the settings are dependent upon other settings being made first. When you complete each setting, your pump will automatically display the screen for the next setting in the sequence. When a screen item flashes, you may change the value of the flashing item by pressing  or .

Turning on the sensor

The sensor must be turned on, started and initialized to report glucose measurements.

- 1 Go to the EDIT SETTINGS screen.

Main > Sensor > Edit Settings

The EDIT SETTINGS screen appears with **Sensor: Off** selected.

- 2 Press **ACT**. The SENSOR ON/OFF screen appears.
- 3 Select **On**, then press **ACT**. The EDIT SETTINGS screen appears showing the sensor is now turned on. You are now ready to turn on the glucose alerts.

Turning on the Glucose Alerts

The Glucose alerts must be turned on if you want the system to send you an alert when the sensor glucose measurements reach or exceed your Glucose Limits. For a Low sensor glucose alert, the pump plays four consecutive tones, in falling pitch, if an audible beep has been selected as the alert type. For a High sensor glucose alert, the pump plays four consecutive tones, in rising pitch, if an audible beep has been selected as the alert type.

- 1 Make sure the EDIT SETTINGS screen is open.
Main > Sensor > Edit Settings
- 2 Select **Glucose Alerts**, then press **ACT**. The GLUCOSE ALERTS screen appears.
- 3 Select **On**, then press **ACT**. The EDIT SETTINGS screen appears showing the glucose alerts are now turned On.

Setting the Glucose Limits

If you turned the Glucose alerts On, then you need to set the high and low Glucose Limits recommended by your healthcare professional. Your recommended glucose limits may vary throughout the day, so your pump allows you to set up to eight pairs.

A Glucose Limit pair includes one High Glucose Limit and one Low Glucose Limit, as shown in the following table:

Glucose Limit pair	Alerts
Low: 4.4 mmol/L (80 mg/dL)	A Low Glucose Limit alert occurs when the sensor glucose measurement reaches or goes below the Low Glucose Limit in this table.
High: 13.3 mmol/L (240 mg/dL)	A High Glucose Limit alert occurs when the sensor glucose measurement reaches or goes above the High Glucose Limit in this table.


Selecting the BG units

You must select the blood glucose units (BG Units) for the system (either mmol/L or mg/dL). All BG measurements will show in the BG unit type you select.

- 1 Make sure the EDIT SETTINGS screen is open.
Main > Sensor > Edit Settings
- 2 Select **BG Units**, then press **ACT**. The BLOOD GLUCOSE UNITS screen shows **mg/dL** and **mmol/L**.
- 3 Select your BG units.
- 4 Press **ACT**. The EDIT SETTINGS screen shows the BG Units type you selected.
You are now ready to set up your glucose limits.

Setting the Glucose Limits start time

If your healthcare professional recommends that you set more than one pair of Glucose Limits, you must enter a start time for each. For example, you may use one pair of Glucose Limits during the day and another pair at night. The example below shows two pairs of Glucose Limits. The first pair begins at midnight, and the second pair begins at 9:00.

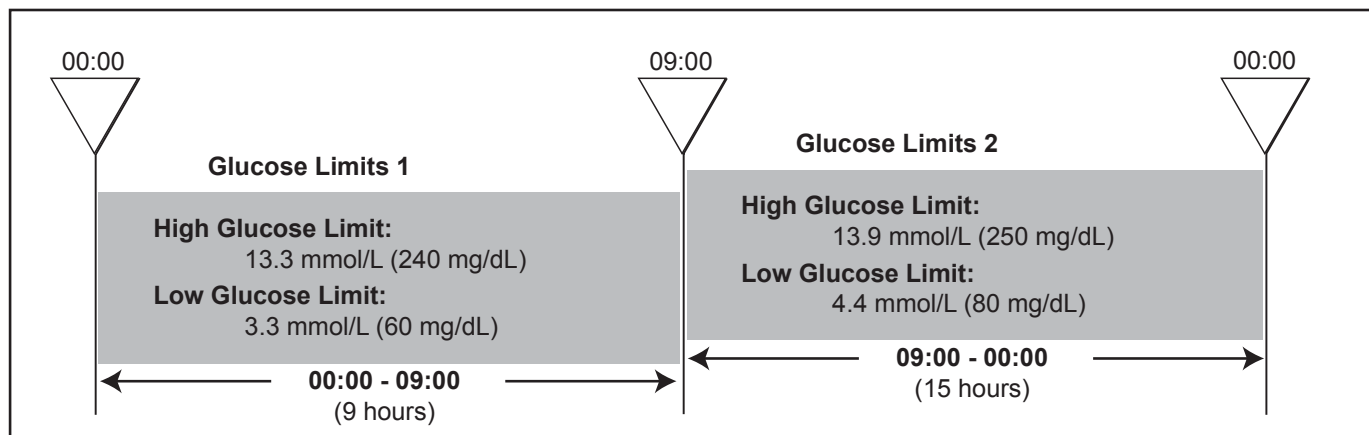

To set up your Glucose Limits:

- 1 Make sure the EDIT SETTINGS screen is open.

Main > Sensor > Edit Settings

- 2 Select **Glucose Limits**, then press **ACT**. The SET GLUCOSE LIMITS 1 screen appears. The Low Glucose Limit flashes 4.4 mmol/L (80 mg/dL).
The start time for your first pair of Glucose Limits is midnight and cannot be changed.
- 3 Select your Low Glucose Limit. The value must be between 2.2 and 21.6 mmol/L (40 and 390 mg/dL). To switch off the Low Glucose Limit, reduce it below 2.2 mmol/L (40 mg/dL) to select **OFF**.
- 4 Press **ACT** to select your High Glucose Limit. The screen flashes 13.3 mmol/L (240 mg/dL). The limit you select must be between 2.8 and 22.2 mmol/L (50 and 400 mg/dL). Your High Glucose Limit must also be at least 0.6 mmol/L (10 mg/dL) above your Low Glucose Limit. To switch off the High Glucose Limit, reduce the High Glucose Limit below 2.8 mmol/L (50 mg/dL) to select **OFF**.
- 5 Press **ACT**. The SET START TIME 2 screen appears. The dashes under the screen name flash. The first pair of Glucose Limits is now set.

If you do not need a second pair of Glucose Limits, press **ESC** and skip to the next section. If you want to set up a second pair of Glucose Limits, follow steps 6 through 11 below.

- 6 To set a second pair of Glucose Limits, select the time of day you want this pair to become active.
- 7 Press **ACT**. The SET GLUCOSE LIMITS 2 screen appears. The default Low Glucose Limit flashes (**OFF**).
- 8 Select your Low Glucose Limit. The value must be between 2.2 and 21.6 mmol/L (40 and 390 mg/dL).
- 9 Press **ACT**. The default High Glucose Limit flashes (**OFF**).
- 10 Select your High Glucose Limit, then press **ACT**. The value must be between 2.8 and 22.2 mmol/L (50 and 400 mg/dL).
- 11 Press **ACT**. The SET START TIME 3 screen appears.

The second pair of Glucose Limits (Glucose Limits 2) is now set. If you do not need any more Glucose Limits pairs, press **ESC**. If you want to set more Glucose Limits, then repeat steps 6 through 11 above for each pair, up to a maximum of eight.

You are ready to set up your High Repeat.

Setting the High Repeat

After you receive and clear a High Glucose, Rise Rate of Change, or High Predictive alert, the alert will repeat until the condition that has caused the alert is resolved. The High Repeat feature allows you to set how frequently you want the alert to repeat after you clear it the first time. For example, if your blood glucose levels usually take one to two hours to fall after a correction bolus, you may want to set the High Repeat to one or two hours to avoid unnecessary alerts.

Example

You set your High Repeat to 15 minutes. You receive an alert at 13:00 and immediately clear it. However, the condition that caused the alert still exists, so you will receive a second alert at 13:15, 15 minutes after you cleared the first alert. If you immediately clear the alert again, then the alert will repeat at 13:30 and continue to repeat every 15 minutes until you resolve the condition that caused the alert. (This example also applies to the Low Repeat.)

To set up your High Repeat:

- 1 Make sure the EDIT SETTINGS screen is open.
Main > Sensor > Edit Settings
- 2 Select **Hi Repeat**, then press **ACT**. The SET HI ALERT REPEAT screen flashes the default High Repeat time of 1 hour (1:00).
- 3 Select your High Repeat time. The time must be between 5 minutes (0:05) and 3 hours (3:00).
- 4 Press **ACT**. The EDIT SETTINGS screen shows the High Repeat value set.

You are now ready to set up your Low Repeat.

Setting the Low Repeat

After you receive and clear a Low Glucose, Fall Rate of Change, or Low Predictive alert, the alert will repeat until the condition that caused the alert is resolved. After you clear a Low Suspend alarm and resume basal delivery, the alarm will repeat if a low glucose condition still exists. The Low Repeat feature allows you to set how frequently you want the above alerts or alarm to repeat after you clear it the first time. The Low Repeat feature works similar to the High Repeat one. See *Setting the High Repeat* section in this chapter for an example.

- 1 Make sure the EDIT SETTINGS screen is open.
Main > Sensor > Edit Settings
- 2 Select **Low Repeat**, then press **ACT**. The SET LO ALERT REPEAT screen flashes the default Lo Alert Repeat time of 20 minutes (0:20).
- 3 Select your Low Repeat time. The time must be between 5 minutes (0:05) and 1 hour (1:00).
- 4 Press **ACT**. The EDIT SETTINGS screen shows the Low Repeat value set. You are now ready to set up the Predictive Alert.

Setting the Predictive glucose alert

The Predictive alerts calculate when you may reach your Low or High Glucose Limits, then send you an alert before you reach those limits. A Predictive alert tells you that **if your sensor glucose measurements keep falling or rising at the current rate**, you may reach your Glucose Limit in the number of minutes you selected. For a LOW PREDICTED alert, the pump plays three consecutive tones, in falling pitch, if an audible beep has been selected as the alert type. For a HIGH PREDICTED alert, the pump plays three consecutive tones, in rising pitch, if an audible beep has been selected as the alert type.

Understanding the Time Sensitivity Settings

Time Sensitivity means the length of time, in minutes, you want to be notified before you reach your High or Low Glucose Limit. If you set the Low time sensitivity setting at 25 minutes and the High time sensitivity setting at 20 minutes, then the Predictive alerts will be sent:

- **25 minutes before** the sensor glucose reaches your Low Glucose Limit
- **20 minutes before** the sensor glucose reaches your High Glucose Limit



To set up the Predictive alerts:

- 1 Make sure the EDIT SETTINGS screen is open.
Main > Sensor > Edit Settings
- 2 Select **Predict Alerts**, then press **ACT**. The PREDICTIVE ALERTS LOW/HIGH screen shows that the alerts are turned **Off**.
- 3 Select **On**, then press **ACT**. The SET TIME SENSITIVITY LOW/HIGH screen appears and the LOW time setting flashes the default time: 15 minutes (0:15).
- 4 Select the Time Sensitivity (how many minutes before you reach your Low Glucose Limit). To turn off the Low Predictive Alert, reduce the Time Sensitivity to **OFF**.
 - Range: **5 to 30 minutes (0:05 to 0:30)**
 - Increments: **5 minutes**
- 5 Press **ACT**. The HIGH time setting flashes the default time: 15 minutes (0:15).
- 6 Select the Time Sensitivity (how many minutes before you reach your High Glucose Limit). To turn off the High Predictive Alert, reduce the Time Sensitivity to **OFF**.
 - Range: **5 to 30 minutes (0:05 to 0:30)**
 - Increments: **5 minutes**
- 7 Press **ACT**. The EDIT SETTINGS screen appears.
You are now ready to set up the rate alerts.

Setting the Rate of Change Alerts

The Rate of Change alerts tell you when your sensor glucose (SG) changes at, or faster than, the per-minute rate pre-selected by you. There are two alerts:

- **FALL RATE** for SG **decreases** at or faster than your pre-selected rate. The pump plays two consecutive tones, in falling pitch, if an audible beep has been selected as the alert type.
- **RISE RATE** for SG **increases** at or faster than your pre-selected rate. The pump plays two consecutive tones, in rising pitch, if an audible beep has been selected as the alert type.

You can set one or both Rate of Change alerts. The rate of change alert default setting is **OFF**. When you press  or , 0.220 mmol/L (4.0 mg/dL) displays. This value changes up or down with each additional press.

The explanation in the following table shows how the Rate alerts work.

If you set the number	Then
Lower than 0.220 mmol/L/min (4.0 mg/dL/min)	<ul style="list-style-type: none">• The system is more sensitive to SG changes• Alerts are more frequent than if you use a higher number
Higher than 0.220 mmol/L/min (4.0 mg/dL/min)	<ul style="list-style-type: none">• The system is less sensitive to SG changes• Alerts are less frequent than if you use a lower number

To set up the Rate of Change alerts:

- 1 Make sure the EDIT SETTINGS screen is open.
Main > Sensor > Edit Settings
- 2 Select **Rate Alerts**, then press **ACT**. The SET FALL RATE LIMIT screen flashes the rate alert default setting: **OFF**.
- 3 Select the Fall Rate Limit.
 - Range: mmol/L/min from **0.065** to **0.275**
 mg/dL/min from **1.1** to **5.0**
 - Increments: mmol/L/min **0.005**
 mg/dL/min **0.1**
- 4 Press **ACT**. The SET RISE RATE LIMIT screen flashes the rate alert default setting: **OFF**.
- 5 Select the Rise Rate Limit. The range and increments are the same as the Fall Rate Limit (see step 3 above).
- 6 Press **ACT**. The Rate of Change alerts will now use the settings you specified. You are now ready to set up the Low Glucose Suspend.

Setting the Low Glucose Suspend

If the Glucose, Predictive, and Fall Rate of Change alerts are set up, your pump will generate an alert to warn you about low or dropping sensor glucose values. However, if you do not treat or are unaware of your low blood glucose, the Low Glucose Suspend feature is set up to stop insulin delivery. When your sensor glucose reaches or drops below your Low Glucose Suspend limit, your pump will automatically stop insulin delivery, and you will receive a **LOW SUSPEND** alarm. See the *LOW SUSPEND* section in the *Troubleshooting, alerts and alarms* chapter for more information. The pump is set at the factory with the Low Suspend feature turned **Off**. Check with your healthcare professional to determine the Low Suspend setting that is best for you.

To set up the Low Glucose Suspend:

- 1 Make sure the EDIT SETTINGS screen is open.
Main > Sensor > Edit Settings
- 2 Select **Low Suspend**, then press **ACT**. The **LOW SUSPEND** screen appears.
- 3 Select **On**, then press **ACT**.
- 4 The **SET LOW SUSPEND LIMIT** screen flashes the default Low Suspend Limit of 3.3 mmol/L (60 mg/dL). Select your limit. You should set the Low Suspend Limit to a value that is less than your Low Glucose Limit. For example, if your Low Glucose Limit is set at 4.2 mmol/L (76 mg/dL), you should set your Low Suspend Limit to a value that is less than 4.2 mmol/L (76 mg/dL). The Low Suspend Limit should be between 3.3 and 6.1 mmol/L (60 to 110 mg/dL).
- 5 Press **ACT**. The EDIT SETTINGS screen shows the Low Suspend Limit value set. You are now ready to set up the Cal Repeat.

Setting the Cal Repeat

After you receive and clear a **METER BG NOW** alert, the device will repeat the alert until you enter a new meter blood glucose measurement. The Cal Repeat feature allows you to set how frequently you want the alert to repeat after you clear it.

- 1 Make sure the EDIT SETTINGS screen is open.
Main > Sensor > Edit Settings
- 2 Select **Cal Repeat**, then press **ACT**. The **SET CAL REPEAT** screen flashes the default Cal Repeat time of 30 minutes (0:30).
- 3 Set your Cal Repeat time between 5 minutes (0:05) and 1 hour (1:00).
- 4 Press **ACT**. The EDIT SETTINGS screen shows the Cal Repeat time you selected. You are now ready to set your Cal Reminder.

Setting the Cal Reminder

The Cal Reminder feature allows you to set a reminder to calibrate your system. For example, if you set your reminder to four hours, then you will receive a METER BG BY alert (Cal Reminder) four hours before the next meter blood glucose (BG) entry is due (eight hours after your last successful sensor calibration).

To set up the Cal Reminder:

- 1 Make sure the EDIT SETTINGS screen is open.
Main > Sensor > Edit Settings
- 2 Select **Cal Reminder**, then press **ACT**. The CAL REMINDER screen shows **On**.
- 3 Press **ACT**. The SET CAL REMINDER screen flashes the Cal Reminder default time of 1 hour (1:00).
- 4 Select your Cal Reminder time between 5 minutes (0:05) and 6 hours (6:00).
- 5 Press **ACT**. The EDIT SETTINGS screen shows the Cal Reminder time you selected. You are now ready to set up the Auto Calibrate feature.

Turning on the Auto Calibration

Auto Calibration feature determines how to calibrate the pump when using a linked meter. When this feature is turned on, the pump uses every linked meter BG reading in the range of 2.2-22.2 mmol/L (40-400 mg/dL) for calibration. When this feature is turned off, the pump asks you every time if you want to use the linked meter BG for calibration.

To turn on the Auto Calibration:

- 1 Make sure the EDIT SETTINGS screen is open.
Main > Sensor > Edit Settings
- 2 Select **Auto Calibrate**, then press **ACT**. The AUTO CALIBRATION screen is displayed.
- 3 Select **On**, then press **ACT**. The EDIT SETTINGS screen appears showing that the Auto Calibration is now turned On. You are now ready to enter the ID number of your transmitter.

Entering the transmitter ID


The transmitter ID (serial number) starts with 2 and is found on the flat side of your transmitter. You must enter the transmitter ID so that the transmitter and the pump can communicate with each other.

Write the Transmitter ID of the transmitter you are using here: _____.

- 1 Make sure the EDIT SETTINGS screen is open.

Main > Sensor > Edit Settings

- 2 Select **Transmtr ID**, then press **ACT**. The SET TRANSMITTER ID screen appears. The screen shows seven dashes for the Transmitter ID, the first of which (on the left) is flashing.
- 3 Select the first number of the Transmitter ID (the number on the screen must match the first number of the Transmitter ID). The transmitter ID is located on the transmitter label, as shown.

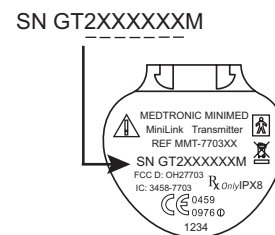

- 4 Press **ACT**. The second dash flashes.
 - 5 Select the next number of the Transmitter ID, then press **ACT**.
 - 6 Repeat step 5 until all seven digits of the Transmitter ID have been selected.
 - 7 A message shows that the Transmitter ID has been changed. After about 30 seconds, the EDIT SETTINGS screen appears. It shows the Transmitter ID you entered.
- You are now ready to set up the Weak Signal.

Setting up the Weak Signal feature

The Weak Signal feature allows you to determine how quickly you will receive a Weak Signal alert when there is no communication between the transmitter and the pump. For example, if you set the Weak Signal value to 30 minutes, you will receive a Weak Signal alert 30 minutes after loss of communication.

To set up the Weak Signal feature:

- 1 Make sure the EDIT SETTINGS screen is open.
Main > Sensor > Edit Settings
- 2 Select **Weak Signal**, then press **ACT**. The SET WEAK SIGNAL screen flashes the default time of 30 minutes (0:30).
- 3 Select your Weak Signal time between 5 minutes (0:05) and 40 minutes (0:40).
- 4 Press **ACT**. The EDIT SETTINGS screen shows the length of time you entered for the Weak Signal option.
- 5 You are now ready to set up the sensor glucose graph timeout.

Setting the sensor glucose graph timeout

The sensor glucose graph timeout is how long the graphs show on the pump screen until the HOME screen automatically appears. The default timeout is two (2) minutes, and you can set it for 2, 4, or 6 minutes.

Using a continuous graph display


You can also select **NONE** to set the sensor graph timeout to display continuously. If you select **NONE**, the sensor graphs will show until another screen or menu is selected, or an alarm/alert occurs. Selecting continuous graph display uses more battery power and decreases battery life. To return to the HOME screen from the sensor glucose graphs, press **ESC** until the HOME screen appears.

To set the sensor glucose graph timeout:

- 1 Make sure the EDIT SETTINGS screen is open.
Main > Sensor > Edit Settings
- 2 Select **Graph Timeout**, then press **ACT**. The SET GRAPH TIMEOUT screen flashes the default time of 2 minutes (0:02).
- 3 Select your Graph Timeout. The time you select must be: 2 minutes; 4 minutes; 6 minutes; or NONE.
- 4 Press **ACT**. The EDIT SETTINGS screen shows the Graph Timeout you selected.
You are now ready to review your sensor glucose monitoring settings.

Reviewing your settings

Make sure that all the settings you have made are correct before you use the system. Use the following procedure to review your settings:

- 1 Open the REVIEW SETTINGS screen.
Main > Sensor > Review Settings
- 2 Press  to scroll through all settings to make sure that they are correct.
- 3 To change any settings, return to the EDIT SETTINGS menu.
Main > Sensor > Edit Settings
- 4 Save your settings after you have set your preferences. This will allow you to restore the settings you have saved if you receive an alarm or error that resets your settings.

Using the Sensor Demo

The Sensor Demo shows you a demonstration of a sensor graph. To view the Sensor Demo screen, you must first turn on the Sensor Demo feature. After the Sensor Demo feature is turned on, the first line of the STATUS and SENSOR STATUS screens will show **SENSOR DEMO: On**.

Turning on the Sensor Demo feature

- 1 Make sure the EDIT SETTINGS screen is open.



Main > Sensor > Edit Settings

The EDIT SETTINGS screen shows **Sensor Demo: Off**.

- 2 Press ACT. The SENSOR DEMO ON/OFF screen appears with **Off** selected.
- 3 Select **On**, then press ACT. The EDIT SETTINGS screen shows that SENSOR DEMO is now turned **On**.

CAUTION: The Sensor Demo feature affects the sensor demonstration graph only. All other features, screens and options are in the normal mode.

Viewing the Sensor Demo screens

- 1 From the HOME screen, press ESC. The screen briefly shows SENSOR DEMO, then the first Sensor Demo graph appears.
- 2 You can move the cursor (flashing vertical line) through the graph to see examples of how the continuous graphs appear. Press  to move the cursor left, press  to move the cursor right.
- 3 You must turn the Sensor Demo off in order to show your actual, continuous sensor glucose graphs. To turn the Sensor Demo off, follow the above procedure and select **Off**.

The transmitter

The Medtronic MiniLink Transmitter (MMT-7703) is a device that takes electronic signals generated by the glucose sensor and sends them by radio frequency to the pump. It has a tester (MMT-7706) and a charger (MMT-7705 or MMT-7715).


Transmitter charger

The transmitter contains a non-replaceable, rechargeable battery that you can recharge as needed with the charger. The charger has a green light that shows the charging status and a red light that communicates any problems during charging. If you see a red light, see *Understanding your transmitter, tester, and charger*, on page 191. The charger needs a AAA or LR-03 battery to operate.

NOTE: A new AAA or LR-03 battery contains enough power to recharge the transmitter more than 40 times. If the battery is installed incorrectly or is low, the charger will not work. Repeat the battery installation steps using a new battery.


Installing a battery in the gray charger

- 1 Push the battery cover in and slide it off (as shown in the following illustration).
- 2 Insert a new Alkaline AAA or LR-03 battery. Make sure the + and - symbols on the battery align with these same symbols shown on the charger.
- 3 Slide the cover back on the charger until it clicks into place.

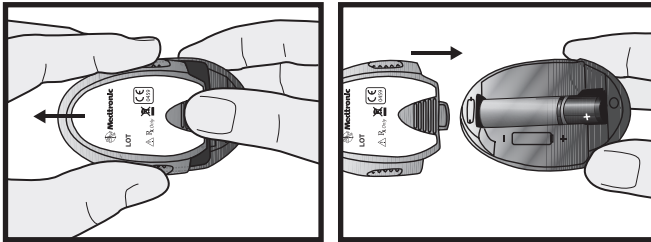

Installing a battery in the blue charger

- 1 Remove the charger battery cap by turning it counter-clockwise 1/4 turn using a coin in the groove of the cap.

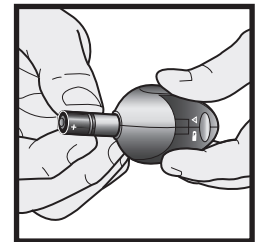



- 2 Insert a new AAA or LR-03 battery with the flat (-) end first. Make sure that you align the small bumps on the battery cap with the small notches in the charger's battery opening. Push in the cap all the way using a coin. Turn the cap clockwise 1/4 turn to close.

Charging the transmitter

Before using the transmitter for the first time, you must fully charge the transmitter battery. This may take up to 8 hours. It is recommended to recharge the transmitter after each sensor use. A fully charged transmitter battery works more than 14 days without recharging. After 14 days of use, the transmitter recharges in less than 2 hours.

To charge the transmitter:

- 1 If a green light on the transmitter is lit or flashing, do not connect it to the charger. The transmitter will not charge with its green light on. Wait until the green light turns off (approximately 30 seconds), and then connect the transmitter to the charger.
- 2 Connect the transmitter to the charger by lining it up, flat side down, with the charger. Push the two components together fully.
- 3 Within 10 seconds after the transmitter is connected, a green light on the charger will flash for 1-2 seconds as the charger powers on. For the rest of the charging time, the charger's green light will continue to flash in a pattern of four flashes with a pause between the four flashes.
- 4 When charging is complete, the green light on the charger will stay on, without flashing, for 15-20 seconds and then turn off.
- 5 After the green charger light turns off, disconnect the transmitter from the charger. The green light on the transmitter will flash for about 5 seconds and then turn off.

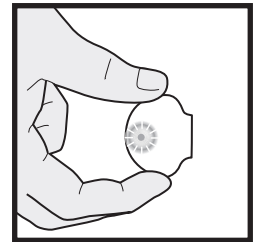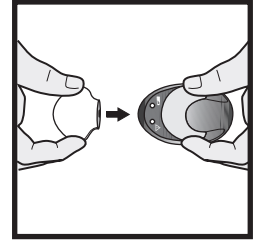

Starting the sensor

To start the sensor working, you must complete the following steps in order:


- Insert a battery into the transmitter's charger.
- Charge the transmitter battery.
- Set up the sensor features.
- Insert the sensor.
- If you are using the Enlite Sensor (MMT-7008), wait five minutes. If you are using a different sensor (MMT-7002/MMT-7003), wait 10-15 minutes.
- Connect the transmitter to the sensor.
- Start the sensor and wait two hours.
- Enter your first meter BG.

Inserting the Enlite™ sensor


This procedure shows you how to insert the Enlite Sensor (MMT-7008) using the Enlite Serter (MMT-7510). If you are using a different sensor (MMT-7002/MMT-7003), use the insertion instructions in the sensor user guide.

CAUTION: Healthcare professionals and caregivers should use universal precautions when handling the sensor.

- 1 Wash your hands.
- 2 Select an insertion site.
- 3 Clean the insertion site with alcohol. Allow the area to air-dry.
- 4 Open the sensor package.
- 5 Remove the pedestal from the package. Place the pedestal on a flat surface.
- 6 Push the Serter down onto the pedestal until the base of the Serter sits flat on the table.

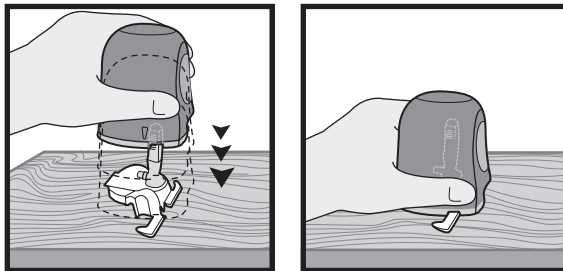


- 7 To detach the Serter from the pedestal, place two fingers on the pedestal arms and slowly pull the Serter straight up. Do not detach the pedestal from the Serter in mid-air, as this might damage the sensor.

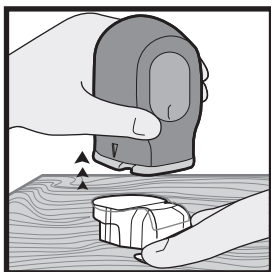

WARNING: Never point a loaded Serter toward any body part where insertion is not desired.


- 8 Place the base of the Serter flat against your insertion site.

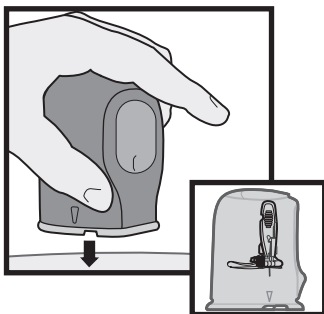

WARNING: The Serter injects sensor upon button release.

NOTE: A second button push is required to remove the Serter from the sensor.

9 To insert the sensor, complete the following four steps:

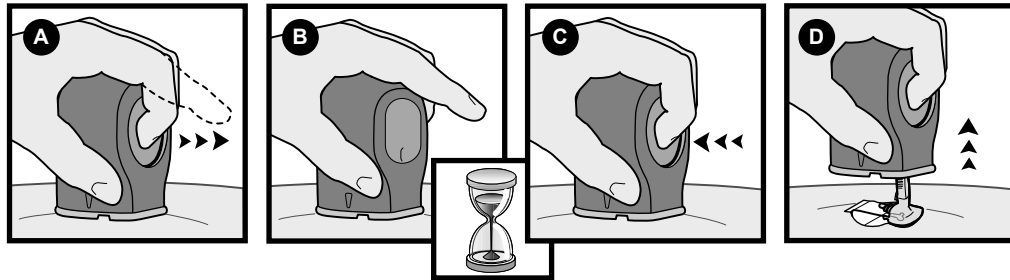

A Press the button in and then release.


B Wait 5 seconds to allow adhesive time to stick to skin.

C Press and hold in the button.

D While holding the Serter button in, lift the Serter away from the skin.


10 If you inserted the sensor into yourself, complete step A. If you are a healthcare professional or caregiver who inserted the sensor into a Patient, complete step B.

Patient:

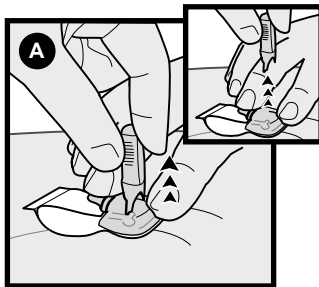

A Gently hold the base of the sensor against skin. Hold the needle housing at the top and slowly pull straight up, away from the sensor.


Healthcare professional or caregiver:

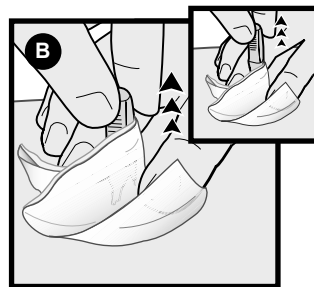

B Wrap sterile gauze around the sensor (as shown in image B). Gently hold the base of the sensor against skin. Hold the needle housing at the top and slowly pull straight up, away from the sensor.

WARNING: If bleeding occurs, apply steady pressure using sterile gauze or clean cloth for up to 3 minutes.



- 11 Straighten the sensor adhesive tab so that it lies flat against the skin.

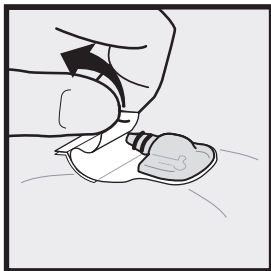

- 12 While holding the sensor in place, gently lift the adhesive tab. Carefully remove the bottom piece of white paper from under the adhesive pad. Press the adhesive against your skin.



- 13 Before connecting the device, apply overtape. **Do not cover the sensor connectors or adhesive tab with overtape.**

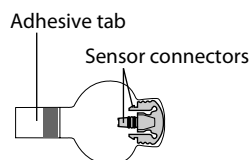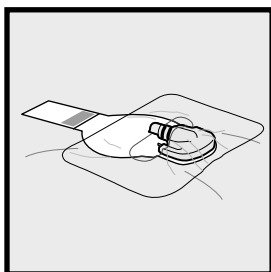


- 14 Follow the directions in the next section for connecting the transmitter to the sensor.

Connecting the transmitter to the sensor

Before connecting the transmitter to the sensor, you must have the transmitter fully charged and set up (see *Charging the transmitter* and *Entering the transmitter ID* sections in this chapter). You must also have the features set up and the sensor inserted (see the instructions in this chapter and those provided with your sensor).


To connect the transmitter to the sensor:

- 1 Connect the transmitter to the sensor, as shown.

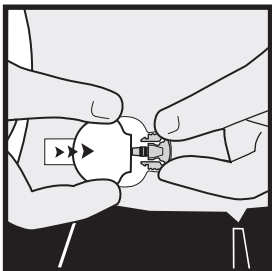


NOTE: Overtape is not shown in the following steps.

- 2 The green light on the transmitter will flash. If the light does not flash, disconnect it from the sensor, wait for several seconds, and reconnect. If the transmitter light still does not flash, charge the transmitter.

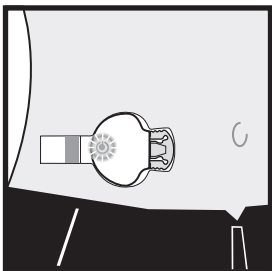

- 3 Cover the transmitter with the Enlite sensor's adhesive tab.

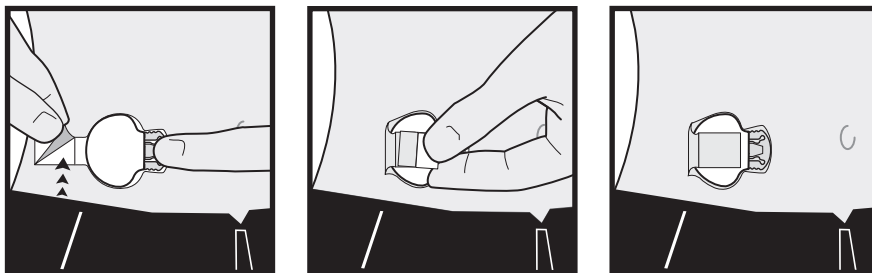

- 4 [Optional]: Apply occlusive dressing over the transmitter and the sensor.
- 5 Follow the instructions in the next section to start your sensor.

Preparing the sensor for communication

You are now ready to use your pump to start communicating with the sensor and start sensor initialization. The pump will start a timer and notify you when the sensor is ready to use.

Use your pump to follow the steps below:

- 1 Go to the LINK SENSOR MENU screen.
Main > Sensor > Link to Sensor
- 2 Select **New Sensor**. Press **ACT**.
- 3 Connect the transmitter to your sensor now if you have not already done so. If the screen times out, do **not** disconnect the sensor. Start from step 1 again.
- 4 After your sensor is connected, press **ACT**. The **SENSOR READY 2 HRS** screen appears.
- 5 The sensor will enter a two hour initialization period. Press any button to continue. After 2 hours your pump will alert you to enter a meter BG to calibrate your sensor.

Calibrating the sensor

Two hours after you use your pump to start the sensor, your pump will alert you to enter a meter BG (**METER BG NOW**). This meter BG entry will be the **first** calibration for your sensor. You have to wait 10-15 minutes after calibration to see the first sensor glucose reading on the pump screen.

After the first calibration, you must calibrate your sensor every 12 hours. If you fail to enter a meter BG reading after 12 hours, your pump displays the **METER BG NOW** alert. Your pump will then stop calculating glucose values. However, about 20 minutes after you have entered a meter BG, your pump will continue calculating glucose values.

Follow these guidelines for best calibration results:

- Calibrate three to four times spread out throughout the day.
- Avoid calibrating when your glucose is changing rapidly, such as after eating or exercise.
- Enter your meter BG reading into the pump immediately after testing your BG. Do not wait to enter it later.
- Always use clean dry fingers when you check your blood glucose.
- Only use your fingertips to obtain blood samples for calibration.

After the transmitter successfully transmits signals to the pump, you may choose to put an occlusive dressing over the transmitter and the sensor.

Entering meter BG

Sensor calibration will be successful only if the BG entry is in the range of 2.2-22 mmol/L (40-400 mg/dL). You should calibrate three to four times spread throughout the day for optimal results.


- 1 Use any of the following methods for your sensor calibrations, including the first calibration after the initialization:
 - a. To enter a BG value automatically through the linked meter, test your BG with a fingerstick. The value will be automatically sent from the meter to your pump. Make sure that your pump and the meter are properly programmed for communication. See the *Meter option* section in the *Using the Bolus Wizard feature* chapter for more information.
 - b. To enter a meter BG value manually, follow this path:

Main > Sensor > Calibrate

Enter the BG value, and then press ACT. Your sensor is now calibrated with the BG value you entered.

- If the Bolus Wizard feature is off, the SENSOR MENU screen appears.
- If the Bolus Wizard is on and the saved BG is below your BG Target range, the Cal Saved screen appears. This screen indicates that your saved BG is below your BG Target set in the Bolus Wizard feature, and the pump beeps three times. See *Setting the BG Targets, on page 83* for more information.
- If the Bolus Wizard is on, the saved BG is above your BG Target range, and the calculated correction estimate is more than the preset scroll rate increment, the Cal Saved screen appears. This screen indicates that your saved BG is above your BG Target set in the Bolus Wizard feature,

and the pump beeps three times. Go to the ENTER BG screen of the Bolus Wizard feature to enter a correction bolus. See the *Using the Bolus Wizard feature* chapter for more information.

- If the information is not correct, select **Cancel**, then press **ACT**. The SENSOR MENU screen reappears with **Calibrate** selected. Repeat the procedure and enter the correct information.
- c. To enter a meter BG through the Capture Event Menu, follow this path:
Main > Capture Event > Enter BG
See the *Entering BG measurements* section in the *Basic programming* chapter for more information.
- d. If at the time of sensor calibration you also need to perform a bolus using the Bolus Wizard, you can use the meter BG for both calibrating the system and bolusing. To enter a meter BG value manually in the Bolus Wizard, press the  button on your pump, or use the MAIN MENU to go to the ENTER BG screen:

Main > Bolus > Use Bolus Wizard

- 2 If you are manually entering your meter BG, use the up and down arrow buttons on your pump to do so. Press **ACT**. For the Bolus Wizard option above, select **Yes** when the following screen displays: BG to update Sensor?

Using your sensor

Status screens

Your STATUS screens tell you what is going on in your pump. In the SENSOR STATUS screen you can check the status of sensor information including when your next calibration will be needed, your sensor's age, and the state of your transmitter battery.

To get to your status screens:

- 1 From the HOME screen, press the **ESC** button twice. This takes you to your pump status screen.
- 2 To see the Sensor Status screen, press **ESC** one more time. This screen will only be available if the Sensor feature has been turned On.



Reading the sensor glucose graphs

NOTE: If your graphs say **DEMO** in the upper left corner of your screen, you must turn off the **SENSOR DEMO** mode to see your actual graphs (see the *Using the Sensor Demo* section in the *Sensor features* chapter).

Your pump shows an updated, continuous glucose measurement. This measurement is generated by data sent from the sensor to the transmitter, and then to the pump every five minutes. The pump converts these measurements to glucose graphs that include the following information:

- The most recent sensor glucose measurement or the reason that one is not displayed.
- The historical sensor glucose measurements or the reason that one is not displayed for the last 3, 6, 12 or 24 hours.
- Arrows that show the relative rate at which the most recent sensor glucose levels have risen or fallen.

The following conditions can prevent a continuous sensor glucose measurement from being taken:





- A Lost Sensor alert.
- A Sensor Error alert.
- A Change Sensor alert.
- A Cal Error alert.
- Find Lost Sensor (15 minutes).
- Weak signal.
- Reconnect old sensor (two hours).
- A meter BG entry is needed for calibration.

Two other conditions also prevent a continuous SG measurement and generate a WARM UP message on your sensor glucose graphs:

- You started a new sensor, and it has entered a two-hour initialization. The sensor is not ready to take glucose measurements.
- You entered a BG for calibration, and your pump is waiting to receive the first valid SG value.

Opening and viewing the graphs


The following table describes how to use the graphs.

If you want to:	Do the following:
Open the graphs	From the HOME screen, press ESC . The last graph you viewed (3-, 6-, 12- or 24-hour) shows the details for your most recent sensor glucose (SG) measurement.
View information in a graph	<p>If there is no SG measurement at the time you select, the reason for the missing measurement will be shown on the right side of the graph (for example, LOST SENSOR).</p> <ul style="list-style-type: none">• When you open the graphs, the cursor (flashing vertical line) is always on the right edge of the graph and the graph shows the most recent SG measurement.• Press  to move the cursor left to see an earlier SG measurement.• Press  to move the cursor right to return to more recent SG measurements.
View different graphs	<ul style="list-style-type: none">• With the cursor at the far right edge of a graph, you can press  to open the next graph.• If the cursor is anywhere to the left of the far right edge, press ESC to move it to the right, then press  to open the next graph.

The graphs


All graphs show High and Low Glucose Limits lines, a continuous SG measurement line, a data section, and the cursor (flashing vertical line).

When you open any graph, the cursor flashes on the right edge of the graph. There are three marks on the cursor at 5, 10, and 15 mmol/L (100, 200, and 300 if mg/dL is selected). The data section shows the most recent SG measurement or the reason why no measurement shows, and the time in the data section matches the time at the top of the screen.

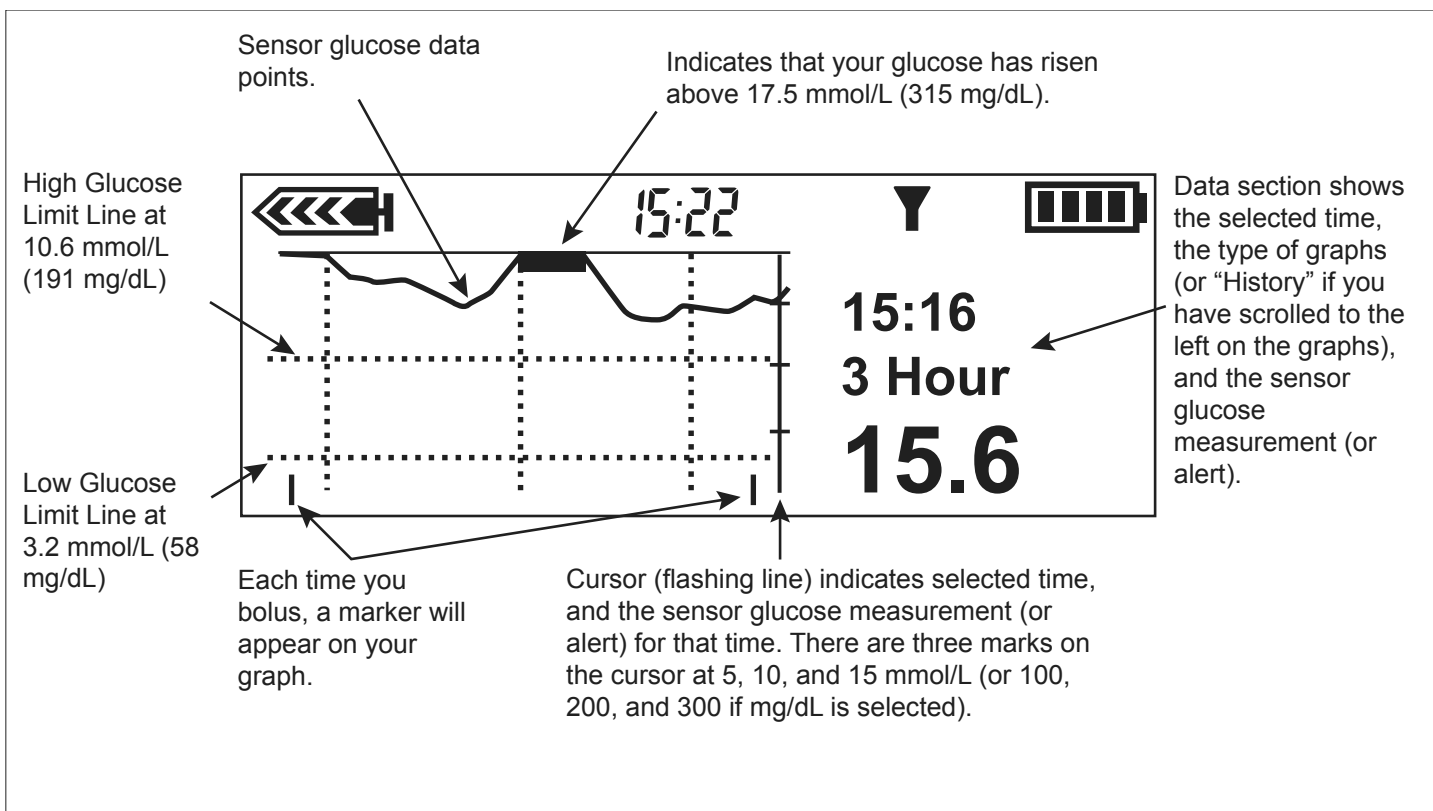

When you move the cursor left to select an earlier SG measurement, the data section of the graph turns black, the word History appears here, and the time in the data section changes to show the time when the SG measurement was taken or the reason why no SG measurement shows for that time.

NOTE: There are many reasons your graphs may not show an SG measurement, but your system does not alert for every one. Therefore, your system's alarm/alert history may not match the number of times your graphs do not show an SG measurement.

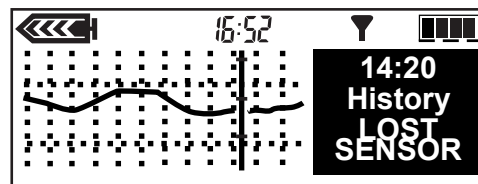


Examples of continuous sensor glucose graphs

The following pages show examples of the different graph types.


3-hour graph

Each part of the graph between the vertical dashed lines represents one hour. In this example, the selected SG measurement was taken at 15:16, and the value at that time was 15.6 mmol/L.

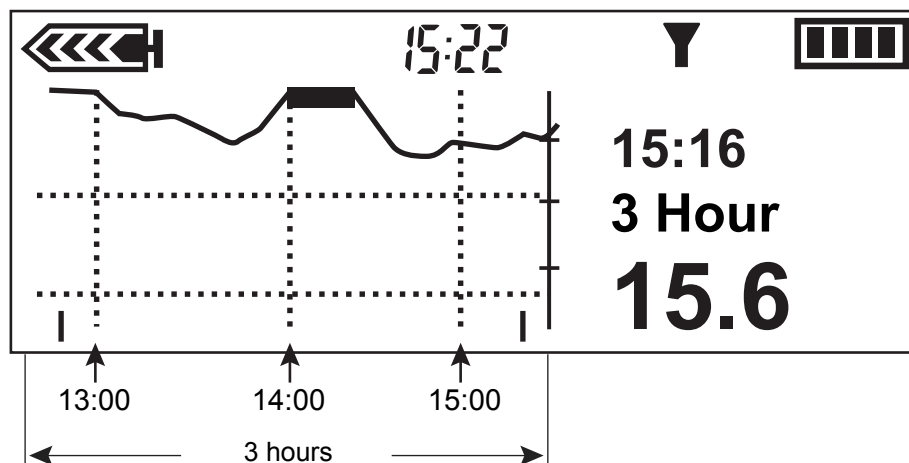


6-hour graph

Each part of the graph between the vertical dashed lines represents one hour. In this example, the selected SG measurement was taken at 1:12, and the measurement at that time was 8.3 mmol/L.

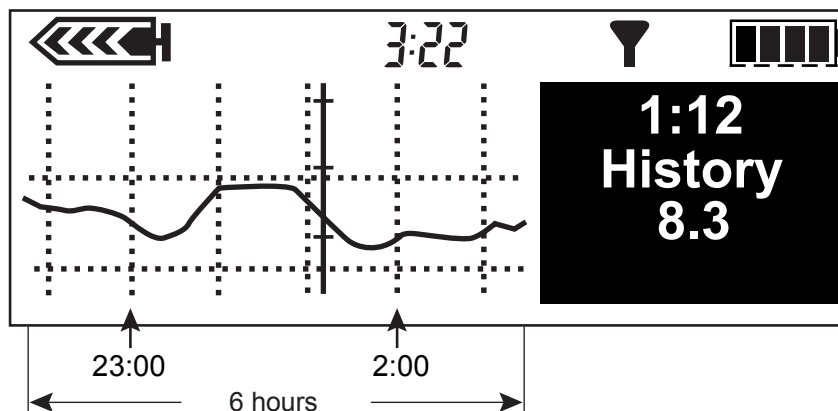


12-hour graph

Each part of the graph between the vertical dashed lines represents one hour. In this example, the reason no measurement is displayed is LOST SENSOR and it occurred at 14:20.

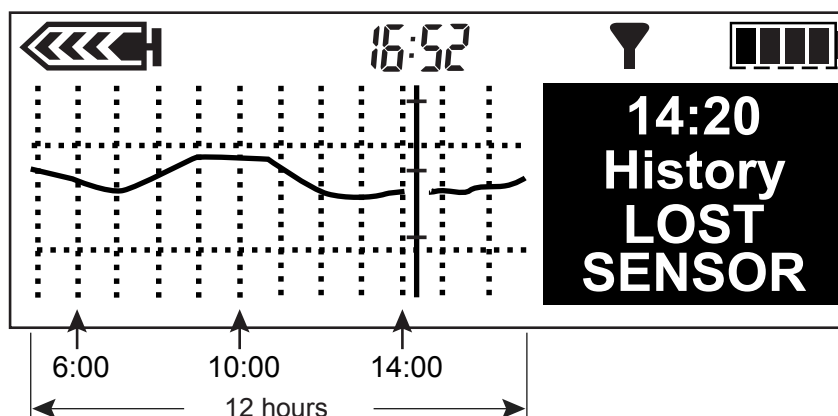

24-hour graph

Each part of the graph between the vertical dashed lines represents 12 hours. The shaded area helps you more easily see the sensor glucose measurements for the previous night. In this example, the selected SG measurement was taken at 11:50 and the value at that time was 11.7 mmol/L.

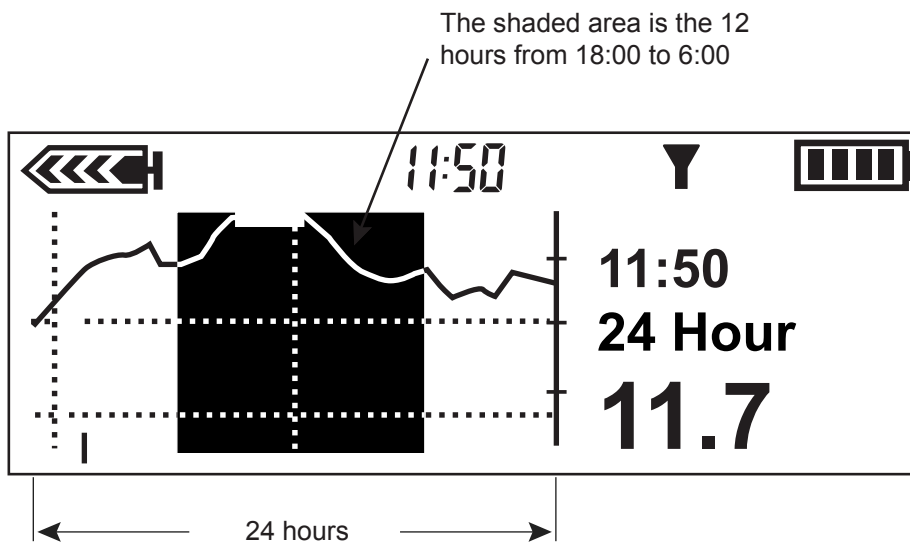


How to check for rapid changes in sensor glucose

If your sensor glucose (SG) rises or falls faster than a certain rate, rapid change arrows automatically appear in the graphs next to the SG measurement, as shown in the following section.


The rapid change arrows

These examples show how the rapid change arrows tell you when your SG rises or falls faster than a certain per-minute rate. The arrow direction shows if your SG measurements are rising or falling. The number of arrows (one or two) shows how rapidly your SG is changing.


- One up arrow \uparrow shows that your SG has been rising at a rate of 0.05 to 0.11 mmol/L (1 to 2 mg/dL) per minute.

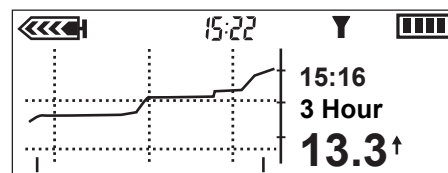


- One down arrow \downarrow shows that your SG has been falling at a rate of 0.05 to 0.11 mmol/L (1 to 2 mg/dL) per minute.

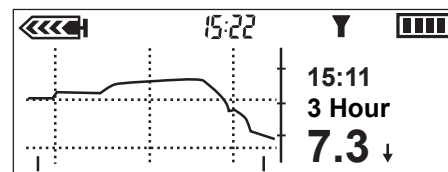

- Two up arrows $\uparrow\uparrow$ show that your SG has been rising at a rate of 0.11 mmol/L (2.0 mg/dL) or more per minute.

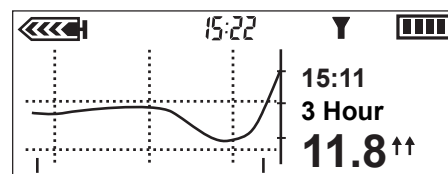

- Two down arrows $\downarrow\downarrow$ show that your SG has been falling at a rate of 0.11 mmol/L (2.0 mg/dL) or more per minute.

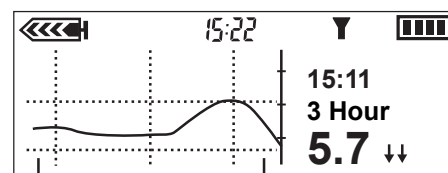

How to make glucose alerts silent

The Alert Silence feature allows you to silence glucose alerts for a preset amount of time. You may use this feature when you are in a business meeting, watching a movie in a theater, and so on. The system records the time and glucose measurement for each alert during this silent period. You can see this alert information in the SENSOR ALERT HISTORY screen. See the *Sensor alert history* section in this chapter for details.

When one or more glucose alerts occur during the time when the Alert Silence feature is turned on, the pump displays the ALERT SILENCE screen without beeping or vibrating. If this alert is not cleared before the end of the silenced period, the pump beeps or vibrates periodically at the end of the preset time until the alert is cleared. Press **ESC**, then **ACT** to clear the alert.

You can make the glucose alerts silent in the following ways:

- **Off** — This means the Alert Silence feature is **Off** and the alerts are **On**. The pump beeps or vibrates when any sensor glucose alert occurs.
- **Hi Alerts** — The pump does not beep or vibrate if a high or rising sensor related alert occurs during the specified time.
- **Lo Alerts** — The pump does not beep or vibrate if a low or falling sensor related alert occurs during the specified time.
- **Hi and Lo Alerts** — The pump does not beep or vibrate if a high or low, rising or falling sensor related alert occurs during the specified time.
- **All Sensor Alerts** — The pump does not beep or vibrate if any sensor glucose alert occurs during the specified time, including the Hi and Lo glucose alerts, Cal Reminder, and Meter BG Now alerts.

The amount of time that you can make glucose alerts silent is:

- Minimum: 30 minutes
- Maximum: 24 hours

To make alerts silent:

- 1 Open the SENSOR MENU.
Main Menu > Sensor
- 2 The SENSOR MENU shows you the current status of the Alert Silence feature.
- 3 Depending on your current settings, one of the following **Alert Silence** options will appear in the SENSOR MENU.
 - Alert Silence: Off
 - Alert Silence: Hi
 - Alert Silence: Lo

- Alert Silence: Hi/Lo
 - Alert Silence: All
- 4 To change the current alert silence settings, select **Alert Silence**, then press **ACT**.
 - 5 Select the alert silence option you want, then press **ACT**. The SET DURATION screen flashes a set of dashes or the length of time that you entered the last time you set an alert to silent.
 - 6 Select your alert silence length between 30 minutes (minimum) and 24 hours (maximum), then press **ACT**. The pump is in Special mode for the set silenced period.
The SENSOR MENU shows the new alert that you selected to be silent.
 - 7 Press the **ESC** button until the MAIN MENU appears.

Calibration history

The CALIBRATE HISTORY screen lists the successful sensor calibration BG values that were entered into the pump more than 15 minutes ago.

To view your Calibration History:

- 1 Go to the CALIBRATE HISTORY screen.
Main > Sensor > Calibration Hist.
- 2 The CALIBRATE HISTORY screen appears. It shows up to 28 calibration values.

Sensor alert history

The SENSOR ALERT HISTORY screen lists all of the sensor alerts that have occurred, displaying up to 36 alerts.

To view your Sensor Alert History:


- 1 Go to the SENSOR ALERT HISTORY screen.
Main > Sensor > Sensor Alert History
- 2 The SENSOR ALERT HISTORY screen appears. The most recent alert is highlighted. Select the alert you want to get more information on, then press **ACT**.
- 3 A new screen will appear with the alert details.

Disconnecting the transmitter and removing the sensor

Disconnecting the sensor from the transmitter

If you are not going to replace the sensor, turn the Sensor features **Off** to avoid getting a LOST SENSOR alert.

- 1 Carefully remove any occlusive dressing from the transmitter and sensor.

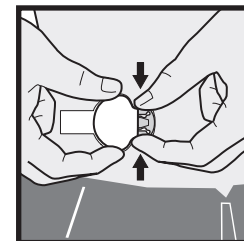

- 2 For the Enlite sensor, remove the adhesive tab from the transmitter.
- 3 Hold the transmitter as shown, and pinch the flexible side arms of the sensor between your thumb and forefinger.
- 4 Gently pull the transmitter away from the sensor.

Removing the sensor

Gently pull the sensor from your body to remove it. Place it in a sharps container.

Storage

If you do not plan to use your transmitter in the next 30 days, connect it to the charger for storage.

Using your system in water

Your pump must not be used in water and needs to be removed if planning water activities.

You should shower, bathe and swim with the transmitter and the sensor by following the guidelines below:

- 1 Disconnect the infusion set from the pump and remove the pump. The pump is not water-tight.
- 2 After the transmitter and sensor are connected, they form a water-tight seal to a depth of 2.4 meters (8 feet) for up to 30 minutes. You can shower and swim without removing them. Avoid taking hot baths as this may significantly reduce the sensor's life.

- 3 Once out of the water, put the pump back on and reconnect the infusion set.
- 4 Check your infusion set tape and the sensor tape to make sure that they are not damaged.

Utilities

Alarm review

You can review alarms and their details in the ALARM HISTORY screen. This screen shows up to 36 past alarms, errors, or LOW RESERVOIR and LOW BATTERY alerts.

- 1 Go to the ALARM HISTORY screen.
Main > Utilities > Alarm > Alarm History
- 2 Scroll through your past alarms.
- 3 In the ALARM HISTORY screen, select the alarm you want to review and press **ACT**. The details for that alarm will appear on the screen.
- 4 Press **ESC** to return to the ALARM HISTORY screen. Select another alarm to review, or exit the menus if you are done.

Setting your alert type

You can select the type of alert your pump uses (for alarms, special conditions and programming). You can select a vibrate (silent) alert, or an audible beep alert. There are three beep types: long, medium and short tones. The factory setting for this feature is beep-medium.

The vibrate alert type is disabled if you use the block feature, and vibrate must be set again once block is turned off. Vibrate uses more battery power than the beep alert type and may shorten battery life. If your alert type is set to vibrate and you get a LOW BATTERY alert, your pump will use the beep alert type instead to conserve battery power.

- 1 Go to the ALERT TYPE screen.
Main > Utilities > Alarm > Alert Type
- 2 Select your alert type and press **ACT**. That alert type is now active. Exit the menus.

Auto-off

The factory setting for this feature is **Off**. This is a safety feature that stops insulin delivery after a defined time period (from 1 to 24 hours). If the pump detects that no buttons have been pressed for the selected amount of time in Auto-off, insulin delivery will stop and an alarm will sound. You may choose to program this feature into your pump based on the number of hours that you usually sleep at night. Discuss what uses and settings are best for you with your healthcare professional.

- 1 Go to the AUTO OFF DURATION screen.

Main > Utilities > Alarm > Auto Off

- 2 Set the number of hours you want to set and press **ACT**.

If you do not want to use the Auto Off feature, make sure the hour is set to zero (0).

- 3 The screen will return to the ALARM MENU. The Auto Off feature is now set. Exit the menus.

LOW RESV WARNING (Low reservoir warning)

Allows you to program the pump to sound an alert before your reservoir is empty. You can select one of these warning types:

- a specified number of units that remain in the reservoir
- a specified maximum amount of time that remains before the reservoir will be empty

The factory setting for this feature is 20 insulin units.

If you use **Time** as the low reservoir warning type and you deliver large boluses, the actual time remaining could be less than the warning time. The warning time is based on the basal insulin delivery rate. The **Time** low reservoir warning type is intended to let you know if you will have enough insulin while you are sleeping.

- 1 Go to the RESV WARNING TYPE screen.

Main > Utilities > Alarm > Low Resv Warning

- 2 Select **Insulin Units** or **Time** and press **ACT**.

WARNING: When the pump detects a low reservoir condition during a bolus or Fill Cannula delivery, the LOW RESERVOIR alert displays after the delivery is finished. Make sure to check the volume of your reservoir to ensure enough insulin is available.

- a. For Insulin units:

In the WARNING UNITS screen enter the number of units you want remaining for the first warning, then press **ACT**.

The pump displays a LOW RESERVOIR alert first when the specified units remain, then again when half that remaining amount is used.

b. For Time:

In the WARNING TIME screen enter the amount of remaining time you want the first warning, then press **ACT**.

The pump displays a LOW RESERVOIR alert first when the specified time remains, then again one hour before empty.

Review daily insulin totals

The DAILY TOTALS screen provides a day-by-day history of the total amount of insulin that delivered for the last 32 days, including today. This screen includes all bolus and basal amounts delivered midnight-to-midnight for each of the last 32 days. The Today line in the DAILY TOTALS screen shows the amount of insulin you delivered so far that day.

The insulin used to fill the tubing or the cannula is not included in the DAILY TOTALS screen. This amount is counted separately and shown in the HISTORY screen of the RESERVOIR + SET MENU.

Question: Why should I review my daily totals?

Answer: Comparing your daily insulin deliveries to your blood glucose records helps you and your healthcare professional identify your optimal daily insulin rate(s).

Question: What is included in the daily totals?

Answer: Daily totals include all basal and bolus insulin deliveries, sensor glucose, meter blood glucose, and AUC information for that day.

Pump data management

The pump data management feature allows you and your healthcare professional to view and manage your basal and bolus insulin delivery, food intake, BG data, SG data, and AUC data with averages. You can view the details of each day individually or you can average the data over a specified number of days (up to 32 days).

- **Daily totals** – The system automatically tracks certain types of information for you on a day-to-day basis. The daily totals are the totals for a single day.
- **Averages** – The system automatically calculates averages of certain types of information for you. There are two basic types of averages:
 - An average over a number of days (you select the number of days)
 - An average for a single day

The insulin delivery, food intake, sensor glucose information, meter blood glucose information, and AUC data includes daily totals and averages, as described in the following sections.

Insulin delivery information

This table lists the different types of insulin delivery and food intake information that the system automatically records. This information is calculated by using all of the basal, bolus, and food data for that day.

An asterisk (*) next to a date means that the data for that day will not be used for the calculations that are displayed in the Daily Average option. Anything, such as a time or date change, that shortens a day at least one hour and twelve minutes (01:12) will cause the asterisk to be displayed.

Item	Description
Carbs	The total number of carbohydrates entered using the Bolus Wizard and Capture Event (Meal marker) features.
Tot Insulin	The total amount of insulin delivered with the pump.
Basal (first line)	The total amount of insulin delivered as basal .
Basal (second line)	The percent of insulin delivered as basal .
Bolus (first line)	The total amount of insulin delivered as bolus .
Bolus (second line)	The percent of insulin delivered as bolus .
Food Only	The total amount of insulin delivered using the Bolus Wizard feature with food amount only .

Item	Description
#Food Only	The number of times the Bolus Wizard feature delivered a food bolus only .
Corr Only	The total amount of insulin delivered using the Bolus Wizard feature with BG correction amount only .
#Corr Only	The number of times the Bolus Wizard feature delivered a BG correction bolus only .
Food+Corr	The total amount of insulin delivered using the Bolus Wizard feature with food and BG correction amount .
#Food+Corr	The number of times the Bolus Wizard feature delivered a food and BG correction bolus .
Manual	The total amount of insulin delivered manually .
#Manual	The number of manual insulin deliveries.

Sensor glucose information

This table lists the different types of sensor glucose information that the system automatically records. The sensor glucose information is calculated by using all of the sensor glucose measurements that were received by the pump during the day.

Item	Description
SG Average	The average sensor glucose measurement.
Above Hi Limit	The percentage of sensor glucose measurements that were above your High Glucose Limit.
Within Limits	The percentage of sensor glucose measurements that were between your Low Glucose Limit and your High Glucose Limit.
Below Low Limit	The percentage of sensor glucose measurements that were below your Low Glucose Limit.
SG Std. Dev.	The standard deviation of the sensor glucose measurements.
#Hi Alert	The number of High Glucose alerts.
#Lo Alert	The number of Low Glucose alerts.
#Predict Hi	The number of High Predictive Glucose alerts.
#Predict Lo	The number of Low Predictive Glucose alerts.


Item	Description
#Rise Rate	The number of Rise Rate alerts.
#Fall Rate	The number of Fall Rate alerts.
Number SGs	The number of SG measurements received from the transmitter.



Meter blood glucose information

The table below lists the different types of meter blood glucose information that the system automatically records. The meter blood glucose information is calculated by using all of the automatic and manual meter blood glucose measurements that were received by the pump for the day, or number of days. This includes measurements that were manually entered into the pump during that time.

Item	Description
BG Average	The average meter blood glucose measurement.
Number BGs	The total number of meter blood glucose measurements taken.
Meter L/H	The lowest and highest blood glucose measurements sent by the meter (L is the lowest, and H is the highest). These can be used for calibration or non-calibration purposes.
Manual L/H	The lowest and highest manually entered blood glucose measurements (L is the lowest, and H is the highest). These can be used for calibration or non-calibration purposes.

To view data for a single day:

- Go to the DAILY TOTALS screen.
Main > Utilities > Daily Totals
- Select the day that you want to view.
 - To view totals for today, select **Today**.
 - To view totals for a previous day, select the date.
- Press **ACT**. The DELIVERY STATS screen shows the basal, bolus, and food intake information for the day you selected in step 2. This day is shown on the title line of the screen.
- Press  to view more information. See the table in the *Insulin delivery information* section in this chapter for descriptions of information shown here.



- 5 To view sensor glucose information for the same day you have selected in step 2, press **ACT** from the DELIVERY STATS screen. The SENSOR STATS screen appears, showing the sensor glucose information for the day you selected. The day you selected is shown on the title line of the screen.
- 6 Press  to view more information. See the table in the *Sensor glucose information* section in this chapter for descriptions of the sensor glucose information.
- 7 To view meter blood glucose measurement information for the same day you have selected in step 2, press **ACT** from the SENSOR STATS screen. The METER STATS screen appears, showing the meter blood glucose information for the day you selected. The day you selected is shown on the title line of the screen.
- 8 Press  to view all of the information in the METER STATS screen. See the table in this section for descriptions of the meter blood glucose information shown here.
- 9 Press **ACT** to toggle between the different STATS screens. Press **ESC** from any STATS screen to go to the DAILY TOTALS screen.


Continue viewing daily totals by following the steps above, or press **ESC** until the HOME screen appears.

The following procedure shows you how to view insulin delivery, food intake, sensor glucose and meter blood glucose data for a selected number of days.

To view data over a number of days:

- 1 Go to the DAILY TOTALS screen.
Main > Utilities > Daily Totals

The DAILY TOTALS screen shows **Daily Average** selected.
- 2 Press **ACT**. The DAYS TO AVERAGE screen flashes the number of days you selected the last time you used this feature.
- 3 Select the number of days that you want to average. For example, if you want to see the averages for the last seven days, select 7.
- 4 Press **ACT**. A message tells you that the pump is calculating the averages.
The DELIVERY STATS screen shows the insulin delivery and food intake information for the number of days you selected. The number of days is shown on the title line of the screen.
- 5 Press  to view more information. See the table in the *Insulin delivery information* section in this chapter for descriptions of the information shown here.
- 6 To view sensor glucose information for the same days you selected in step 3, press **ACT**. The SENSOR STATS screen appears, showing the sensor glucose information.
- 7 Press  to view more information. See the table in the *Sensor glucose information* section in this chapter for descriptions of the information shown here.

- 8 To view meter blood glucose measurement information for the same days you selected in step 3, press **ACT**. The METER STATS screen shows the meter blood glucose information.
- 9 Press  to view more information. See the table in the *Meter blood glucose information* section in this chapter for descriptions of the information shown here.

AUC calculation

The AUC (Area Under the Curve) feature measures how much and how long your sensor glucose measurements are outside the pre-defined AUC Limits. Check with your healthcare professional to determine the AUC settings that are best for you. The AUC Limits are different than the Glucose Limits. The AUC Limits are used to analyze the sensor glucose measurement data that is stored in your system. The closer the AUC values are to zero, the closer your SG values are to your AUC Limits. If your AUC values are zero, then your SG values are within your AUC Limits.

The default AUC Limits are:

- Low AUC Limit: 3.9 mmol/L (70 mg/dL)
- High AUC Limit: 10.0 mmol/L (180 mg/dL)

To set the AUC Limits:

- 1 Go to the AUC MENU.
Main > Utilities > Daily Totals > Calculate AUC
- 2 Select **AUC Limit**, then press **ACT**. The SET AUC LIMITS screen appears. The default LOW AUC Limit flashes 3.9 mmol/L (70 mg/dL).
- 3 Select your Low AUC Limit. The number must be between 2.2 and 22.2 mmol/L (40 and 400 mg/dL).
- 4 Press **ACT**. The default High AUC Limit flashes 10.0 mmol/L (180 mg/dL).
- 5 Select your High AUC Limit. The number must be between 2.2 and 22.2 mmol/L (40 and 400 mg/dL). It can be equal to, but not lower than, the AUC Low Limit.

6 Press **ACT**. The AUC MENU screen shows the AUC Limit you selected.

To view AUC data for a single day:

1 Go to the AUC MENU.

Main > Utilities > Daily Totals > Calculate AUC

2 Select the day that you want to view.

- To view totals for today, select **Today**.
- To view totals for a previous day, select the date.

3 Press **ACT**. AUC High and AUC Low values appear for the selected day.

- **AUC High** is the average of AUC for sensor glucose values above the AUC High Limit for that day.
- **AUC Low** is the average of AUC for sensor glucose values below the AUC Low Limit for that day.

To view AUC data over a number of days:

1 Go to the AUC MENU.

Main > Utilities > Daily Totals > Calculate AUC

2 Select **AUC #Days**, then press **ACT**.

3 The DAYS TO AVERAGE flashes the number of days to average. Select the number of days you want. For example, if you want to see the averages for the last seven days, select 7.

4 Press **ACT**. A message tells you that the pump is calculating averages.

5 The Day Average screen shows the data for the number of days you selected. The number of days is shown on the screen. For example, if you have selected seven days, the screen title will be 7Day Average. AUC High and AUC Low appear below the title.

- **AUC High** is the average value of AUC for sensor glucose values above the AUC High Limit for the number of days selected.
- **AUC Low** is the average value of AUC for sensor glucose values below the AUC Low Limit for the number of days selected.

Personal reminders

Alarm clock

The alarm clock is a feature that allows you to set a maximum of eight daily reminders for various events. The factory setting for this feature is off. The alarm clock can be useful to remind you when to check your blood glucose, eat, bolus, etc. When the alarm clock goes off, the ALARM CLOCK screen appears. Press **ESC**, then **ACT** to clear the reminder.

- 1 Go to the ALARM OPTION screen.
Main > Utilities > Alarm Clock
- 2 Select **On/Set**. Press **ACT**.
- 3 Select **Add Alarm**. Press **ACT**.
- 4 Enter the hour (flashing). Press **ACT**.
- 5 Enter the minutes (flashing). Press **ACT**.
- 6 Repeat steps 3 through 5 to program additional alarm times. Exit the menus when you are done.

Remote control option

The factory setting for this feature is off. You may want to explore the remote option after you have become completely familiar with the basic functions of your pump. It is important that you consult with your healthcare professional before using this feature. Remote controls can be purchased from Medtronic Diabetes.

Refer to the remote control user guide for operating instructions.

NOTE: *The use of RF (radio frequency) devices with the pump reduces battery life.*

To use the remote control, these pump settings must be programmed:

- Remote Options = On
- Remote control ID code entered in pump (code is on back of remote)
- Easy Bolus = On

Turn on remote control option

WARNING: If there is a **LOW BATTERY** condition, the pump will not receive signals from the remote. To ensure the pump communicates with the remote control, make sure the pump does not have a low battery. (Replacing the low battery with a new battery will restore remote control function.)

- 1 Go to the REMOTE OPTION screen.
Main > Utilities > Connect Devices > Remotes
- 2 Select **On**, then press **ACT**.
- 3 The REMOTE ID MENU screen appears. Add, delete or review your remote ID as described in the next section. Exit the menus if you are done.

Add, delete, review remote control IDs

Each remote control has its own unique ID. Up to three different remote control IDs can be programmed in your pump. You must turn on the remote option to add, delete or review the remote control IDs programmed in your pump.

- 1 **Adding a remote ID**
 - a. Select **Add ID** from the REMOTE ID MENU screen, then press **ACT**.
 - b. Enter each of the six ID numbers found on the back of the remote control. Press **ACT** after each entry. After you set the last number of the ID, the screen will return to the REMOTE ID MENU.
- 2 **Deleting a remote ID**
 - a. Select **Delete ID** from the REMOTE ID MENU screen, then press **ACT**.
 - b. Select the remote ID that you want to delete, then press **ACT**. The selected ID is now deleted.
- 3 **Reviewing remote IDs**
 - a. Select **Review ID** from the REMOTE ID MENU screen, then press **ACT**.
 - b. The programmed IDs will show in the REVIEW REMOTE ID screen.
- 4 Exit the menus when you are done.

Other Devices option

The factory setting for this feature is **Off**.

Other devices may be available to be linked to your pump. If available, refer to the appropriate device user guide for programming and operating instructions.

NOTE: The use of RF (radio frequency) devices with the pump reduces battery life.

Block feature

Block restricts access to pump programming. The factory setting for this feature is off. Block is an important safety feature if the pump user requires someone else to maintain complete control of pump operation. When block is on, the remote control is used to deliver a bolus and suspend/resume the pump. Direct pump programming is limited to suspend, block, and selftest. You can, however, still view status-type screens such as STATUS, BOLUS HISTORY, BASAL REVIEW, and DAILY TOTALS. Discuss what uses and settings are best for you with your healthcare professional. (You can order the remote control from Medtronic Diabetes.)

Turning block on

NOTE: The vibrate alert-type is disabled when block is on.

- 1 Go to the BLOCK OPTION screen.
Main > Utilities > Block
- 2 Select **On**, then press **ACT**. The HOME screen will appear with an open circle. The Block option is now on and the pump is in Special mode. Exit the menus.



Example 1: Block

Nicholas is a very active toddler who wears a Paradigm pump. His parents don't want to worry that he will play with the pump and accidentally change his programmed settings. They simply activated the Block feature, and now, except for the Suspend and Self-Test, no other features are active when using the pump buttons. When Nicholas needs a bolus, his parents and caregivers simply program it with the Remote Control.

Example 2: Block

Oscar is an elderly man with diabetes who needs assistance with all of his daily living activities. He needs his caregiver or family member to help him with his pump as well. To be sure that Oscar does not change any pump settings, his family programmed his Paradigm pump with the Block feature turned on. They use the Remote Control to give him his boluses when he needs them.

Lock keypad feature

Lock keypad prevents accidental pump keypad presses. You can only press the  to view the STATUS screen, and the  button to turn on the backlight. The remote control can be used to give a bolus or put the pump into Suspend.


A locked keypad is automatically unlocked during the following:

- Battery insertion
- Alarms
- Alerts

Locking the keypad

- 1 Go to the UTILITIES MENU.
- 2 Select **Lock Keypad**, then press **ACT**.
- 3 Press **ACT** again to lock the keypad. The KEYPAD LOCKED screen appears with instructions on how to unlock the keypad.

Unlocking the keypad

Press the  button and the  button at the same time. Keypad Unlocked screen appears.

Selftest

Selftest is a safety utility that allows you to check if your pump is operating properly. This self-diagnostic feature can be used for maintenance or to check your pump if it operates unusually. During selftest, your pump will automatically run internal tests, including a check for proper operation of the beep and vibrate modes. The selftest is additional to the routine tests that run independently while the pump operates.

Contact your local help line or representative if any of the tests do not occur as described here.

NOTE: *If the pump detects a condition such as low battery, the selftest will not finish. A message will appear to show the condition that caused the test to stop.*

- 1 Go to the UTILITIES MENU.
Main > Utilities > Selftest
- 2 Select **Selftest**, then press **ACT**.

- 3 Periodically, you will hear beeps as different mechanisms in the pump are being tested. As part of the selftest, the pump will do these tests:
 - a. Screen Test:
The screen will appear all black.
 - b. Selftest:
The pump will count down from 10.
 - c. Tone Test:
You should hear beeps.
 - d. Vibrate Test:
You will feel vibrations.
- 4 After the selftest is finished, TEST COMPLETE screen appears. The screen will return to the UTILITIES MENU, then to the HOME screen.

User settings

The user settings function allows you to save, restore, and clear all pump settings. You can also view a listing of the dates and times of all recent user settings operations you have done. The Save Settings feature lets you keep a set of pump settings that you can restore to your pump if it is cleared or you need to go back to these settings for any reason.

When you clear your pump, the pump settings are restored to the factory defaults, and you must either use Restore Settings, if you have saved a set of pump settings, or reprogram all your settings before you can use your pump again. The pump does not clear the internal pump memory.

WARNING: Do not clear your pump settings while it is connected to your body.

CAUTION: Do not clear your pump settings unless directed by your healthcare professional or a Medtronic Diabetes representative. If you clear your pump settings, it will be necessary to reprogram all your personal pump settings as directed by your healthcare professional. Additionally, you will have to rewind your pump.

Saving the settings

Do these steps to save your current pump settings:

- 1 Go to the UTILITIES MENU screen and select **User Settings**.
Main > Utilities > User Settings
- 2 Hold down **↶B** and press **ACT**.
- 3 The USER SETTINGS screen is displayed with **Save Settings** highlighted. Press **ACT**.
- 4 If this is the first time you have saved pump settings, go to step 5. If you have previously saved pump settings, a message displays, indicating the date of your last pump settings save. Read the instructions on the screen, then press **ACT** to save your current settings. You can press **ESC** if you want to cancel the save.
- 5 The SETTINGS SAVED message displays to confirm that your current pump settings have been saved. Exit the menus.

Restoring the settings

Do these steps to restore the most recent pump settings you have saved to your pump:

- 1 Go to the UTILITIES MENU screen and select **User Settings**.
Main > Utilities > User Settings
- 2 Hold down **↶B** and press **ACT**.
- 3 The USER SETTINGS screen is displayed. Select **Restore Settings**, then press **ACT**.
- 4 A message displays, giving you the option of restoring the pump settings that were saved on the given date and erasing the ones currently on your pump. Read the instructions on the screen then press **ACT** to restore the settings. You can press **ESC** if you want to cancel the restore.
- 5 The SETTINGS RESTORED message displays to confirm that your current pump settings have been replaced with the pump settings you had saved on the given date. Exit the menus and check your pump settings to verify the restore.

Clearing the settings

Take the following steps only if you want to clear your pump to factory default settings.

WARNING: Do not clear your pump settings while it is connected to your body.

CAUTION: Do not clear your pump settings unless directed by your healthcare professional or a Medtronic Diabetes representative. If you clear your pump settings, it will be necessary to reprogram all your personal pump settings as directed by your healthcare professional. Additionally, you will have to rewind your pump.

- 1 Go to the UTILITIES MENU screen and select **User Settings**.
Main > Utilities > User Settings
- 2 Hold down **⬅B** and press **ACT**.
- 3 The USER SETTINGS screen is displayed. Select **Clear Settings**, then press **ACT**.
- 4 A CONFIRM message displays, giving you the option of clearing your pump settings or not. Select **YES** and press **ACT** to clear the settings. You can press **ESC** if you want to cancel.
- 5 You will see a RESET screen, and then the pump will go through various screens while it restarts. After the pump clears all of your settings, the screen will go to the TIME/DATE SETUP screen.
- 6 Reset the time and date as described in *Setting the time and date* section in the *Basic programming* chapter.
- 7 After you set the time and date, you must rewind your pump. Refer to the *Rewinding your pump* section in the *Starting on insulin* chapter for instructions. Remember, all your settings have been cleared and you must either restore or reprogram all your settings.

History

If you want to view a listing of the dates and times of all recent User Settings operations you have done on your pump, such as saves and restores, take the following steps:

- 1 Go to the UTILITIES MENU screen and select **User Settings**.
Main > Utilities > User Settings
- 2 Hold down **⬅B** and press **ACT**.
- 3 The USER SETTINGS screen is displayed. Select **History**, then press **ACT**.
- 4 The SETTINGS HISTORY screen displays, giving you a list of all the dates and times of your most recent User Setting operations. Scroll down to view the entire history. When you are done looking at the history, press **ESC** to exit the menu.

Therapy management software

Medtronic provides software to optimize your pump use. Contact your local representative for more information.

CareLink® Personal Software

Your pump provides continuous glucose values so that you can track glucose concentration patterns and possibly identify episodes of low and high blood glucose. It also stores glucose data so that it can be analyzed to track patterns or downloaded to CareLink Personal for analysis of historical glucose values.

CareLink Personal is a Web-based system designed to help you manage your diabetes. It has many key features:

- It copies (uploads) data from your devices: insulin pump and supported blood glucose meters.
- Device data is stored on a centralized database.
- There is an online logbook where you can record self-reported information, such as how many carbohydrates you consumed.
- Uploaded data and other information stored on the system can be viewed through several different types of treatment reports.
- Access to data and personal information on this system is secure.

Reports

Reports show you data gathered from your device(s), and from what you directly entered into the system. Examples of the report data you can get include:

- Pump operation, such as alarm settings and maximum allowed bolus amount.
- Insulin delivery rates, types and times.
- Blood glucose reading levels and times.
- Sensor glucose reading levels and times.
- Bolus Wizard or online logbook entry data.

Reports are displayed in PDF format. They can be viewed online, saved, or printed.

You can share your reports with your healthcare professional. These reports can help you and your healthcare professional to discover trends and other information. This can lead to improved therapy management for greater control.

Logbook

The logbook allows you to enter the following self-reported information:

- Carbohydrates consumed
- Exercise activity
- Urine ketone results
- HbA1c results
- Infusion set changes

The data from your logbook is used in the reports you create. You can see how, for example, your carbohydrate intake compares to blood glucose levels for the same day or time. You can also use it as an online diary of your diabetes management.

Insulin pump therapy follow-up

We hope that you are now comfortable using the pump and your blood glucose values have improved through insulin pump therapy. Diabetes management requires much more than blood glucose control. You need to take care of your complete physical and mental health. This includes seeking treatment for any condition both directly related to and not related to diabetes. The following recommendations apply to general diabetes as well as insulin pump therapy follow-up. Remember, your healthcare professional is your best resource for successful diabetes management.

Recommended follow-up

Everyday

- Check blood glucose four to six times a day and always before bed.
- Test before driving and have a fast-acting carbohydrate with you when you drive.
- If your blood glucose is above 13.9 mmol/L (250 mg/dL) twice in a row, take an injection and change the infusion set.

Every month

- Review DKA prevention guidelines.
- Check 03:00 blood glucose at least once during the month.
- Check two-hour post-meal blood glucose for all meals on a given day.

Every three months

- Visit your healthcare professional, even if you feel well and your blood glucose values are within target range.
- Review your blood glucose log and insulin pump settings with your healthcare professional.
- Make sure you have an HbA1c test done.

Laboratory tests

- Test for HbA1c four or more times a year
- Test for cholesterol, HDL, LDL, triglyceride yearly
- Test for microalbuminuria yearly

Every visit

- Blood pressure check
- Foot exam
- Review goals for blood glucose, meal plan and exercise

Annually

- Dilated eye exam by a qualified ophthalmologist
- Annual flu shot
- Regular dental visits
- Nerve function tests
- EKG test over age 35
- Prostate exam for men, breast exam for women
- Diabetes education review
- Replace Glucagon Emergency Kit (new prescription from your healthcare professional)

Troubleshooting, alerts and alarms

This chapter is designed to help you understand the messages that your pump is giving you when you get an alarm or an alert message. The procedures in the beginning of this chapter are to be used when you get specific alarms or cover a likely condition that might happen. A list of alarms is provided at the end of this chapter.

NOTE: *It is recommended that you read your warranty statement included with your pump for information on what is covered during your warranty period.*

My pump has a NO DELIVERY alarm

When a NO DELIVERY alarm occurs, it means that insulin delivery has stopped. Your pump is not broken, but it has detected that something is preventing insulin from being delivered. Do the following steps:

- 1 Check to see if your infusion set has become dislodged or kinked.
- 2 Check your blood glucose and follow your safety guidelines for treating high glucose if necessary.
- 3 Clear the alarm by pressing **ESC** and **ACT**. A screen appears with two choices: **Resume** and **Rewind**.
- 4 Make sure that there is insulin in your reservoir.
 - a. If you identified and fixed a kink in your tubing, select **Resume**.
 - b. If the reservoir is empty, select **Rewind** and change your reservoir and infusion set.
- 5 Monitor your blood glucose closely.
- 6 If you have followed these steps and receive another NO DELIVERY alarm, call your local representative.

What happens if I leave the pump battery out too long?

If you leave the pump battery out too long (more than five minutes) you may receive a BATT OUT LIMIT alarm message when you install the new battery. Do the following steps:

- 1 Set your pump clock to the correct time, date, and year.
- 2 Check to make sure that all your settings, such as basal rate, are set as desired. If need be, reapply your last saved settings to the pump by using the **Restore Settings** option under **User Settings** in the UTILITIES MENU (you can only use this option if you have previously saved your pump settings). See the *Restoring the settings* section in the *Utilities* chapter.
- 3 Check the ALARM HISTORY screen and the STATUS screen for any alarms and/or alerts that may still need attention.

If you leave the battery out for more than three or four days, you may receive an A17 and A47 alarms when you install a new battery. Do the following steps:

- 1 Clear the A17 and A47 alarms.
- 2 Set your pump clock to the correct time, date and year.
- 3 Check to make sure that all your settings, such as basal rate, are set as desired. If need be, reapply your last saved settings to the pump by using the **Restore Settings** option under **User Settings** in the UTILITIES MENU (you can only use this option if you have previously saved your pump settings). See the *Restoring the settings* section in the *Utilities* chapter.
- 4 Check the ALARM HISTORY screen and the STATUS screen for any alarms and/or alerts.

Why doesn't my pump battery last very long?

A short battery life does not mean that something is wrong with your pump. Battery life in your pump is variable and based on the conditions below:

- The brand of battery you use (we recommend Energizer).
- The storage and/or handling of the battery before use (avoid high or low temperatures).
- The usage of pump in cold temperatures; this may shorten the battery life.
- The usage of your pump, such as how often the buttons are pushed, the number of alerts/alarms and set changes.
- The amount of insulin the pump is delivering.
- The usage of some features. The backlight, vibrate, sensor, remote control and meter options decrease battery life.

What is a CHECK SETTINGS alarm?

This alarm occurs after an E alarm or after you clear your pump. It is advising you to make sure that all your settings are correct. A CHECK SETTINGS alarm occurs after any of these actions:

- All user settings were cleared (set back to their defaults) because there was an E-error alarm,
- The Clear Settings function was performed, or
- After you rewind when practicing without a reservoir when you first get your pump. In this case, it is just a reminder to make sure all your settings are programmed before you begin using the pump with insulin.

My screen appears distorted

The screen may appear distorted or have a rainbow-like appearance if you are wearing polarized sunglasses, are in bright sunlight, or in extreme high or low temperatures. If your screen appears distorted:

- Take off your sunglasses.
- Move into the shade.
- Make sure your pump is not in direct heat (next to a heater) or cold (worn on the outside of your clothing on a very cold day).
- **Do not return the pump:** this is a normal property of this type of screen on any device.

The pump is asking me to rewind

This is normal after any of the following:




- 1 Any E-alarms,
- 2 The Clear Settings function,
- 3 Or, a NO DELIVERY alarm (during the Reservoir + Set sequence).

My bolus stopped

The Bolus Stopped error can occur if the battery cap is loose or the pump was bumped or dropped during a bolus. It can also happen if the pump receives a static shock. As a safety measure, the pump stops the bolus when this happens.

- 1 If you dropped your pump, visually inspect it to make sure that it is not damaged in any way.
- 2 Review your bolus history and reprogram the remaining bolus, if needed.

My pump buttons are not acting right during a bolus

If the ,  or  buttons are pressed and held down while a bolus is being delivered, the screen will freeze on that amount. Once the button is released, the units will ramp up to the amount delivered so far. Pressing and holding down the button will not stop the delivery of a bolus.

My pump will not display my blood glucose reading from my meter

- 1 Make sure you are using the correct meter (a blood glucose meter supported by MWT1 technology). Your pump will communicate with this meter only.
- 2 Make sure your meter's RF (radio frequency) feature is on and working correctly.
- 3 Make sure the meter option in your pump is set to **On** and you have set the meter ID number correctly on the pump.
- 4 Make sure your pump does not have a low battery alert condition.
- 5 Make sure the meter is within 1.2 meters (4 feet) of the pump without anything in between such as another person, a wall, etc.
- 6 Make sure there is no RF (radio frequency) interference from other electronic devices that could prevent communication. These devices can include some cell phones, cordless phones, televisions, computers, radios, other Paradigm pumps, meters and pump remote controls. To restore communication, simply move away from these other types of devices, or turn them off.
- 7 Your pump will not show another reading. Make sure the pump is idle and the HOME screen is blank.
- 8 If your pump still does not receive your blood glucose reading from the meter, use the up/down buttons to manually enter your blood glucose (in the ENTER BG screen).

I dropped my pump

Take care to protect your pump from being dropped.

- 1 Check that all connections are still tightly in place.
- 2 Check the LCD, keypad and pump case for cracks or damage.
- 3 Check infusion set, including the tubing connector and tubing for cracks or damage.
- 4 Review the status screen, basal rates and other pump settings.
- 5 Perform the Selftest procedure located in the UTILITIES MENU.
- 6 Call your local Medtronic Diabetes representative for assistance.

I submerged my pump in water

Your pump is designed to resist accidental contact with water. Do not submerge in water during bathing, swimming, or other water activities.

- 1 Pat the outside of the case until dry.
- 2 Open the reservoir compartment and check the compartment and reservoir for water. If wet, dry it completely within 10 minutes of exposure to water. Exposure to liquids, including water or insulin can corrode the mechanism.
- 3 Dry the reservoir completely—do NOT place a wet reservoir in the pump.
- 4 Do not use hot air to dry your pump. This may damage your pump's internal electronics.
- 5 Check the battery compartment and the battery—if wet, let them dry completely before using the pump.
- 6 Perform a selftest.

I cannot get to the USER SETTINGS screen

If you press **ACT** when you have **User Settings** highlighted, the following message appears: This feature is not normally accessible. To access, consult your User Guide.


To get to the USER SETTINGS screen:

- 1 Go to the UTILITIES MENU screen.
Main > Utilities
- 2 Select **User Settings**. Hold down **⬅B** and press **ACT**.
- 3 See the *User settings* section in the *Utilities* chapter for information about the menu options.

Alerts

Your pump has a sophisticated network of safety checks and systems. If it detects an unusual condition that requires your immediate attention, it beeps or vibrates periodically to alert you. The pump goes into Special mode (an open circle displays), and the backlight illuminates.

Question: Why are alerts important?

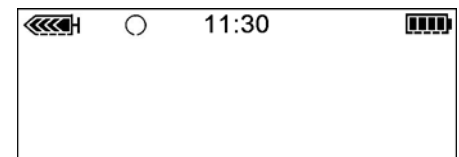

Answer: Your pump monitors activities and will alert you if a Special mode is active. Some alerts are a normal part of pump therapy, such as an active temporary basal. There are alerts that indicate a condition that is outside normal pump activity. For example, your pump notifies you with an alert when you need to replace the reservoir (LOW RESERVOIR) or replace your pump battery (LOW BATTERY).

It is important that you understand how your pump behaves when you do not clear an alarm or an alert:

- Alarm — If you do not clear an alarm, the pump does not display any subsequent alarms or alerts until the first alarm is cleared. After you clear the first alarm, the pump displays any subsequent alarms, followed by any alerts in order of their importance.
- Alert — If you do not clear an alert, the pump does not display any subsequent alerts until the first alert is cleared. However, the pump displays an alarm even if you do not clear the alert. After you clear the first alert or alarm, the pump displays any subsequent alarms or alerts in order of their importance.

Check your pump and test your blood sugar. Follow the guidelines prescribed by your healthcare professional to manage your diabetes as needed.

Example:

John's pump displays a BUTTON ERROR alarm, but he does not clear it. During the time that his pump displays the BUTTON ERROR alarm, his pump receives a LOW RESERVOIR alert. A few minutes later, his pump receives a NO DELIVERY alarm. John's pump does not display the LOW RESERVOIR alert or the NO DELIVERY alarm until he clears the first alarm. After John clears the BUTTON ERROR alarm, his pump displays the NO DELIVERY alarm. After he clears the NO DELIVERY alarm, his pump displays the LOW RESERVOIR alert.

What to do

When your pump beeps or vibrates notifying you that an alert condition exists:

- 1 Read and follow the instructions on the screen. Press **ESC**, **ACT** to clear the alert message.
- 2 Check the STATUS screen to determine what caused the alert.
- 3 If the condition is due to a low battery, replace the battery.
- 4 If the condition is due to a low reservoir, monitor the reservoir volume frequently and change the reservoir when appropriate. Make sure you have a new reservoir, infusion set and vial of insulin with you.

Pump alert conditions

Listed below are alerts that indicate a condition that is outside normal pump activity.

LOW RESERVOIR

You can program the pump to sound an alert when either a specified number of units remains or a specified amount of time remains before the reservoir will be empty.

LOW BATTERY

If you get this alert, **do not go to sleep** without replacing the battery. The backlight, the remote, and meter functions are disabled during a LOW BATTERY condition. If the alert type is set to Vibrate, the pump will change to the audio alert Beep-medium. Clear (**ESC**, **ACT**) this alert before you replace your battery.

Sensor alert conditions

Listed below are the alerts that you may encounter while using the sensor feature of your pump, and how to resolve the alert condition.

WEAK SIGNAL

Reason: Occurs when the pump does not receive data from the transmitter for a pre-defined period of time (as set in Weak Signal).

Action: Move the pump closer to the transmitter or move the transmitter and the pump to the same side of your body.

ALERT SILENCE

Reason: This alert screen displays when one or more glucose alerts occur during the time that the Alert Silence feature is turned on.

Action: Press **ESC**, then **ACT** to clear the alert. Check the SENSOR ALERT HISTORY screen for the recorded glucose alert information. See the *Sensor alert history* section in the *Using your sensor* chapter for details.

MISSED BOLUS

Reason: You have not delivered a bolus during the time period you specified in the Missed Bolus Reminder feature.

Action: Set a meal bolus if needed.

LOST SENSOR

Reason: The pump has not received a signal from the transmitter. **Do not disconnect the transmitter from the sensor.**

Action:

- 1 Make sure that the sensor is inserted correctly.
- 2 Check the REVIEW SETTINGS screen to make sure that the transmitter ID entered into the pump matches the ID on your transmitter:

Main > Sensor > Review Settings

- 3 Check the transmitter and sensor connection. Touch the inserted sensor at the back of the assembly to prevent movement and push the transmitter firmly:
 - a. If you hear a click, wait 20 seconds and look for a green light on the transmitter flash for 10 seconds to confirm a good connection. If you see the green light, the alert was due to the transmitter and the sensor not being connected.
 - b. If you hear a click, but do not see a green light on the transmitter, check to make sure that the transmitter is charged.
 - c. If you do not hear a click when you check the connection, the alert was due to a transmission problem. Bring the pump closer to the sensor and transmitter.
- 4 Use the **Find Lost Sensor** function to find your sensor (see the *Troubleshooting sensor features* section in this chapter):

Main > Sensor > Link to Sensor > Find Lost Sensor

LOW TRANSMTR

Reason: Occurs when the transmitter battery is close to running out of power. This alert will repeat daily at noon while this condition exists. The transmitter will continue sending sensor signals for several hours and may last for several days until the battery becomes depleted. Be prepared to recharge your transmitter immediately when its battery becomes depleted.

Action: Recharge your transmitter as soon as possible. See the *Charging the transmitter* section in the *Sensor features* chapter.

CHARGE TRANSMTR

Reason: The transmitter battery is depleted.

Action: Recharge your transmitter immediately. See the *Charging the transmitter* section in the *Sensor features* chapter.

CHANGE SENSOR

Reason: You may receive this alert after you receive two Cal Errors in a row, without the Cal Errors, or when initializing your sensor.

Action: If the alert resulted from two Cal Errors in a row, replace the sensor. If the alert happened without two Cal Errors, use the transmitter's tester to make sure your transmitter is working properly. If you received this alert during initialization, you may be able to resolve this alert without replacing your sensor. Call your local representative for further assistance.

SENSOR END

Reason: The sensor has reached the end of its life.

Action: Replace your sensor. The sensor has a maximum life of 144 hours (6 days). The 144-hour life span of the sensor begins when the pump receives the first METER BG NOW alert.

CAL ERROR

Reason: An error occurred when entering a new meter BG measurement to calibrate the system. Some possible causes are:

- An incorrect blood glucose number was entered from the meter into the pump.
- The entered blood glucose measurement was not current.
- Your blood glucose is rising or falling rapidly.
- The sensor needs more time to stabilize after being inserted.
- The sensor is no longer reading the sensor glucose correctly.

Action: Follow the guidelines in the *Calibrating the sensor* section in the *Sensor features* chapter to enter a new meter BG for calibration. If you receive a Cal Error on your second calibration, a CHANGE SENSOR alert will occur. Call your local representative if you have questions.

METER BG NOW

Reason: A meter BG measurement is needed immediately to calibrate your sensor so that you can continue receiving sensor glucose readings.

Action: Take and enter a meter blood glucose measurement. Follow the guidelines in the *Calibrating the sensor* section in the *Sensor features* chapter.

METER BG BY

Reason: A meter BG measurement must be entered by the time that is shown to calibrate your sensor. The METER BG BY alert is also known as the CAL REMINDER alert.

Action: Enter a blood glucose measurement to avoid a METER BG NOW alert.

LOW X.X MMOL/L (XX MG/DL)

Reason: The sensor glucose number is equal to or lower than your low glucose limit. If you do not set a Low Glucose Limit and do not set the glucose alerts ON, then you will not get a Low sensor glucose alert. The pump plays four consecutive tones, in falling pitch, if an audible beep has been selected as the alert type.

HIGH XX.X MMOL/L (XXX MG/DL)

Reason: The sensor glucose number is higher than or equal to your high glucose limit. If you do not set a High Glucose Limit and do not set the glucose alerts ON, then you will not get a High sensor glucose alert. The pump plays four consecutive tones, in rising pitch, if an audible beep has been selected as the alert type.

LOW PREDICTED

Reason: The sensor glucose measurements may **reach or go below** your Low Glucose Limit in the length of time you selected for the Low Predictive alert. The pump plays three consecutive tones, in falling pitch, if an audible beep has been selected as the alert type.

HIGH PREDICTED

Reason: The sensor glucose measurements may **reach or go above** your High Glucose Limit in the length of time you selected for the High Predictive alert. The pump plays three consecutive tones, in rising pitch, if an audible beep has been selected as the alert type.

RISE RATE

Reason: The sensor glucose measurements are rising at a rate that is equal to or faster than the Set Rise Rate Limit you selected for the alert. The pump plays two consecutive tones, in rising pitch, if an audible beep has been selected as the alert type.

FALL RATE

Reason: The sensor glucose measurements are falling at a rate that is equal to or faster than the Set Fall Rate Limit you selected for the alert. The pump plays two consecutive tones, in falling pitch, if an audible beep has been selected as the alert type.

SENSOR ERROR

Reason: The sensor signals are either too high or too low.

Action: You do not need to change the sensor. Clear the alert. If the alert persists, test your transmitter with the tester.

Alarms

Your pump has a sophisticated network of safety checks and systems. If the safety network detects anything unusual, your pump notifies you of conditions that require your immediate attention. The backlight illuminates the pump screen and the alarm message displays on the screen.

It is important that you understand how your pump behaves when you do not clear an alarm or an alert. See the *Alerts* section in this chapter for an explanation and an example.

Check your pump and test your blood sugar. Follow the guidelines prescribed by your healthcare professional to manage your diabetes as needed.

NOTE: The STATUS screen shows any alarms and alerts that are active.


Question: Why are alarms important?

Answer: Your pump monitors activities and notifies you if there is an unusual pump status or your attention is required. When an attention alarm is active, INSULIN DELIVERY IS STOPPED and immediate operator interaction is required.

If the vibrate mode is on, all alarms and alerts start as vibrations and then change to beeps. For your safety, if there is no response after 10 minutes, the beeps change to a siren. The siren continues every minute until the alarm is cleared.

What to do

When an alarm is triggered, the pump goes into Attention mode and an alarm message shows on the screen. The pump then defaults to the HOME screen. Do these steps when you get an alarm:

- 1 **View the alarm:** From the HOME screen, press any button to see the alarm message.
- 2 **Read all of the alarm text.** There are instructions on how to fix the alarm condition. (Press  to read more text, if available.)
- 3 **Clear the alarm.** Press ESC then ACT after you read the alarm instructions.
- 4 The HOME screen appears.
- 5 **Follow the instructions** that appeared with the alarm to fix the alarm condition.
- 6 **Check your settings.** Check your time, date, basal and other settings to make sure they are correct.

Alarm conditions

Alarms put the pump in Attention mode.



A (ALARM)

This alarm shows the letter *A* followed by two numbers. A-alarms cause all insulin delivery to stop. Your pump settings are retained. If this alarm repeats often, call your local Medtronic Diabetes representative for assistance.



AUTO OFF

Alerts you that no buttons were pressed during the time limit you set for the AUTO OFF DURATION feature, and so insulin delivery has been stopped.



BATT OUT LIMIT

Occurs if the battery has been out of the pump for more than five minutes. Verify that the pump time and date are correct. If not correct, go to the UTILITIES MENU and reset the date and time.



BOLUS STOPPED

If this alarm occurs, it is very important to check bolus history to review how much of the bolus was actually delivered. If necessary, reprogram a bolus with the amount not yet delivered.



BUTTON ERROR

Occurs if a button has been continually pressed for more than three minutes.

CHECK SETTINGS

When this alarm is active, you should check and/or reprogram your pump settings, including the time/date.

E (ERROR)

An E-error alarm will show the letter *E* followed by two numbers. E-alarms cause all insulin delivery to stop, the pump resets, and all your settings are cleared. After receiving this alarm, note the error number then call your local Medtronic Diabetes representative for assistance.

EMPTY RESERVOIR

There is no insulin in the reservoir. Change the reservoir immediately.

FAILED BATT TEST

The pump tests the voltage of each battery installed. This test ensures a battery with low voltage is not used. If the battery does not have enough voltage, this alarm will occur. The pump will not function and the battery must be replaced. (Always make sure that you install a NEW battery into the pump.)

FINISH LOADING

You have not completed filling the infusion set with insulin. Clear the alarm. This resumes basal delivery. See the *Filling the cannula* section in the *Starting on insulin* chapter to fill the infusion set cannula with insulin. If your infusion set has a needle cannula, press **ESC** to skip this step.

MAX FILL REACHED

If you use more than 30 units of insulin to fill the infusion set tubing, this message appears. Press **ESC**, **ACT** to clear the message. See the *Filling the tubing* section in the *Starting on insulin* chapter for more information.

MAX DELIVERY

This alarm alerts you when you have taken more insulin than expected based on Maximum bolus and Maximum basal rates.

MOTOR ERROR

Insulin delivery has stopped. This alarm occurs if your pump detects a motor error.

NO DELIVERY

Insulin delivery has stopped. This alarm will occur if your pump detects a blockage or if the reservoir is empty. Refer to the first page of the Troubleshooting chapter for additional information.

NO RESERVOIR

The reservoir is not inserted correctly or no reservoir has been inserted.

OFF NO POWER

The battery is dead. Replace battery immediately. Follow the directions on the screen. Check to make sure that the time is correct on the screen. Reset the time if necessary.

RESET

The Reset alarm triggers when pump settings are cleared because of one of these reasons:

- Pump settings were cleared (Clear Settings function) and settings have not been reprogrammed.
- A download attempt from the PC is incomplete. (The download function is applicable to the optional software feature. Refer to the software user guide for more information.)

WEAK BATTERY

The pump tests the voltage of each battery installed. If the battery voltage is less than full strength, this alarm may occur. The pump will operate normally, but the battery life will be shorter than expected. Always make sure to install a new battery in the pump.



LOW SUSPEND

When your sensor glucose measurement is below the Low Suspend Limit, the LOW SUSPEND alarm screen appears, and all insulin delivery stops. A siren starts immediately to get your attention or to get help from others if you are unable to respond to the alarm. The pump vibrates and sounds the siren until you or someone else pushes a button on the pump.

If you are unable to respond to the alarm:

After two minutes of the LOW SUSPEND alarm, a new screen appears with the following message: **I have diabetes, call for emergency assistance. To clear ESC then ACT.**

After two hours, the pump resumes delivering insulin at the **standard or pattern** basal rate that is programmed to be active for that time. **This cancels any programmed temporary basal insulin delivery.** The following message appears as the insulin delivery resumes: **I have diabetes, call for emergency assistance. Low Suspend canceled, check Status. ESC-ACT to cont.**

If you are able to respond to the alarm:

- 1 Press **ESC** then **ACT** to clear the alarm.
 - If you clear the alarm within two minutes, the LOW SUSPEND/CHECK BG screen appears.
 - If you clear the alarm after two minutes, the pump returns to the HOME screen. Press **ACT**, and the LOW SUSPEND/CHECK BG screen appears.
- 2 You have two options: to keep the insulin delivery suspended (**Suspend**) or to resume insulin delivery (**Restart Basal**). See the following sections for details.

To resume insulin delivery after clearing the alarm:

- 1 Select **Restart Basal** in the LOW SUSPEND/CHECK BG screen, then press **ACT**.
- 2 The pump resumes delivering insulin at the basal rate (standard, pattern or temporary) that is programmed to be active for that time. The following message appears as insulin delivery resumes: **Basal delivery is resumed. Check BG. Press any key to cont.**
- 3 Check your BG and treat as appropriate.
- 4 Press any button to continue using your pump.

To keep your insulin delivery suspended after clearing the alarm:

- 1 Press **ACT** with the **Suspend** option highlighted in the LOW SUSPEND/CHECK BG screen.
- 2 The SUSPEND screen appears, and the pump soon returns to the HOME screen.
- 3 Check your BG and treat as appropriate. After your blood glucose has returned to an acceptable measurement, see the following section to restart delivering your pump's programmed basal.

NOTE: Two hours after you choose to keep your insulin delivery suspended, the pump automatically resumes delivering insulin at the basal rate (standard, pattern or temporary) that is programmed to be active for that time. The pump returns to the HOME screen as insulin delivery resumes. When you press **ACT** from the HOME screen, the following message appears: **Basal delivery is resumed. Check BG.** Press any key to cont.

To resume insulin delivery after you suspend your pump:


- 1 Press **ACT** from the HOME screen when you are ready to restart your basal delivery.
- 2 The LOW SUSPEND/CHECK BG screen appears. Select **Restart Basal**, then press **ACT**.
- 3 The pump resumes delivering insulin at the basal rate (standard, pattern or temporary) that is programmed to be active for that time. The following message appears as insulin delivery resumes: **Basal delivery is resumed. Check BG. Press any key to cont.**
- 4 Press any button to continue using your pump.

Testing the transmitter

The tester acts like a sensor. If you get sensor-related alerts, use the tester to make sure that the transmitter is working.

Connecting the tester


- 1 Hold the transmitter and the tester as shown. Line up the flat side of the tester with the flat side of the transmitter.

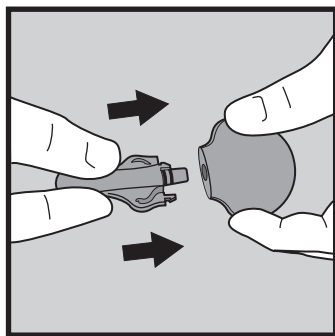

- 2 Push the tester into the transmitter until the flexible side arms of the tester click into the notches on both sides of the transmitter.

- 3 Within 20 seconds, the green light on the transmitter will flash for about 10 seconds when properly connected. The Sensor feature must be turned **On** on your pump. Go to the RECONNECT OLD SENSOR screen to start your sensor.


Main > Sensor > Link to Sensor > Reconnect Old Sensor

You should see the sensor icon  on the pump screen, which means that the transmitter and the sensor are communicating.

- 4 Within the next seven to eight minutes, go to the SENSOR STATUS screen to look for the Sensor ISIG value:
 - a. The sensor ISIG value on this screen should be between 24.00-29.00 nA. This range of ISIG values means that the transmitter electronics are working properly. Since the transmitter is sending signals correctly, the sensor must have caused the alert. Remove and discard the sensor. Insert a new sensor in a new site.
 - b. If you see a Sensor ISIG value that is less than 24.00 nA or more than 29.00 nA, contact your local representative. It may be time to replace your transmitter.

Disconnecting the tester

- 1 Hold the transmitter body as shown and pinch the side arms of the tester.

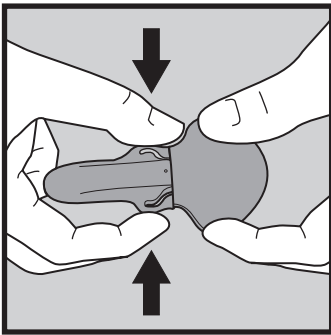

- 2 With the tester arms pinched, gently pull the transmitter away from the tester. To save transmitter battery life, do not leave the tester connected after testing.

Troubleshooting sensor features

Reconnect old sensor

You should only use this feature if you have disconnected the sensor from the transmitter and have to reconnect them.

To Reconnect old sensor:

- 1 Go to the RECONNECT OLD SENSOR screen.
Main > Sensor > Link to Sensor > Reconnect Old Sensor
- 2 Press **ACT** in the RECONNECT OLD SENSOR screen.
- 3 The SENSOR READY 2 HRS screen appears. Press any button to continue. An alert will sound in two hours when a meter BG is required for calibration.

Find lost sensor


If you receive a Lost Sensor alert:

- 1 Relocate the pump closer to the sensor, and go to the LINK SENSOR MENU.
Main > Sensor > Link to Sensor
- 2 Select **Find Lost Sensor**, then press **ACT**.
- 3 The SENSOR READY 15 MIN screen appears. Press any button to continue. An alert will sound in 15 minutes when a meter BG is required for calibration.

Understanding your transmitter, tester, and charger

The charger has a green light that shows the charging status and a red light that communicates any problems during charging. See the following graphic for the location of these lights:


Question: Why did the flashing green charger light turn off and a longer flashing red charger light turn on during charging?

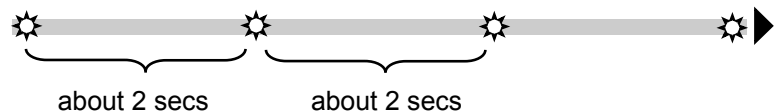

Answer: The transmitter battery is very low. Leave the transmitter on the charger for eight hours to completely recharge. If the red light is still flashing after eight hours, call your local representative. It may be time to replace your transmitter.


Question: Why do I see quick flashing red lights on the charger?

Answer: Your charger battery is low. Make sure that the transmitter is not connected to the charger. Replace the charger battery with a new AAA or LR-03 battery.

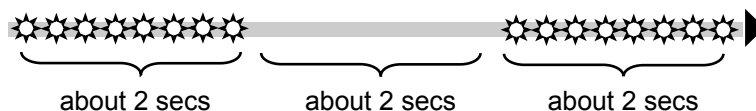

Question: Why do I see a mix of quick and long flashing red lights on the charger?

Answer: Your charger and transmitter batteries are very low. Replace the charger's AAA or LR-03 battery. If you now get the pattern for very low transmitter battery, leave the transmitter on the charger for eight hours to recharge. If the red light is still flashing after eight hours, call your local representative. It may be time to replace your transmitter.

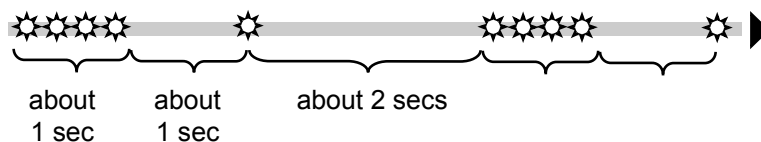

Question: I had my transmitter on the charger for a day. Will this damage my transmitter?

Answer: It will not damage the transmitter. You cannot overcharge it.

Question: What should I do if the green light on the transmitter does not flash when connected to the sensor?

Answer: Is the sensor inserted in your body? If it is not inserted, the transmitter will not flash green or send signals to a monitoring device.

If the sensor is inserted in your body, disconnect the transmitter from the sensor, wait several seconds and then reconnect it. If the green light still does not flash, charge the transmitter.

Question: Why did I not see the green light on the transmitter flash after connecting it to the tester?

Answer: Check the connection. If you still do not see a green light flash, fully recharge the transmitter battery. Test the transmitter with the tester. If you still do not see a green light flash, call your local representative. It may be time to replace your transmitter.

Maintenance

It is recommended that you read your warranty statement included in this user guide for information about what is covered during your warranty period.

Battery

The Paradigm pump uses a single 1.5 V AAA alkaline battery, size E92, type LR03. As a safety measure, Medtronic Diabetes has designed the pump to only accept a NEW battery. If you insert a used battery, the FAILED BATT TEST alarm may be triggered. Refer to the *Install battery* section in *The basics* chapter for instructions.

The use of cold batteries causes erratic pump behavior. To prevent this, do not use batteries that have been in cold storage (such as in the refrigerator or your car in the winter). It takes several hours for these batteries to warm to room temperature.

Certain features on the pump use a lot of battery power. Your battery will need to be replaced more frequently if you use these features:

- Remote control
- Meter
- Backlight
- Vibrate alert type setting

CAUTION: It is recommended that you use a AAA Energizer battery. Do not use a rechargeable or carbon zinc battery in the pump. Do not remove the battery unless you are changing your battery (installing a NEW battery). Replace it within five minutes. If not replaced within five minutes, the screen may display an alarm message. Follow the instructions in the message and make sure the time and date are set correctly. Check that the battery is inserted correctly. If the battery has been installed backwards, remove the battery and install it properly.

Storage

If you have to remove and store your pump, it is recommended that you store it with the battery in place. Keep a record of your current basal rates. To preserve battery life, reset the basal rates to 0 (zero), turn off the remote and meter options, and set the Auto-off to dashes or zeroes.

Cleaning your pump

- 1 Use only a damp cloth and mild detergent mixed with water to clean the outside of your pump.
- 2 Rinse the pump with a clean cloth, dampened with clean water.
- 3 Dry with a clean cloth.
- 4 Never use organic solvents, such as lighter fluid, nail polish remover, or paint thinner to clean your pump.
- 5 Keep the reservoir compartment and battery compartment dry and away from moisture.
- 6 Do not use any lubricants with your pump.
- 7 Use a 70 percent alcohol wipe to disinfect your pump.
- 8 Use a dry clean cotton tip to remove any battery residue from the battery cap.
- 9 Use a dry clean cloth to remove any battery residue from the battery compartment opening.

Cleaning the transmitter

CAUTION: The charger and the tester are NOT water-tight. Do NOT immerse in water.

Do NOT discard the transmitter in a medical waste container or otherwise subject it to incineration. Transmitter contains a battery which may explode upon incineration.

- 1 Wash your hands thoroughly.
- 2 Dampen a clean cloth with mild liquid soap and warm water. Wipe the outside of the transmitter.
- 3 Rinse the transmitter under warm tap water but do not get water inside the connector. If you get water inside the connector, shake the water out and allow it to air dry.
- 4 Using an antibacterial hand-sanitizer (available at a local drugstore) on a clean, dry cloth, wipe the transmitter's surface. Do NOT get any hand-sanitizer inside the connector. Repeated exposure to the hand-sanitizer could damage the connectors and affect the transmitter's performance as a result. If you get hand-sanitizer inside the connector, allow it to air dry.
- 5 Place the transmitter on a clean, dry cloth and air dry for 2-3 minutes.

Cleaning the Serter

This procedure describes how to clean the Enlite Serter (MMT-7510). If you are using a different sensor and insertion device, follow the directions in the user guide for the insertion device you are using.

- 1 Dampen a clean cloth with mild liquid soap and warm water. Wipe the Serter.
- 2 Rinse the Serter with warm tap water.
- 3 Wipe the Serter with 70% isopropyl alcohol.
- 4 Place the Serter on a clean, dry cloth and let air dry.

Pump specifications

This section provides detailed information on specifications related to your pump. The safety features of your pump are individually listed and described.

Alarms and error messages

- Indicators: audible tone (beep) or vibration (silent).
- All alarms and errors show messages on the pump's screen and provide instructions on what to do. Unresolved alarms will escalate to siren mode for added safety.

Alarm history

Maximum records shown: 36

Audio frequency

Name	Frequency (+/- 25%)
Alarm	1850 Hz
Escalated alarm	2.5 kHz, followed by 3.2 kHz
Alert	1.8 kHz, followed by 2.1 kHz and 1.8 kHz
Low Glucose Suspend alarm	2.5 kHz, followed by 3.2 kHz
High Glucose alert	2.5 kHz, followed by 2.79 kHz, 2.98 kHz and 3.2 kHz
Low Glucose alert	3.2 kHz, followed by 2.98 kHz, 2.79 kHz and 2.5 kHz
Predictive High Glucose alert	2.5 kHz, followed by 2.98 kHz and 3.2 kHz
Predictive Low Glucose alert	3.2 kHz, followed by 2.98 kHz and 2.5 kHz
Rising rate of change alert	2.5 kHz, followed by 3.2 kHz
Falling rate of change alert	3.2 kHz, followed by 2.5 kHz

Backlight

- LCD (Liquid Crystal Display) type
- Time-out: 30 seconds

Basal

- Delivery: 0.025-35 units/hour (maximum units: 35/hour)
- Factory maximum setting: 2.0 units/hour
- 3 patterns maximum, each with 48 rates
- Increments:
 - 0.025 units for basal amounts in the range of 0.025 to 0.975 units
 - 0.05 units for basal amounts in the range of 1 to 9.95 units
 - 0.1 units for basal amounts of 10.0 units or larger

BG Target

Maximum targets: 8

- Range: 3.3-13.9 mmol/L (60-250 mg/dL)
- Warning limits: less than 5.0 or greater than 7.8 mmol/L (less than 90 or greater than 140 mg/dL)

Bolus delivery

- Insulin delivered/stroke:
 - 0.025 units for bolus amounts in the range of 0.025 to 0.975 units
 - 0.05 units for bolus amounts larger than 0.975 units
- Fluid delivered/stroke: 0.25 µL (microliter) for 0.025 unit pump stroke; 0.5 µL for 0.05 unit pump stroke

Maximum bolus

- Range: 0.0-75.0 units
- Warning limit: greater than 25.0 units

Bolus history

Maximum records shown: 24

Bolus units

Increments: Based on selected scroll rate

Bolus Wizard feature

See the *Bolus Wizard feature specifications* section in this chapter for information.

Carb ratios

Maximum ratio settings	Range	Warning limits
8	1-200 grams/u	less than 5 or greater than 50 grams/u
	0.075-15.0 u/exch	less than 0.3 or greater than 3.0 u/exch

Carb units

The food entry when using Bolus Wizard feature:

- grams: 0-300 (increments: 1 gram)
- exchanges: 0.0-20 (increments: 0.5 exch)

Daily totals

Maximum records shown: 32 days of data, maximum display: 999.975 units/day. Daily delivery accuracy: 0.025 units.


Default screen

The HOME screen. For most screens, if no buttons are pressed for 30 seconds, the pump returns to this screen.

Delivery accuracy

- Delivery accuracy: +/-5 percent.
- All normal boluses will be delivered within a maximum of 28 minutes.
- Displacement (piston) accuracy: +/-2 percent.

The following is a representative delivery accuracy curve. All infusion sets have +/-5 percent overall delivery accuracy error.

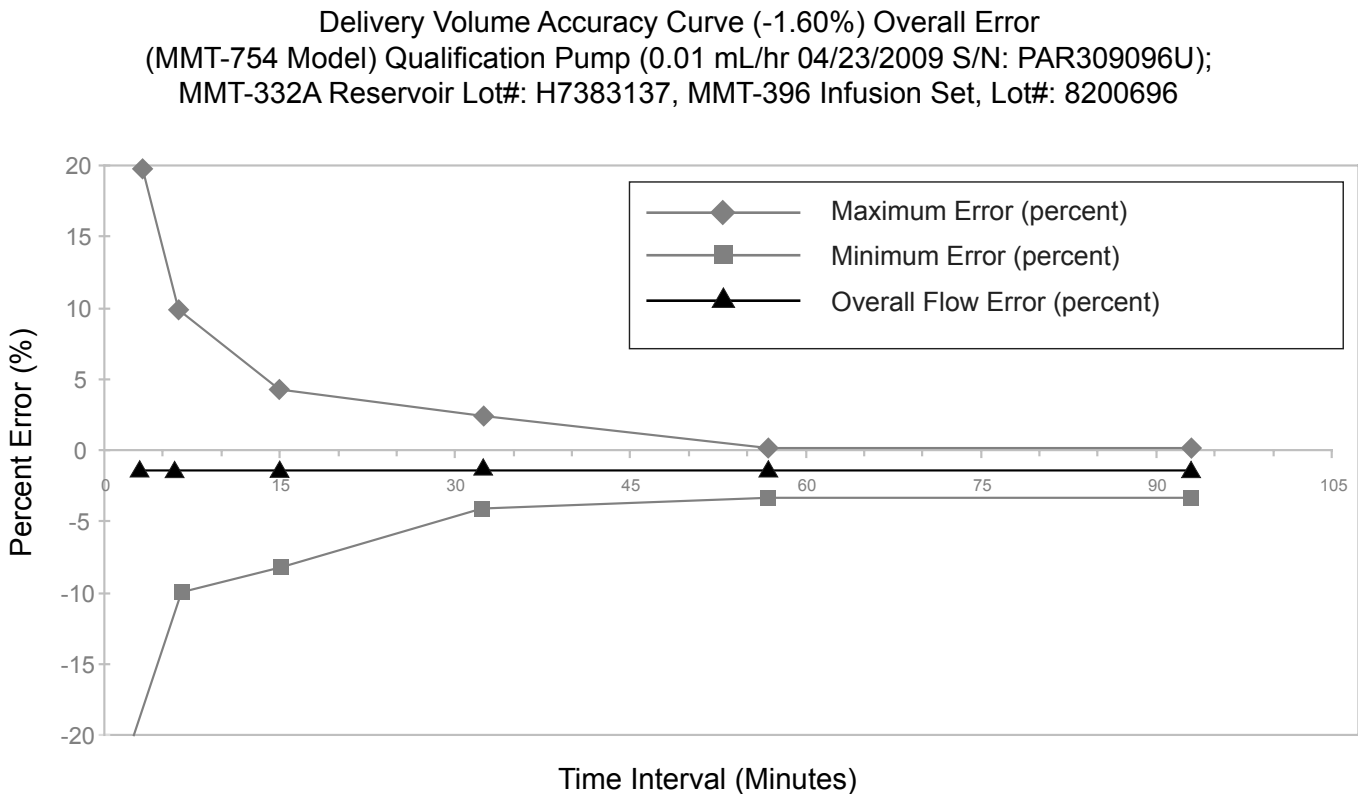

Pump motor

The pump motor has a unique, patented design with an integrated safety check system. The system delivers in precise increments.

Dual Wave bolus

Delivers a Normal bolus followed by a Square Wave bolus (limited by maximum bolus).

Easy bolus

- Features programming using audible tones (or vibrate pulses) in user-defined increments.
 - Beep mode range: 0 to maximum bolus
 - Vibrate mode range: 0 to 20 steps or maximum bolus, whichever comes first
- Factory default: off
- Default step increment: 0.1 unit
- Step size \leq maximum bolus.
- Adjustable step size: 0.1 to 2.0 units per step.
- Accessible from the remote control or pump buttons.

Filling the infusion set

- Filling the cannula: up to 10 units (limited by maximum bolus)
- Filling the tubing: warning at 30 units, then at each 10 units thereafter.
- Fill rate: 1 to 5 units/second.
- Insulin used to fill the infusion set is not counted in daily totals but is recorded separately in the Reservoir + Set history.

Infusion pressure

Maximum infusion pressure and occlusion pressure: 13.7 PSI (94.46 kPa).

Insulin sensitivity

- Maximum settings: 8
- Factory default: 2.8 mmol/L/u (50 mg/dL/u)
- Range: 0.5-22.2 mmol/L/u (10-400 mg/dL/u)
- Warning limits: less than 1.1 or greater than 5.6 mmol/L/u (less than 20 or greater than 100 mg/dL/u)

Low resv (reservoir) warning

Values are based on displayed amount, not actual amount.

time:	2-24-hours, and @ 1:00 hour before empty	08:00 hours (default when time is selected)
units:	5-50 units, and @ half of the remaining amount	20 units (factory default)

Meter value

The BG measurement received from the meter. Appears in the ENTER BG screen during bolus programming. Appears on the screen when the pump is idle at the HOME screen.

- Expiration: 12 minutes
- Range: 1.1-33.3 mmol/L (20-600 mg/dL)
- Maximum meter ID entries: 3

Normal bolus

Range 0.025 - 75.0 units of insulin (limited by maximum bolus setting).

Occlusion detection

When occlusion is detected, the NO DELIVERY alarm occurs. The occlusion alarm is triggered by an average of 2.77 units of missed insulin. The Paradigm pump is intended for use with U100 insulin. This table shows occlusion detection for three different situations when using U100 insulin.

Rate	Minimum time before alarm	Typical time before alarm	Maximum time before alarm
bolus delivery (10 units at 2 u/min)	23 seconds	49 seconds	67 seconds
basal delivery (1.0 u/h)	2.25 hours	3.13 hours	3.9 hours
basal delivery (0.025 u/h)	66 hours	122.73 hours	189 hours

Percent temp basal

Default value: 100% of basal programming

Power supply

The pump is powered by a standard 1.5 V AAA alkaline battery, size E92, type LR03 (Energizer brand recommended).

Reservoir + Set history

Maximum records shown: 20 (fill tubing and fill cannula)

Program safety checks

Maximum infusion with single fault condition: 0.0 units

Pump size

The dimensions of the pump are approximately:

- 554 pump:
 - 5.1 x 8.3 (7.1 at the battery cap) x 2.0 cm
 - 2.0 x 3.3 (2.8 at the battery cap) x 0.82 inches
- 754 pump:
 - 5.1 x 9.4 (8.9 at the battery cap) x 2.1 cm
 - 2.0 x 3.7 (3.5 at the battery cap) x 0.84 inches

Pump weight

- 554 pump: approximately 95 grams (with battery installed).
- 754 pump: approximately 102 grams (with battery installed).

Remote control

Uses radio signals to allow users to program Normal boluses or to suspend/resume their pumps.

Reservoir

- The user-filled reservoir is made from impact-resistant, insulin-compatible polypropylene.
- 554 pump volume: up to 176 units of U100 insulin
- 754 pump volume: up to 300 units of U100 insulin

Square Wave bolus

Delivers bolus insulin over a duration of 30 minutes up to eight hours (limited by the Max bolus setting).

Temporary (temp) basal rate

Allows you to temporarily change the current basal rate for a duration of 30 minutes up to 24 hours (limited by Maximum basal setting). The temp basal rate can be set to either Percent of basal or Insulin rate.

Time and date screen

12-hour or 24-hour formats. Pump users set the time/date, including the year, month and day. The date is included in the STATUS screen. Time always appears on the top of the screen.

Environmental conditions

- Pump operating temperature range: from 3° C (37° F) to 40° C (104° F)
- Air pressure range: from 700 hPa to 1060 hPa (from 10.2 psi to 15.4 psi)
- Operating humidity range of the case: 20% to 95%. This requirement exceeds EN 60601-2-24:1998, Section 10.2.1b (20% to 90%).

Status screen

Item	When	What
Active Insulin	(if any)	(Active Ins:) the amount of active insulin, displayed to the next 0.025 U
Alarm Clock*	(appears if scheduled to alarm)	Time set is displayed.
Alert Silence	(if active)	Hi, Lo, Hi/Lo, All with time left
Auto-off	(appears if enabled)	X HR
Basal pattern information	(if active)	Pattern A or B
Battery Status	(always appears)	Normal, Low, Off
BG meter value (most recent BG value received)	(appears if BG meter is enabled)	XX.X mmol/L (XXX mg/dL) time and date received

Item	When	What
BG Reminder* (Only if enabled)	Time remaining before BG Reminder is set to go off H:MM h (if less than 1 hour, 0:XXh where XX is minutes remaining)	
Block	(if active)	ON
Current date	(always appears)	
Current temp basal information	(if active)	rate (units per hour), duration, time remaining
Last alarm/alert	(if cleared within 24 hours)	(Last:) with the date, time and abbreviation for the most recently cleared alarm/alert
Last bolus information	(appears if at least one bolus was started)	type and units delivered delivery time and date (S = Square, N = Normal, DN = Normal part of dual wave bolus, DS = Square part of dual wave bolus)
Meter: Disabled	(appears if enabled but battery is low or empty)	
Meter: On	(appears if enabled)	
Missed Bolus Reminder	(if this feature is turned on)	Bolus Before and the time that the reminder is due
Pump model number	(always appears)	
Remote: Disabled	(appears if enabled but battery is low or empty)	
Remote: On	(appears if enabled)	
Reservoir started	(always appears)	date, time, units left, time left
Serial number	(always appears)	
Sensor Demo	(if Demo Mode is selected for the sensor graphs)	SENSOR DEMO: On
Software version	(always appears)	
Standard basal delivery data	(always appears)	current basal rate (basal 1, basal 2, etc.)
Status of pump	(such as Rewind, Suspended, Low Reservoir, Set Time, and so on)	

Item	When	What
Time	(always appears)	
*If all of the Reminders and the Alarm Clock are enabled, only the reminder that is nearest to alarming will appear in the STATUS screen.		

Bolus Wizard feature specifications

There are four different formulas the Bolus Wizard feature uses to estimate a bolus, depending on your current BG. The following formulas apply only when the carb units are in grams.

- 1 If your current BG is greater than your High BG Target, the Bolus Wizard feature subtracts active insulin from the BG correction estimate, then adds this to the food estimate to get the total bolus estimate. However, if the result of subtracting active insulin from BG correction estimate is a negative number (less than zero), the total bolus estimate is based only on the food estimate.

$$\begin{array}{ccccccc} & \text{(food estimate)} & & \text{(correction estimate)} & & & \\ \text{total bolus} & & & & & & \\ \text{estimate} = & \frac{A}{B} & + & \frac{C - D}{E} & - & \text{active insulin} & \end{array}$$

where: A = food (grams)
 B = carb ratio
 C = current BG
 D = High BG Target
 E = insullin sensitivity

- 2 If your current BG is less than your Low BG Target, the Bolus Wizard feature adds the BG correction estimate to the food estimate to get the total bolus estimate.

$$\begin{array}{ccccccc} & \text{(food estimate)} & & \text{(correction estimate)} & & & \\ \text{total bolus} & & & & & & \\ \text{estimate} = & \frac{A}{B} & + & \frac{C - D}{E} & & & \end{array}$$

where: A = food (grams)
 B = carb ratio
 C = current BG
 D = Low BG Target
 E = insullin sensitivity

- 3 If your current BG is between or equal to your High or Low BG Target, the total bolus estimate is based only on the food estimate.

(food estimate)

$$\text{total bolus estimate} = \frac{\text{food (grams)}}{\text{carb ratio}}$$


- 4 If you do not enter a BG, the total bolus estimate is based only on the food estimate.

NOTES:

- *If a Dual Wave bolus is less than the estimate due to the max bolus limit or a user change, the square (sq) portion is reduced first.*
- *Based on the Active Insulin Time setting you choose, the pump keeps track of how much insulin is still active in your body. This is shown as Active Ins. or Act. Ins. in the STATUS, SET BOLUS, SET NORMAL BOLUS, SET SQUARE BOLUS, SET DUAL BOLUS TOTAL, and SET EASY BOLUS screens. This prevents stacking of insulin, and lowers the chances of hypoglycemia.*
- *The active insulin amount is also shown in the Bolus Wizard's ESTIMATE DETAILS screen but is calculated differently than in the above screens. The active insulin amount calculated in the ESTIMATE DETAILS screen includes the insulin that has already been delivered and insulin that is going to be delivered by the active Square Bolus. To show this difference in calculation, the active insulin appears with an asterisk (*Active Insulin) in this screen.*



- The Bolus Wizard feature may utilize your current BG measurement, carbohydrate consumption, and active insulin to calculate your estimated bolus.
- Active Insulin Curves

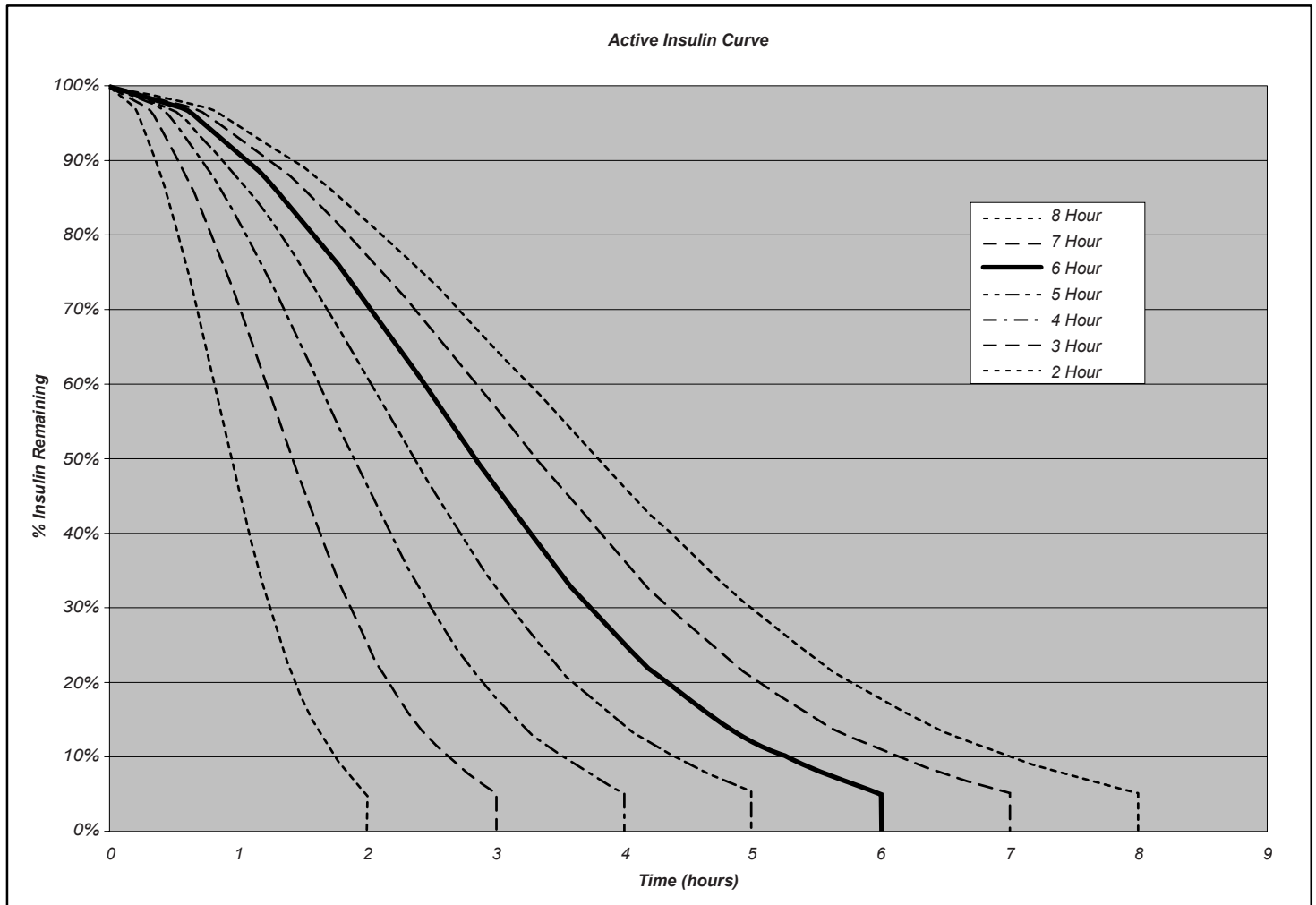

Graph adapted from Mudaliar and colleagues, *Diabetes Care*, Volume 22, Number 9, Sept. 1999, page 1501.

Bolus Wizard feature examples

Settings:

- Carb ratio: 30 grams/unit
- BG Target: 5.0-6.7 mmol/L (90-120 mg/dL)
- Insulin Sensitivity: 2.2 mmol/L/units (40 mg/dL/unit)
- Active Insulin Time: 6 Hours

1 No active insulin from previous bolus delivery. User inputs 60 grams of carbs and does not enter BG.

$$\text{Food estimate:} \quad \frac{60 \text{ g}}{30 \text{ g/u}} = 2 \text{ units}$$

$$\text{Bolus estimate:} \quad 2 \text{ units}$$

2 No active insulin from previous bolus delivery. User inputs 60 grams of carbs and a BG of 11.1 mmol/L (200 mg/dL).

$$\text{Food estimate:} \quad \frac{60 \text{ g}}{30 \text{ g/u}} = 2 \text{ units}$$

$$\text{Correction estimate:} \quad \frac{11.1 \text{ mmol/L} - 6.7 \text{ mmol/L}}{2.2 \text{ mmol/L/u}} = 2 \text{ units}$$

or

$$\frac{200 \text{ mg/dL} - 120 \text{ mg/dL}}{40 \text{ mg/dL/u}} = 2 \text{ units}$$

$$\text{Bolus estimate} = \text{Food estimate} + \text{Correction estimate}$$

$$\text{Bolus estimate} = 2 \text{ units} + 2 \text{ units} = 4 \text{ units}$$

- 3** No active insulin from previous bolus delivery. User inputs 60 grams of carbs and a BG of 3.9 mmol/L (70 mg/dL).

Food estimate: $\frac{60 \text{ g}}{30 \text{ g/u}} = 2 \text{ units}$

Correction estimate: $\frac{3.9 \text{ mmol/L} - 5.0 \text{ mmol/L}}{2.2 \text{ mmol/L/u}} = -0.5 \text{ units}$
or

$\frac{70 \text{ mg/dL} - 90 \text{ mg/dL}}{40 \text{ mg/dL/u}} = -0.5 \text{ units}$

Bolus estimate = Food estimate + Correction estimate

Bolus estimate = 2 units + (-0.5) units = 1.5 units

- 4** No active insulin from previous bolus delivery. User inputs 60 grams of carbs and a BG of 5.6 mmol/L (100 mg/dL).

Food estimate: $\frac{60 \text{ g}}{30 \text{ g/u}} = 2 \text{ units}$

Correction estimate: Correction is 0 because the current blood glucose reading is between the Low and High BG Targets.

Bolus estimate = Food estimate + Correction estimate

Bolus estimate = 2 units + 0 = 2 units

- 5** Previous bolus activity results in a calculation of 1.5 units unabsorbed (active) insulin. User inputs 60 grams of carbs and a BG of 11.1 mmol/L (200 mg/dL).

$$\text{Food estimate: } \frac{60 \text{ g}}{30 \text{ g/u}} = 2 \text{ units}$$

Active insulin: 1.5 units

$$\text{Correction estimate: } \frac{11.1 \text{ mmol/L} - 6.7 \text{ mmol/L}}{2.2 \text{ mmol/L/u}} - \text{Active insulin} = 0.5 \text{ units}$$

or

$$\frac{200 \text{ mg/dL} - 120 \text{ mg/dL}}{40 \text{ mg/dL/u}} - \text{Active insulin} = 0.5 \text{ units}$$

Bolus estimate = Food estimate + Correction estimate

$$\text{Bolus estimate} = 2 \text{ units} + 0.5 \text{ units} = 2.5 \text{ units}$$

- 6** Previous bolus activity results in a calculation of 3.5 units unabsorbed (active) insulin. User inputs 60 grams of carbs and a BG of 11.1 mmol/L (200 mg/dL).

$$\text{Food estimate: } \frac{60 \text{ g}}{30 \text{ g/u}} = 2 \text{ units}$$

Active insulin: 3.5 units

$$\text{Correction estimate: } \frac{11.1 \text{ mmol/L} - 6.7 \text{ mmol/L}}{2.2 \text{ mmol/L/u}} - \text{Active insulin} = -1.5 \text{ units}^*$$

or

$$\frac{200 \text{ mg/dL} - 120 \text{ mg/dL}}{40 \text{ mg/dL/u}} - \text{Active insulin} = -1.5 \text{ units}^*$$

*This negative number indicates that active insulin is sufficient to cover the correction that is needed. Thus, correction is 0 units. Active insulin does not reduce the food portion of the estimate.

Bolus estimate = Food estimate + Correction estimate

$$\text{Bolus estimate} = 2 \text{ units} + 0 = 2 \text{ units}$$

Insulin delivery default settings

Menu	Item	Default Setting	Limits	Increments
Bolus Menu:	Bolus Wizard feature:	Off		
	Easy bolus:	Off		
	Easy bolus step:	0.1 u	2.0 u	0.1 u
	Scroll Rate:	0.10 u		0.025 u; 0.05 u; 0.10 u
	Dual/Square bolus:	Off		
	Maximum bolus:	10.0 u	0.0-75.0 u (per single bolus)	
	BG Reminder:	Off	0:00 - 5:00	0:30
Basal Menu:	Patterns:	Off		
	Maximum basal rate:	2.00 u/h	0.000-35.0 u/h	0.025 u (for 0.025-0.975 u/h); 0.05 u (for 1.00-9.95 u/h); 0.1 u (for rates of 10 u/h or more)
	Basal rate:	0.000 u/h		0.025 u (for 0.025-0.975 u/h); 0.05 u (for 1.00-9.95 u/h); 0.1 u (for rates of 10 u/h or more)
	Temp basal type:	U/H	Max basal rate setting	0.025 u (for 0.025 - 0.975 u/h); 0.05 u (for 1.00 - 9.95 u/h); 0.1 u (for rates of 10 u/h or more); (or 1%)
Utilities Menu:	Lock Keypad:	Off		
	(Alarm) History:	(no defaults)		
	Alert type:	audio, beep-med		
	Auto-off:	Off		

Menu	Item	Default Setting	Limits	Increments
	Capture Option:	Off		
	Low reservoir warning:	(20) insulin units	if units: 1st @ 5-50 u; 2nd @ half amount remaining (if time: 1st @ 2:00-24:00; 2nd @ 1:00 hour before empty)	20 u (0:30)
	(Time/Date) Time:	12 a.m. (midnight)		
	(Time/Date) Date:	1/1/07		
	(Time/Date) Time format:	12-hour		
	Block:	Off		
	Alarm clock:	Off		
	Remote option:	Off		
	Meter option:	Off		
	User Settings:	(no defaults)		
	Language:	English		

Bolus Wizard feature default settings

Item	Default Setting	Limits	Increments	Warning Limits
carb units:	grams	—	—	—
ins to carb (or exch) ratio:	15 grams/u or 1 unit/exch	1-200 g/u or 0.075- 15.0 u/exch	0.1 g/u for 0-9.9 g/u; 1 g/u for ratios of 10.0 g/u or more (or 0.001 u/exch for 0.075- 0.099 u/exch; 0.01 u/exch for 0.10-9.99 u/exch; 0.1 u/exch for 10.0- 15.0 u/exch)	5-50 g/u or 0.3- 3.0 u/exch
(insulin) sensitivity:	2.8 mmol/L/ u or 50 mg/dL/u	0.5-22.2 mmol/L/u or 10-400 mg/dL/u	0.1 mmol/L/u or 1 mg/dL/ u	1.1-5.6 mmol/L/ u or 20- 100 mg/dL/u
BG Target:	5.6- 5.6 mmol/L or 100- 100 mg/dL	3.3-13.9 mmol/L or 60-250 mg/dL	0.1 mmol/L or 1 mg/dL	5.0-7.8 mmol/L or 90-140 mg/dL
Active Ins Time:	6 hours	2-8 hours	1 hour	—

Sensor features default settings

Feature	Options	Default Setting	Limits, Range of Values
Sensor	<ul style="list-style-type: none"> On Off 	Off	—
Glucose Alerts	<ul style="list-style-type: none"> On Off 	Off	Set the limits using the Glucose Limits feature.
Glucose Limits	<ul style="list-style-type: none"> Maximum of 8 Glucose Limits. One pair of Glucose Limits (Glucose Limits 1) is already set up for you, with a default High Glucose Limit and Low Glucose Limit. Glucose Limits 2 through 8 can be set as needed. 		

Feature	Options	Default Setting	Limits, Range of Values
Glucose Limits 1	<ul style="list-style-type: none"> Low: 4.4 mmol/L (80 mg/dL) High: 13.3 mmol/L (240 mg/dL) 	Low: On 4.4 mmol/L (80 mg/dL) High: On 13.3 mmol/L (240 mg/dL)	Low Glucose Limit range: ^a 2.2 - 21.6 mmol/L (40 - 390 mg/dL) High Glucose Limit range: ^b 2.8 - 22.2 mmol/L (50 - 400 mg/dL)
Glucose Limits 2 through Glucose Limits 8	<ul style="list-style-type: none"> Low: Off, or a valid range from 2.2 to 21.6 mmol/L (40 to 390 mg/dL) High: Off, or a valid range from 2.8 to 22.2 mmol/L (50 to 400 mg/dL) 	Low Off High Off	See the Glucose Limits 1 range.
High Repeat	<ul style="list-style-type: none"> On Off (Set to On when High Glucose, High Predictive, or Rising Rate of Change Alert is on)	On 1 hour	Range: 0:05-3:00 (5 minutes to 3 hours)
Low Repeat	<ul style="list-style-type: none"> On Off (Set to On when Low Glucose, Low Predictive, or Falling Rate of Change Alert is on.)	On 20 minutes	Range: 0:05 - 1:00 (5 minutes to 1 hour)
Predictive Alert	<ul style="list-style-type: none"> On Off 	Off Default time sensitivity: Low 15 minutes High 15 minutes	Range: Off 0:05 - 0:30 (5 to 30 minutes)

Feature	Options	Default Setting	Limits, Range of Values
Rate of Change Alerts	<ul style="list-style-type: none"> On Off 	Fall Rate Limit Off Rise Rate Limit Off	Fall Rate Limit range: 0.065-0.275 mmol/L/min (1.1-5.0 mg/dL/min) Rise Rate Limit range: 0.065-0.275 mmol/L/min (1.1-5.0 mg/dL/min)
Low Glucose Suspend	<ul style="list-style-type: none"> On Off 	Off	Range: 3.3-6.1 mmol/L (60-110 mg/dL)
Cal Repeat	(Always on - cannot be turned off)	On 30 minutes	Range: 0:05 - 1:00 (5 minutes to 1 hour)
Cal Reminder	<ul style="list-style-type: none"> On Off 	On 1 hour	Range: 0:05 - 6:00 (5 minutes to 6 hours)
Auto Calibration	<ul style="list-style-type: none"> On Off 	Off	—
Blood Glucose Units (BG Units)	<ul style="list-style-type: none"> mg/dL mmol/L 	mg/dL	—
Weak Signal	(Always on - cannot be turned off)	30 minutes	Range: 0:05 - 0:40 (5 to 40 minutes)
Sensor glucose graph timeout	<ul style="list-style-type: none"> 2 minutes 4 minutes 6 minutes NONE^c 	2 minutes	—

Feature	Options	Default Setting	Limits, Range of Values
Alert Silence	<ul style="list-style-type: none"> Off Hi Alerts Lo Alerts Hi and Lo Alerts All Sensor Alerts 	Off	If any option other than Off is selected in the ALERT SILENCE screen, the SET DURATION screen displays. The value range is from 00:30 to 24:00.
AUC Limits	Always on	Low Limit 3.9 mmol/L (70 mg/dL) High Limit 10 mmol/L (180 mg/dL)	Low Limit range: ^d 2.2-22.2 mmol/L (40-400 mg/dL) High Limit range: ^e 2.2-22.2 mmol/L (40-400 mg/dL)

- a. Your Low Glucose Limit must be at least 0.6 mmol/L (10 mg/dL) below your High Glucose Limit. For example, if your High Glucose limit is 11.1 mmol/L (200 mg/dL), then your Low Glucose limit cannot be higher than 10.5 mmol/L (190 mg/dL).
- b. Your High Glucose Limit must be at least 0.6 mmol/L (10 mg/dL) above your Low Glucose Limit.
- c. If you select this option the sensor glucose graphs will show continuously on the pump screen.
- d. The AUC Low Limit can be equal to, but not higher than the AUC High Limit.
- e. The AUC High Limit can be equal to, but not lower than the AUC Low Limit.

Guidance and manufacturer's declaration


Guidance and Manufacturer's Declaration - Electromagnetic Emissions		
The Paradigm 554/754 pump system (consisting of the MMT-554/754 pump and MMT-7703 transmitter) is intended for use in the electromagnetic environment specified below. The customer or the user of the Paradigm 554/754 pump system should make sure that it is used in such an environment.		
Emissions Test	Compliance	Electromagnetic Environment - Guidance
RF emissions CISPR 11	Group 1	The Paradigm 554/754 pump system uses RF energy only for system communication functions. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11 Harmonic emissions IEC 61000-3-2 Voltage fluctuations/flicker emissions IEC 61000-3-3	Class B Not applicable Not applicable	The Paradigm 554/754 pump system is suitable for use in all establishments other than domestic and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes. NOTE: The preceding statement is required by IEC 60601-1-2 for Group 1, Class B devices. However, since the Paradigm 554/754 pump system is battery powered, its emissions will not be affected by the establishment power supply and there is no evidence of any issues associated with the use of the system in domestic establishments.

Guidance and Manufacturer's Declaration - Electromagnetic Immunity			
<p>The Paradigm 554/754 pump system (consisting of the MMT-554/754 pump and MMT-7703 transmitter) is intended for use in the electromagnetic environment specified below. The customer or the user of the Paradigm 554/754 pump system should make sure that it is used in such an environment.</p>			
Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment - Guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±6 kV contact ±8 kV air	±30 kV air (<5% relative humidity)	The Paradigm 554/754 pump system should not be affected by electrostatic discharge that might occur under normal conditions of use.
Electrical fast transient/burst IEC 61000-4-4	±2 kV for power supply lines ±1 kV for input/output lines	Not applicable Not applicable	Not applicable Not applicable
Surge IEC 61000-4-5	±1 kV line(s) to line(s) ±2 kV line(s) to earth	Not applicable Not applicable	Not applicable Not applicable
Voltage dips, short interruptions and voltage variations on power supply lines IEC 61000-4-11	<5% U_T (>95% dip in U_T) for 0.5 cycle 40% U_T (60% dip in U_T) for 5 cycles 70% U_T (30% dip in U_T) for 25 cycles <5% U_T (>95% dip in U_T) for 5 seconds	Not applicable Not applicable Not applicable Not applicable	Not applicable Not applicable Not applicable Not applicable
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
<p>NOTE: U_T is the a.c. mains voltage prior to application of the test level.</p>			

Guidance and Manufacturer's Declaration - Electromagnetic Immunity

The Paradigm 554/754 pump is intended for use in the electromagnetic environment specified below. The customer or user of the Paradigm 554/754 pump system should make sure that it is used in such an environment.

Immunity Test	IEC 60601 Level	Compliance Level	Electromagnetic Environment Guidance
Conducted RF IEC 61000-4-6	10 Vrms 150 kHz to 80 MHz	Not applicable	<p>Portable and mobile RF communications equipment should be used no closer to any part of the Paradigm 554/754 pump, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.</p> <p>Recommended separation distance</p> <p>Not applicable</p>

Guidance and Manufacturer's Declaration - Electromagnetic Immunity			
The Paradigm 554/754 pump is intended for use in the electromagnetic environment specified below. The customer or user of the Paradigm 554/754 pump system should make sure that it is used in such an environment.			
Immunity Test	IEC 60601 Level	Compliance Level	Electromagnetic Environment Guidance
Radiated RF IEC 61000-4-3	10 V/m 80 MHz to 6.0 GHz	10 V/m	<p>$d = 0.35 \sqrt{P}$ 80 MHz to 800 MHz</p> <p>$d = 0.70 \sqrt{P}$ 800 MHz to 6.0 GHz</p> <p>Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m).</p> <p>Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey^a, should be less than the compliance level in each frequency range^b.</p> <p>Interference may occur in the vicinity of equipment marked with the following symbol:</p> 

Guidance and Manufacturer's Declaration - Electromagnetic Immunity

The Paradigm 554/754 pump is intended for use in the electromagnetic environment specified below. The customer or user of the Paradigm 554/754 pump system should make sure that it is used in such an environment.

Immunity Test	IEC 60601 Level	Compliance Level	Electromagnetic Environment Guidance
<p>NOTE: At 80 MHz and 800 MHz, the higher frequency range applies.</p> <p>NOTE: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption, and reflection from structures, objects and people.</p> <p>^aField strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcasts and TV broadcast cannot be predicted theoretically with accuracy. To access the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the Paradigm 554/754 pump is used exceeds the applicable RF compliance level above, the Paradigm 554/754 pump should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the Paradigm 554/754 pump.</p> <p>^bOver the frequency range 150 kHz to 80 MHz, field strengths should be less than $[V_1]$ V/m.</p>			

Recommended separation distances between portable and mobile RF communications equipment and the Paradigm 554/754 pump system

This section provides information on the recommended separation distance between portable and mobile RF communications equipment and the Paradigm 554/754 pump. The Paradigm 554/754 pump is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. Paradigm 554/754 pump users can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the Paradigm 554/754 pump as recommended below, according to the maximum output power of the communications equipment.













	Separation distance according to the frequency of transmitter (m)	
Rated maximum output power of transmitter (W)	80MHz to 800MHz $d = 0.35 \sqrt{P}$	800MHz to 6.0GHz $d = 0.70 \sqrt{P}$
0.01	0.035	0.07
0.1	0.11	0.11
1	0.35	0.7
10	1.1	2.2
100	3.5	7



For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where p is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

Icon table

Follow instructions for use	
Attention: Read all warnings and precautions in instructions for use.	
Manufacturer	
Date of manufacture (year - month)	
Batch code	LOT
Use by (year - month)	
Catalogue number	REF
Device serial number	SN
Storage temperature range	
Magnetic Resonance (MR) unsafe	
Fragile product	
Type BF equipment (Protection from electrical shock)	
Conforms to IEC60601-1 sub-clause 44.6 and IEC60529 standard	IPX7
Recycle	
Radio communication device	
Signifies compliance to Australian EMC and Radio communications requirements	

Signifies compliance with Industry Canada EMC and Radio communications requirements	IC
European conformity. This symbol means the device fully complies with MDD93/42/EEC (NB 0459).	
European conformity. This symbol means the device complies with R&TTE1999/5/EC directives (NB 0976). Applies to 868.35 MHz device (MMT-554WW, MMT-754WW).	

Glossary

A

A1C (HbA1c) test - A 2-3 month average of blood glucose values expressed in percent. The normal range varies with different labs and is expressed in percent (such as 4 - 6%).

Accept - Pressing the **ACT** button to approve the selection or setting.

Active insulin - Bolus insulin that has been delivered to your body, but has not yet been used.

Adverse reaction - An unexpected, unpleasant or dangerous reaction to a sensor when it is inserted into the body. An adverse reaction may be sudden or may develop over time.

Alarm - Audible or vibrating (silent) notice that indicates the pump is in Attention mode and immediate attention is required. Alarms are prefixed in the alarm history with the letter A.

Alarm clock - Feature you can set to go off at specified times of the day.

ALARM HISTORY - Screen that displays the last 36 alarms/errors that have occurred on your pump.

Alarm icon - a solid circle that shows at the top of the screen and the pump beeps or vibrates periodically until the condition is cleared (see Attention mode).

Alert - Audible or vibrating (silent) indicator that notifies you the pump needs attention soon or that you should be reminded of something. Insulin delivery continues as programmed.



Alert icon - an open circle that appears at the top of the screen and the pump beeps or vibrates periodically to remind you of the condition (see Special mode).

Attention mode - Operating mode that stops all current insulin delivery. This mode indicates an alarm or condition

exists that requires immediate attention.

AUC Limits - The AUC Limits feature allows you to see how long your sensor glucose measurements are outside your target range. Your target range is the area between your Low AUC Limit and your High AUC Limit. When you set the AUC Limits, you select the Low AUC Limit and High AUC Limit you want for your target range.

B

Backlight - Pump screen light. Turns on whenever the  button is pressed or in combination with the  button if not on the HOME screen. The backlight also appears during alarms and alerts, except for LOW BATTERY alert and OFF NO POWER alarm.

Basal patterns - The user can program up to three different basal rate patterns into the pump for the delivery of basal

insulin: Standard, A, and B. For each pattern, there is an option of setting up to 48 basal rates.

Basal rate - The pump setting that provides a continuous infusion of insulin to keep the blood glucose stable between meals and during the night. Basal insulin mimics pancreatic insulin delivery—which meets all the body's non-food related insulin needs.

BASAL REVIEW screen - shows the basal rates programmed in the pump, with the 24-hour total for each rate.

BG - Blood Glucose

BG Reminder - Feature that you can set to remind you to check your blood glucose after a bolus.

BG Target - Normal blood glucose level

BG Units - Blood glucose units used by the pump (mg/dL or mmol/L).

Block - Feature that restricts access to all programming except suspend, selftest and the delivery of a bolus with the remote control.

Blood glucose (BG) - A form of digested sugar in the bloodstream. Glucose is the major source of energy for living cells, and is carried to each cell through the blood stream. Cells cannot use

glucose without the help of insulin.

Bolus - A dose of insulin given to cover an expected rise in blood glucose (such as the rise after a meal) or to lower a high blood glucose down to target range.

BOLUS HISTORY - This screen displays the last 24 boluses delivered by your pump.

Bolus Wizard® feature - Calculates the bolus amount based on personal information of the pump user.

C

Cal Reminder - The value of the Cal Reminder is the amount of time before the current calibration value expires when the user wants to be reminded to calibrate by having the pump issue a METER BG BY alert. For example, if the Cal Reminder is set to two hours, the METER BG BY alert will occur two hours before the calibration is required.

Cal Repeat - Once a METER BG NOW alert occurs, the pump will not repeat the alert until after this selected Cal Repeat time.

Calibrate - Check, adjust, or set to a standard (calibrate your pump).

Cannula - A short, thin, and flexible tubing at the end of the infusion set that is inserted into the subcutaneous tissue to deliver insulin.

Carb ratio (Carbohydrate ratio) - Used when counting carbohydrates in grams. The amount of carbohydrates covered by one unit of insulin. (Also see exch ratio.)

Carb units - The food entry when using the Bolus Wizard feature. Entered as (carbohydrate) grams or exchanges.

Carbohydrate - One of the three main sources of energy or calories in food. Carbohydrates are mainly sugars found in fruits, milk, and starches that the body breaks down into glucose.

CH - Carbohydrate

Contraindication - A condition that makes a specific treatment or procedure NOT advisable. A contraindication is literally contra- (against) an indication, against something that is indicated as advisable or necessary.

Correction bolus - The amount of insulin needed to return a high blood glucose level back down to target range.

D

Daily totals - Shows the total insulin delivered (basal and

bolus) in the last 24-hours.
Maximum records: 32 days

Default - A setting or value for your pump that is assigned automatically by your system. Some default settings cannot be changed; other default settings remain in effect until you change them.


Diabetic Ketoacidosis (DKA) - A serious condition that occurs when the insulin levels are low, blood glucose level are elevated, and the body uses fat for energy. This process produces ketones which upset the body's acid-base balance leading to a potentially life threatening situation.

DKA - Diabetic Ketoacidosis

Dual Wave® bolus - Combination of a Normal bolus that is delivered immediately, then followed by a Square Wave bolus. The Square Wave portion is delivered evenly over a period of time.

Duration - Amount of time it takes to administer a bolus or basal delivery. Also, length of time for an action or condition.

E

Easy bolus™ - Method of delivery for a Normal bolus using the Easy Bolus button .


Electromagnetic compatibility - The condition that exists when

systems and devices that use electromagnetic energy operate properly without causing or suffering accidental electromagnetic interference to or from other electromagnetic devices.

Electrostatic discharge - The rapid, spontaneous transfer of electrostatic charge induced by a high electrostatic field. The charge usually flows through a spark, known as static discharge, between two objects at different electrostatic levels as they approach each other (for example, when people touch each other after dragging their feet on a carpet).

Exch ratio - (Exchange ratio)
Used when counting carbohydrates as exchanges. The amount of insulin required to cover one (1) carbohydrate exchange. (Also see carb ratio.)

Exchange system - A system of estimating carbohydrate and calorie intake that groups by type. Each exchange group lists similar types of food in serving sizes. A person can exchange, trade, or substitute a food serving in a particular group for another food serving in the same group.

Express bolus - Method of delivery for any bolus type using the **EXPRESS BOLUS** button .

F

Fats - One of the three main classes of foods and a source of calories or energy in foods. Fats help the body use some vitamins and keep the skin healthy.

Fiber - A substance found in foods that comes from plants. Fiber may help in the digestive process and is thought to lower cholesterol and help control blood glucose. Fiber grams are not absorbed, and can be subtracted from the Total Carbohydrate amount for foods or meals with more than 5 grams of fiber.

Fingerstick - One method used to test blood glucose by taking a sample of blood for testing from a finger using a lancet or automatic finger puncture device.

Food bolus - A dose of insulin given to cover the expected rise in blood glucose that occurs after eating.

G

Gastroparesis - A condition of the digestive system that slows down the emptying of food from the stomach.

Generate - create a result (your pump will generate a graph).

Glucagon - A hormone that raises the level of glucose in the blood. The alpha cells of the pancreas make glucagon when the body needs to put more sugar into the blood. Glucagon primarily causes sugar stored in the liver to be released into the blood stream.

Glucagon Emergency Kit - An injectable form of glucagon that can be used to treat severe low blood glucose levels (hypoglycemia). A Glucagon Emergency Kit can be obtained in a pharmacy with a prescription.

Glucose tablet - A simple sugar in tablet form, that is taken orally to treat low blood glucose (hypoglycemia). Check with your local suppliers to purchase glucose tablets.

H

HbA1c - Glycosylated hemoglobin

HDL - High-density lipoprotein A complex of lipids and proteins in approximately equal amounts that functions as a transporter of cholesterol in the blood.

High Glucose - The pump displays an alert if the sensor indicates that your sensor glucose is at or above this value. You have the option to turn this feature on or off.

High Repeat - For the High Glucose, High Predictive, and Rising Rate of Change alerts, allows the user to set the delay between the first alert and any subsequent alerts. This will allow the user to avoid getting an alert every five minutes until the condition is corrected.

Hold - Press and continue to press a pump button.

Hyperglycemia - A condition in which blood glucose values are elevated above target. Symptoms may include excessive thirst, frequent urination, dry mouth, headache, fatigue, and blurred vision.

Hypoglycemia - A condition in which blood glucose values are low below 3.9 mmol/L (70 mg/dL). Symptoms may include behavioral changes, pale complexion, hunger, sweating, sudden weakness, headache, confusion, drowsiness, non-responsiveness to questions, and the most serious: seizure, convulsion or loss of consciousness.

I

Idle - The pump is at the HOME screen.

Increment - A small increase or decrease in a measurement.

Indication - A condition that makes a specific medical treatment or procedure advisable or necessary.

Infusion set - Flexible tubing with a reservoir connector and an infusion site. This tubing delivers insulin from the pump to the body.

Infusion site - The end of the infusion set held to the body with a tape. It consists of a cannula and an introducer needle.

Initialize - To set a starting value and/or prepare your pump, transmitter or sensor for use.

Insulin - A hormone that helps the body use glucose (sugar) for energy. The beta cells of the pancreas make the insulin.

Insulin sensitivity - The amount by which your blood glucose (BG) level is reduced by one unit of insulin. (Bolus Wizard feature data)

Insulin to carbohydrate ratio - The amount of insulin required to cover a given number of carbohydrates. This ratio helps you estimate what size bolus you should take when eating carbohydrates. Your healthcare professional will determine your insulin to carbohydrate ratio.

Interstitial fluid - The fluid between the cells in the human body.

Introducer needle - This needle allows the insertion of a cannula or a sensor into the subcutaneous tissue. It is removed and discarded after insertion leaving only the cannula or the sensor in the body.

K

Ketone - A chemical created by the human body when there is not enough insulin in the blood.

kPa (kilopascal) - A unit of measurement for quantifying force. Used to measure atmospheric pressure. Equivalent to 10,000 dynes per square centimeter.

L

LDL (Low-density lipoprotein) - A complex of lipids and proteins, with greater amounts of lipid than protein, that transports cholesterol in the blood.

Link - To turn on and setup the meter option that allows the pump to receive BG readings from a meter that communicates with your pump.

Log - A record of measurements, alarms or other actions by your pump.

Low Glucose - The pump displays an alert if the sensor indicates that your sensor glucose is at or below this value. You have the option to turn this feature on or off.

Low Glucose Suspend - This feature is set up to stop insulin delivery when your sensor glucose reaches or drops below your Low Suspend limit.

Low Repeat - For the Low Glucose, Low Predictive, and Falling Rate of Change alerts, allows you to set the delay between the first alert and any subsequent alerts. This will allow the user to avoid getting an alert every five minutes until the condition is corrected.

Low resv warning - Programmable warning that notifies you with an alert when either a specified number of units remain in the reservoir or a specified amount of time remains before the reservoir will be empty.

M

Manual bolus - Selectable item available in the BOLUS MENU when the Bolus Wizard feature is active. It is one method of programming a bolus without the Bolus Wizard feature. (see *Set bolus*)

Marker - A measurement of a specific event - such as a BG measurements, exercise, or the amount of insulin used.

Max basal rate - The maximum amount of basal insulin that the pump will deliver at one time. (set by the user)

Max bolus - The maximum amount of bolus insulin that the pump will deliver at one time. (set by the user)

Meter - An optional blood glucose meter powered by MWT1 technology. Your pump can be programmed to receive your BG readings from this meter.

Meter option - Feature that allows the pump to receive BG readings from a meter powered by MWT1 technology.

Mineral - Minerals and vitamins are essential micronutrients, needed in very small amounts to keep our bodies functioning properly.

MWT1 technology - MWT1 is the wireless Radio Frequency (RF) technology that is used to transmit information from the meter to the pump. You can program your pump to automatically receive your BG reading from this meter.

N

Normal bolus - An immediate delivery of a specified unit amount of insulin.

Normal mode - Regular operating mode. No special features are active, no alert or alarm conditions exist. Insulin delivery is normal during this mode.

Now portion - The normal portion of a Dual Wave bolus. The now portion delivers immediately and is then followed by the Square portion.

Nutrition - The process by which the body draws nutrients from food and uses them to make or mend cells.

Nutrition label - Used to determine the nutritional content of foods, as well as calories, carbohydrates, fats, proteins, vitamin and mineral content.

Nutritional content - The breakdown of nutrients found in a particular food.

O

Occlusive dressing - A bandage that seals a wound from air or bacteria.

P

Pattern A/B - Basal pattern that supports activities that are not a part of your day-to-day routine, but are normal in your lifestyle. Such activities could be a sport that you do once a

week or a change in your sleep pattern over the weekend, extended periods of higher or lower activity, or menses.

Post-prandial reading - Blood glucose reading taken after a meal

Predictive alerts - These alerts can be used so that you receive alarms before the sensor glucose measurements reach your glucose limits. This alert predicts when the sensor glucose measurements will reach your glucose limits by using information about your glucose limits and recent sensor glucose measurements.

Press - To push and release a button.

Protein - One of the three main sources of energy or calories in food. Proteins are made of amino acids, which are called the building blocks of the cells. The cells need proteins to grow and to mend themselves. Protein is found in many foods such as meat, fish, poultry, and eggs.

PSI - Pound-force per square inch.

Pump S/N - Pump S/N is the serial number of the pump currently in use.

R

Rate of Change alerts - These alerts can be used to set up the system so that you receive

alerts if the sensor glucose measurements change too rapidly for your needs. You can set the limit for these alerts so that you receive an alert whenever the sensor glucose measurements fall or rise faster than the limit.

Reservoir - The syringe that holds insulin.

Reset - To set again, or change the information (reset your glucose values).

Resume - Restarts basal delivery after the pump is suspended.

Rewind - The pump drive moves back to its starting position to prepare the pump for a new reservoir.

RF - Radio frequency

S

Scroll - Press the up or down arrow buttons to move through the screen text.

Select - Pressing the up or down arrow buttons to highlight a desired screen item.

Sen-serter - The Sen-serter is indicated as an aid for insertion of the Medtronic Diabetes glucose sensor.

Sensitive - Easily irritated (sensitive skin), or able to measure very slight differences (a sensitive instrument).

Sensor Age - Sensor age is the amount of time, in days and hours, since the sensor was first inserted.

Sensor features - Optional continuous glucose monitoring capabilities.

Sensor glucose (SG) - Glucose levels in the fluid under your skin measured by the sensor.

Sensor icon - An icon on your pump that tells you whether the Sensor feature is On or Off.

Set - To enter or establish a value for your pump (set your BG Reminder).

Set bolus - Selectable item available in the BOLUS MENU when the Bolus Wizard feature is inactive. One method of programming a bolus without the Bolus Wizard feature. (See *Manual bolus*.)

S/N# - Serial number.

Sodium - Is an essential element for your body. As table salt, it is used freely in foods as a flavor enhancer and a preservative. Sodium is listed on food labels and is expressed in milligrams.

Special mode - Operating mode that indicates one or more special functions is active or a condition exists that requires attention.

Square Wave® bolus - Immediate bolus delivered evenly over a

specified time period (30 minutes to 8 hours).

Square Wave portion - (Sq) The second part of a Dual Wave bolus. The Square Wave portion delivers evenly over a period of time after the NOW portion delivers.

Stabilize - To make steady or keep from changing (stabilize your blood glucose).

Standard pattern - Your normal basal that supports your usual day-to-day activity. When the Patterns feature is off, the pump uses your standard basal pattern.

Status - The condition of a part of your system (battery status, alarm/alert status).

Status screen - Displays the current operations of the pump, including active functions, the most recent basal and bolus deliveries, reservoir information, and battery condition.

Step - Measurement of insulin that you set and use for Easy Bolus and other bolus deliveries.

Sterile - Free from substances or germs that could cause infection.

Sterile gauze - A thin medical bandage with loose weave that is free from substances or germs that could cause infection.

Submerge - To place under water or cover completely with water.

Suspend - Function that stops all current insulin deliveries. The basal delivery is paused until restarted.

T

Temp - Temporary

Temp basal - (Tmp basal)
Temporary one-time basal insulin with a specified amount and duration. Used to support insulin needs for special activities or conditions that are not part of the normal daily routine.

Timeout - If you do not enter information or do something in a specific length of time, your pump automatically stops what it is doing and goes to the HOME screen.

Transmitter - The device that sends a signal to your pump with information from your sensor.

Transmtr Batt - The status of the transmitter battery. Possible values are **Good**, **Low**, or **Bad**.

Transmtr ID - The serial number of the transmitter currently in use.

Transmtr Ver - The software version of the transmitter currently in use.

W

Weak Signal - The pump displays an alert if it has not received data from the sensor for an amount of time that you set.

Index

Numerics

3-hour graph 145
 6-hour graph 146
 12-hour graph 146
 24-hour graph 147

A

accessories 2
 activity guard 3
 CareLink USB 3
 holster 3
 leather case 3
 meter 2
 pump clip 3
 remote control 2
 sensor 3
 transmitter 2
 active insulin 84
 Active insulin time
 setting 85
 adverse reactions 7
 alarm clock 162
 Alarm conditions 185
 alarm frequency 197
 ALARM HISTORY screen 153
 Alarm types
 A (alarm) 185
 Auto off 185
 Batt out limit 185

Bolus stopped 185
 Button error 185
 Check settings 175, 186
 E (error) 186
 Empty reservoir 186
 Failed batt test 186
 Finish loading 186
 Low suspend 188
 Max delivery 187
 Max fill reached 186
 Motor error 187
 No delivery 173, 187
 No reservoir 187
 Off No power 187
 Reset 187
 Weak battery 187
 Alarms
 responding to 184
 review 153
 alert frequency 197
 Alert Silence 149
 default settings 218
 options 149
 Alert Silence alert 179
 alert type
 Setting your 153
 alerts 179
 Alert Silence 179
 Cal Error 181, 182
 Change Sensor 181
 Charge Transmtr 181
 Fall Rate 183

for alert conditions 179
 low reservoir alert 179
 High Predicted 183
 High XX.X MMOL/L 182
 Lost Sensor 180
 Low Predicted 182
 Low Transmtr 180
 Low X.X MMOL/L 182
 Meter BG By 182
 Meter BG Now 182
 Missed Bolus 180
 Rise Rate 183
 sensor 179
 Sensor End 181
 Sensor Error 183
 Weak Signal 179
 Area Under the Curve 160
 areas
 inserting infusion set 67
 Attention mode 32
 AUC calculation 160
 AUC Limits 160
 default settings 218
 entering 160
 Auto Calibration 126
 default settings 217
 auto-off 154

B

backlight 31, 227
 bad transmitter battery 191

- BASAL MENU 34
- basal rate
 - daily 56
 - maximum 57
 - patterns 105
 - programming 55
 - stop delivery of 58
 - temp 108
- bathing 151
- battery, charger
 - low 192
- battery, pump 193, 203
 - install 26
 - removing 174
 - short life 174
- beep/vibrate, alert type 31
- BG, entering 50
- BG Reminder 40, 48, 88, 97, 101
- BG Targets 83
- BG units 82
 - default settings 217
- bleeding 5, 7
- Block feature 164
 - turn on 164
- BLOCK OPTION screen 164
- blood glucose
 - checking 138
 - entering incorrect number 181
 - rapid change 181
- blood glucose
 - measurements 50, 125
 - non-calibration meter, examples of 50
- blood glucose units
 - default settings 217
 - setting 118
- bolus
 - basic 39
 - details 45
 - history 44

- manual 96
 - maximum limit 47
 - stop delivery of 58
- bolus delivery
 - Bolus Wizard 87
 - Easy Bolus 101
- bolus delivery, methods for
 - normal
 - remote control 2
- BOLUS MENU 33
- bolus types
 - Dual Wave 95
 - Normal 40
 - Square Wave 95
- Bolus Wizard
 - about 75
 - entering meter BG 139
 - estimate details 46
 - HIGH BG warning 79
 - how it works 78
 - LOW BG warning 79
 - MAX BOLUS EXCEEDED warning 79
 - maximum delivery 79
 - on/off 80
 - review settings 85
 - settings 75, 76, 77, 78
 - Active insulin time 78
 - BG Target range 78
 - BG units 76
 - Carb ratios 76
 - Carb units 76
 - Insulin sensitivity 77
 - specifications 207
 - your BG reading 75
- bruising 7
- buttons, pump 25, 27

C

- Cal Error alert 181, 182

- Cal Reminder 182
 - default settings 217
 - setting 126
- CAL REMINDER screen 126
- Cal Repeat
 - default settings 217
 - setting 125
- CALIBRATE HISTORY 150
- calibrating
 - frequency 137
- calibrating sensor 137
 - frequency 137
- calibrating your system
 - Capture Event feature 49
- calibration
 - successful 150
 - values 150
 - when to 138
- Calibration history 150
- CAPTURE EVENT 34
- Capturing Events 49
- Carb ratios
 - exchanges 81
 - grams 81
- carb units 80
- CARB UNITS screen 52
- carbohydrate counting 19
- CareLink Personal
 - Software 49, 169
- CareLink USB 3
- Change Sensor alert 181
- Charge Transmtr alert 181
- charger 2
 - battery power 129
 - low AAA battery 192
 - transmitter 129
 - water-tightness 194
- CHECK BG message 49
- checking
 - blood glucose 138
- cleaning
 - transmitter 194
- cleaning the pump 194

- Clear settings 168
- connecting transmitter to
 - sensor 136
 - tester 189
- consumables 2
- continuous glucose graphs
 - reading 142
- continuous graph display 128
- continuous sensor glucose
 - levels 144
- correction bolus 40, 87, 88
- correction estimate 139, 207

D

- Daily totals 156
- DAILY TOTALS screen
 - Review 155
- default settings
 - Alert Silence 218
 - AUC Limits 218
 - Auto Calibration 217
 - BG Units 217
 - Cal Reminder 217
 - Cal Repeat 217
 - Glucose Alerts 215
 - Glucose Limits 215
 - Glucose Limits 1 216
 - Glucose Limits 2 through 8 216
 - High Glucose Limit 216
 - High Repeat 216
 - insulin delivery 213
 - Low Glucose Limit 216
 - Low Repeat 216
 - Predictive Alerts 216, 217
 - sensor 215
 - sensor features 215
 - sensor glucose graph
 - timeout 217
 - Weak Signal 217
- DELIVERY STATS screen 158, 159

- discarding
 - transmitter 194
- disconnecting
 - sensor from transmitter 151
 - tester from transmitter 190
 - transmitter from sensor 151
- downloading pump data 3
- Dual Wave bolus programming
 - with Bolus Wizard 99
 - without Bolus Wizard 96
- dual/square option
 - on-off 96

E

- Easy bolus 101
 - cancel 103
 - deliver 102
 - step value 102
- electromagnetic emissions 219
- Electromagnetic
 - environment 224
- electromagnetic immunity 220
- emergency kit 1
- emissions 219
- ENTER BG screen 50
- ENTER FOOD screen 52
- entering
 - carbohydrate information 52
 - exercise information 52
 - insulin injection
 - information 51
 - meter BG 137, 138, 139
 - alert 137
 - Bolus Wizard 139
 - frequency 137
 - manually 138
 - Other markers 53
- Entering meter BG
 - measurements 50
- ESC 3
 - cancels programming 34

- Examples of continuous sensor
 - glucose graphs 145
- exercising
 - precaution 7
- EXPRESS BOLUS button 28

F

- Fall Rate alert 183
- Fall Rate of Change
 - setting alert frequency 122
- FILL TUBING screen 66
- filling
 - cannula 70
 - reservoir 61
 - tubing 66
- FIND LOST SENSOR 191
- food bolus 40, 76, 78, 87, 88
- food entry 75, 88
- food estimate 207, 208
- food label 20
- frequency
 - alarm 197
 - alert 197

G

- Glucose Alerts
 - default settings 215
 - turning on 118
- GLUCOSE ALERTS screen 118
- Glucose Limits
 - default settings 215
- Glucose Limits 1
 - default settings 216
- Glucose Limits 2 through 8
 - default settings 216

- glucose monitoring settings
 - entering
 - 118, 119, 120, 122, 124, 126, 160
 - selecting BG units 119
 - setting AUC limits 160
 - setting Glucose Limits 118
 - setting Glucose Limits, guidelines 120
 - setting Glucose Limits, setting start time 120
 - setting Glucose Limits, switching limits off 120
 - setting Predictive glucose alert 122
 - setting Rate of Change alerts 124
 - setting the Cal Reminder 126
 - turning on Glucose Alerts 118
- graphs 144
 - continuous glucose 142
 - continuous sensor glucose levels 144
 - differentiating between current and earlier events in 143
 - opening 143, 144
- green light
 - transmitter 190

H

- high BG levels
 - Bolus Wizard and 79
- HIGH BG warning 79
- High Glucose Limit
 - default settings 216
 - setting 118, 182
- High Predicted alert 183
- High Predictive alert 123

- High Repeat
 - default settings 216
 - setting 121
- High Repeat alert
 - setting 121
- High sensor glucose alert 182
- history
 - alarms 153
 - bolus delivery 44
 - calibration 150
 - delivered insulin 70
 - sensor alert 150
 - user settings 168
- HOME screen 29

I

- infection 5, 7
- infusion set 2, 61
 - changing your 64
 - inserting 67, 68
 - best sites 68
 - rotating sites 68
- initialization 137, 142
- insulin delivery
 - stopped 125
- Insulin sensitivity 82
- interference
 - transmitter 5
- irritation
 - occlusive dressing 5, 7
 - skin 5
- ISIG 190

L

- Language setting 39
- list
 - calibration values 150
 - sensor alert 150
- Lock keypad feature 165

- Lost Sensor alert 180
- low battery
 - affects meter option 86
 - affects remote control 163
 - affects selftest 165
 - affects your alert type 153
- low BG levels
 - Bolus Wizard and 79
- LOW BG warning 79
- Low Glucose
 - setting alert frequency 122
- Low Glucose Limit
 - default settings 216
 - setting 118, 182
- Low Glucose Suspend
 - setting 125
- Low Predicted alert 182
- Low Predictive
 - setting alert frequency 122
- Low Repeat
 - default settings 216
 - setting 122
- low reservoir alert 179
- Low resv warning feature 154
- Low sensor glucose alert 182
- low transmitter battery 180
- Low Transmtr alert 180
- LR-03 battery
 - low 192

M

- MAIN MENU 33
- manually entering meter BG 138
- Max basal rate 57
- MAX BOLUS EXCEEDED
 - warning 79
- Maximum bolus limit 47
- menus 33
- meter 2, 75
 - entering BG 137, 138
- Meter BG By alert 182

METER BG NOW alert 125
 Meter BG Now alert 182
 meter IDs
 add, delete, review 87
 meter option 86
 rules 86
 METER OPTION screen 87
 METER STATS screen 159, 160
 MiniLink transmitter 2
 Missed Bolus alert 180
 Missed Bolus Reminder 104
 modes
 Attention 32
 Normal 32
 Special 32
 MWT1 2
 MWT1 technology 75, 138

N

no delivery 173
 Normal bolus
 with Bolus Wizard 87
 without Bolus Wizard 40
 Normal mode 32
 Now 100
 Now portion
 of the Dual Wave
 bolus 97, 100
 nutrition 19

O

occlusive dressing
 irritation 5, 7
 reaction 5, 7
 Other markers 53

P

pain 7
 Patterns
 basal 105
 on/off 106
 pattern A, B 105
 Program a basal pattern 106
 Select a basal pattern 107
 standard basal pattern 105
 Personal reminders 162
 personal settings 75
 practicing 65
 precaution
 sensor 7
 Predictive Alerts 216, 217
 PREDICTIVE ALERTS LOW/
 HIGH screen 123
 programming sensor 117
 pump
 battery compartment 25
 buttons 25, 27
 if you remove your 35
 screen 25, 29
 using in water 177
 warranty vi

R

Rate of Change alerts
 sensitivity and alarm
 frequency 124
 reaction to occlusive
 dressing 5, 7
 reading the continuous glucose
 graphs 142
 RECONNECT OLD SENSOR 191
 Reminder, alarm clock 162
 remote control 2
 remote control IDs
 add, delete, review 163

Remote control option
 turn on 163
 removing
 sensor 151
 removing your pump 35
 reservoir
 filling the 61
 inserting in your pump 65
 low reservoir warning 154
 plunger rod 62
 removing the 64
 transfer guard 62
 tubing connector 25
 RESERVOIR + SET 34
 reservoirs 2
 Restore settings 167
 REVIEW SETTINGS screen 128
 reviewing
 sensor settings 128
 rewinding your pump 64
 during practice 65
 RF communications 224
 RF features
 meter 86
 other devices 163
 remote control 162
 use on aircraft 7
 Rise Rate alert 183

S

Save settings 167
 screens
 DELIVERY STATS 158, 159
 METER STATS 159, 160
 SENSOR STATS 159
 scroll bar 31
 Scroll rate
 setting 47
 Selftest 165, 176

- ul style="list-style-type: none;">
- sensor
 - alerts 150, 179
 - history 150
 - list 150
 - calibrating 137, 182
 - connecting with
 - transmitter 136
 - default settings 215
 - disconnected 180
 - disconnecting from
 - transmitter 151
 - information 141
 - ISIG 190
 - life 181
 - precaution 7
 - programming 117
 - reconnecting 191
 - removing 151
 - setting up 117
 - signals 183
 - too high 183
 - too low 183
 - stabilizing 181
 - starting 132
 - status screens 141
 - Transmitter ID 126
 - using in water 151
 - warnings 5
- SENSOR ALERT HISTORY 150
- sensor alerts
 - Alert Silence 179
 - Cal Error 181, 182
 - Change Sensor 181
 - Charge Transmtr 181
 - Fall Rate 183
 - High Predicted 183
 - High XX.X MMOL/L 182
 - Lost Sensor 180
 - Low Predicted 182
 - Low Transmtr 180
 - Low X.X MMOL/L 182
 - Meter BG By 182
 - Meter BG Now 182
 - Missed Bolus 180
 - Rise Rate 183
 - Sensor End 181
 - Sensor Error 183
 - Weak signal 179
- Sensor Demo 129
- SENSOR DEMO ON/OFF
 - screen 129
- Sensor End alert 181
- Sensor Error alert 183
- Sensor function 33
- sensor glucose graph
 - timeout 128
 - default settings 217
- sensor glucose graphs
 - continuous display 128
 - WARM UP 142
- sensor glucose levels
 - graphs 144
- sensor graph
 - example 129
- SENSOR ON/OFF screen 117
- SENSOR STATS screen 159
- SENSOR STATUS 141
- SET AUC LIMITS screen 160
- SET CAL REMINDER screen 126
- SET CAL REPEAT screen 125
- SET FALL RATE LIMIT screen 124
- SET GLUCOSE LIMITS 1
 - screen 120
- SET GRAPH TIMEOUT screen 128
- SET HI ALERT REPEAT
 - screen 121
- SET LO ALERT REPEAT
 - screen 122
- SET RISE RATE LIMIT screen 124
- SET TIME SENSITIVITY LOW/
 - HIGH screen 123
- SET TRANSMITTER ID screen 127
- SET WEAK SIGNAL screen 127
- setting
 - Cal Reminder 126
 - Cal Repeat 125
 - High Repeat 121
 - Low Glucose Suspend 125
 - Low Repeat 122
 - sensor glucose graph
 - timeout 128
- setting up sensor 117
- settings
 - entering
 - 118, 120, 122, 124, 126, 160
 - glucose monitoring
 - settings, setting AUC limits 160
 - glucose monitoring
 - settings, setting Glucose Limits, guidelines 120
 - glucose monitoring
 - settings, setting Glucose Limits, setting start time 120
 - glucose monitoring
 - settings, setting Predictive glucose alert 122
 - glucose monitoring
 - settings, setting Rate of Change alerts 124
 - glucose monitoring
 - settings, setting the Cal Reminder 126
 - glucose monitoring
 - settings, turning on Glucose Alerts 118
 - glucose
 - monitoring 118, 120, 122, 124, 126, 160
 - setting AUC limits 160
 - setting Glucose Limits 118
 - setting Glucose Limits, guidelines 120
 - setting Glucose Limits, setting start time 120

- setting Predictive
 - glucose alert 122
- setting Rate of Change
 - alerts 124
- setting the Cal
 - Reminder 126
- turning on Glucose
 - Alerts 118
- sensor glucose graph
 - timeout 128
- settings, reviewing 128
- showering 151
- silent alerts 149
- skin
 - bleeding 5, 7
 - bruising 7
 - infection 5, 7
 - redness 7
 - swelling 5, 7
- Special mode 32
- Square portion
 - of the Dual Wave bolus 100
- Square Wave bolus programming
 - with Bolus Wizard 99
 - without Bolus Wizard 96
- starting sensor
 - procedure 137
 - steps before 132
- STATUS screen 34
- Step value setup 102
- supplies, ordering 3
- Suspend function 33, 58
- swelling 5, 7
- swimming 151
- system initialization 137

T

- Temp basal
 - canceling 113
 - delivering 112

- rates 108
- selecting type 112
- verify delivery 113
- temp basal types
 - Insulin rate 109
 - Percent of basal 110
- tenderness 7
- terms and symbols 3
- test plug. See tester 189
- tester 2
 - connecting to
 - transmitter 189
 - disconnecting from
 - transmitter 190
 - transmitter 129
 - water-tightness 194
- therapy follow-up
 - Recommended 171
- time and date, setting 37
- transmitter
 - bad battery 192
 - battery 129
 - charger 129
 - cleaning 194
 - connecting with sensor 136
 - depleted battery 181
 - discarding 194
 - disconnected 180
 - disconnecting from
 - sensor 151
 - end of life 192
 - green light 190
 - interference 5
 - low battery 180
 - tester 2, 129
 - troubleshooting 129
 - using in water 151
 - warning 5
- Transmitter ID
 - setting 126
- troubleshooting 191
- tubing
 - filling 66

U

- user safety
 - adverse reactions 7
- User settings 166
- User Settings screen 177
- using the Capture Event
 - feature 49
- Utilities 153
- UTILITIES MENU 34

V

- viewing
 - calibration values 150
 - sensor alert history 150

W

- WARM UP 142
- warnings
 - sensor 5
 - transmitter 5
- warranty vi
- water, using your system in 151
- water-tightness 151
 - charger 194
 - tester 194
- Weak Signal 179
 - default settings 217
- Weak Signal alert 179
- Weak Signal feature 127
 - description 127

