




START HERE

1 


Wash your hands. Press **O**. Select **Options**.

2 


Select **Reservoir & Tubing**.

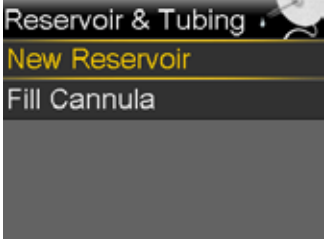3 


Select **New Reservoir**.

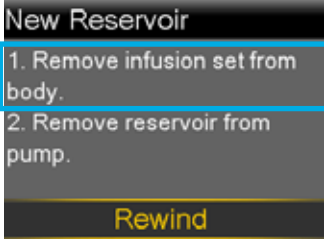4 


Remove the infusion set you have been using by loosening the adhesive and pulling away from body.

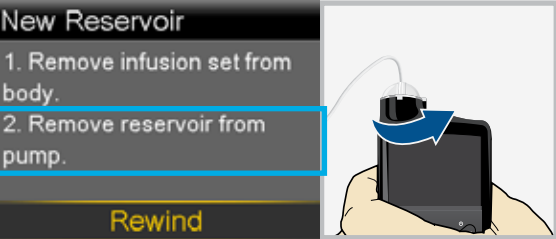
5 

Remove the used reservoir from the pump.

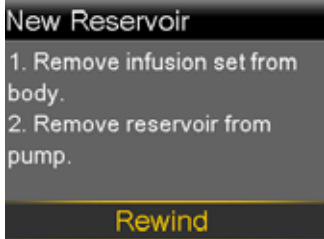6 


Select **Rewind**.



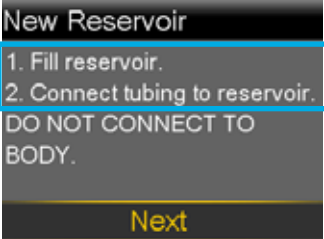


FILL RESERVOIR & CONNECT TO THE INFUSION SET TUBING


Follow the next steps to fill reservoir with insulin and connect to the infusion set tubing.

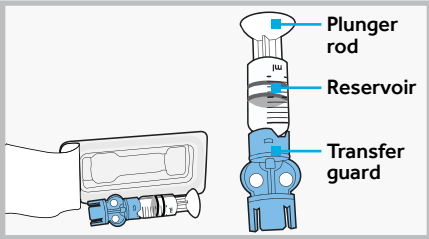1 


Remove from package. Make sure insulin vial is at room temperature to reduce the risk of air bubbles.

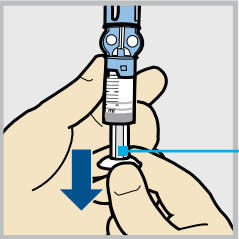2 

Pull plunger down to the amount that you plan to fill with insulin.

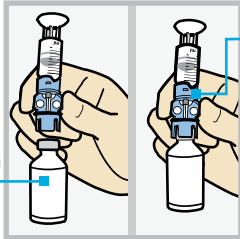
3 


Wipe vial with alcohol. Place vial on table. Firmly press the blue transfer guard onto vial.

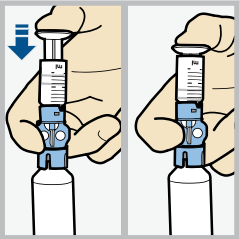
4 


Push and hold plunger down.

WARNING If insulin or any liquid gets inside the tubing connector, it can temporarily block the vents that allow the pump to properly fill the infusion set. This may result in the delivery of too little or too much insulin, which could cause hyperglycaemia or hypoglycaemia.

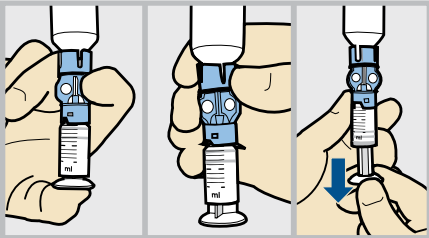- 5** 

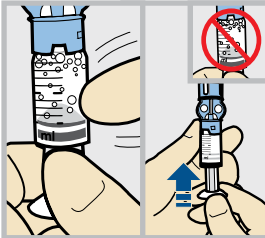With your thumb still on the plunger, flip over so vial is on top. Release thumb and pull plunger down to fill with insulin.
- 6** 

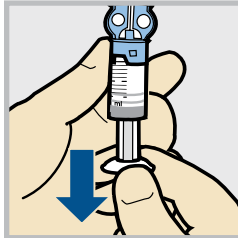
Tap the reservoir to move air bubbles to top of reservoir. Push plunger up to move air into vial.
- 7** 

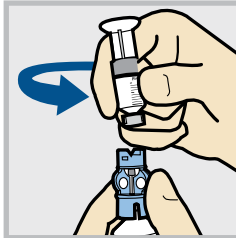
If needed, pull plunger back down to amount of insulin needed for 2-3 days.
- 8** 


To avoid getting insulin on the top of the reservoir, turn vial over so it is upright. Hold transfer guard and turn reservoir counter-clockwise and remove from transfer guard.


CONNECT RESERVOIR TO INFUSION SET

You will place the tubing connector onto the end of the infusion set to the filled reservoir.

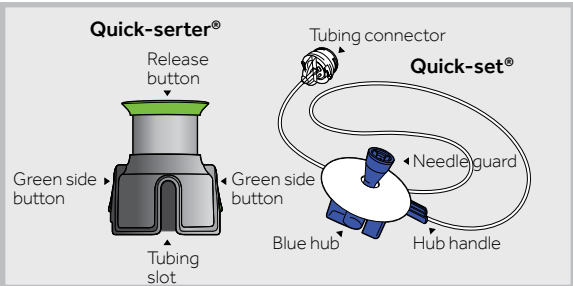- 1** 

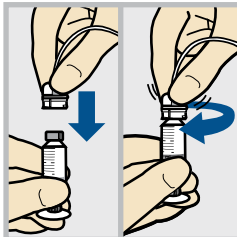Remove infusion set from package. Remove the paper that holds the tubing together.
- 2** 

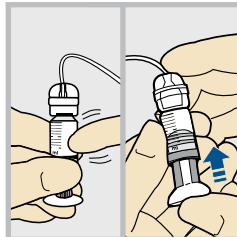Gently push the tubing connector onto reservoir. Turn clockwise until locked. You will hear a click.
- 3** 

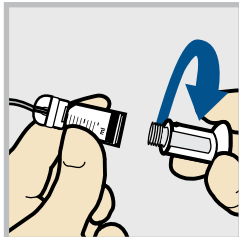

If you see air bubbles, tap reservoir to move them to top. Push plunger just a bit to move them into tubing.
- 4** 

Twist plunger counter-clockwise to loosen and remove.

THE BACKLIGHT MAY HAVE TURNED OFF

Press any button to turn the screen back on.



Press  to open the menu. If the pump is locked, you will need to unlock the pump after pressing .



Select **Load Reservoir** from the menu.



Select **Next**.

Continued on next page

PLACE RESERVOIR INTO PUMP


Now place the filled reservoir into the reservoir compartment of the pump.

1 New Reservoir


3. Place reservoir into pump and lock.
DO NOT CONNECT TO BODY.

Next

Place reservoir into pump.

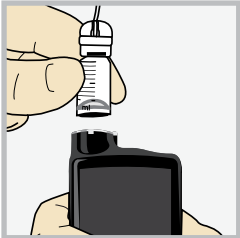

2

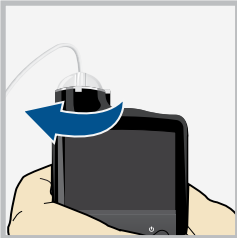

Turn clockwise until you feel reservoir lock into place.

3 New Reservoir

3. Place reservoir into pump and lock.
DO NOT CONNECT TO BODY.

Next

Select **Next**.

LOAD RESERVOIR AND FILL TUBING

Follow these steps to load the reservoir and fill the tubing.

1 Load Reservoir

Select Load and hold until complete.
DO NOT CONNECT TO BODY.

Load Next

Select **Load** and keep holding ○.

2 Load Reservoir

Complete

DO NOT CONNECT TO BODY.

Load Next

When you see this screen, select **Next**.


3 Fill Tubing

DO NOT CONNECT TO BODY.
Hold Fill until drops appear.
Then select Next.

0.0 u

Fill Next

Select **Fill** and keep holding ○ until you see drops at the end of tubing, then let go.

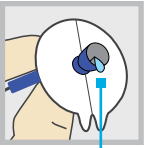

Drops at end of tubing

4 Fill Tubing

DO NOT CONNECT TO BODY.
Hold Fill until drops appear.
Then select Next.

11.3 u

Fill Next

After you see drops, press > and select **Next**.

INSERT INFUSION SET


Next, follow the steps to insert the infusion set into your body.

Fill Cannula?

1. Insert infusion set into body.
2. Select Fill to fill cannula or Done if not needed.


Fill Done

1

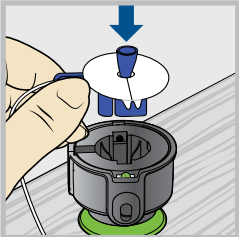


Place the MiniMed Quick-setter insertion device onto a sturdy flat surface with handle facing down. Place **blue** hub intoserter, placing the handle in the tubing slot.

2

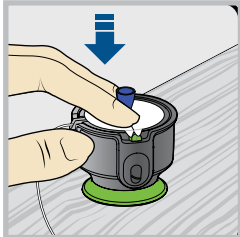


Use two fingers to seat the infusion set inside the inserter securely and gently push down.

3

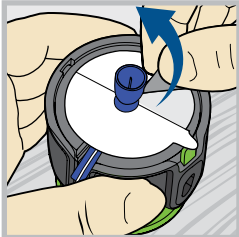


Peel the paper from the adhesive on both sides of the needle guard.


Continued on next page

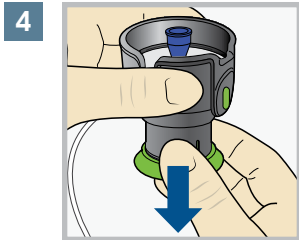


4 Pull green handle down until it clicks. The two **green** side buttons will also click into place. Do not press the **green** side buttons at this time.

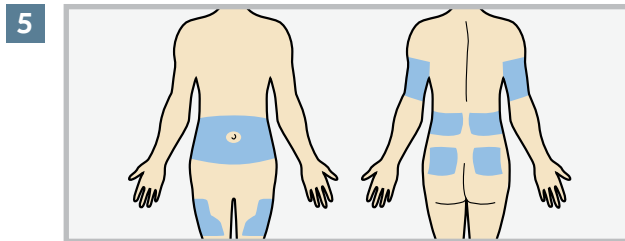


5 The best areas of the body to insert the infusion set are shown in the shaded areas of this drawing. Choose areas that are away from the natural bend of your waist or belt line of your pants. Avoid areas of hard or scar tissue. Consult with your healthcare provider on appropriate areas to insert your infusion set. Wipe with alcohol or other antiseptic.

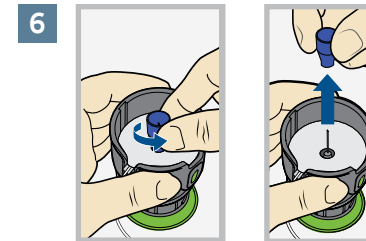


6 Turn to loosen needle guard and pull.

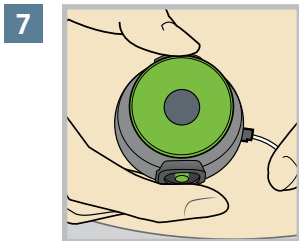


7 Hold serter against cleaned site.

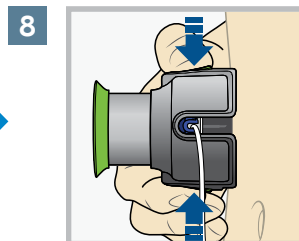


8 Press the two **green** side buttons at the same time. If they are not pressed at the same time the infusion set will not insert properly.

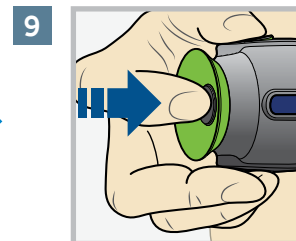


9 Press down the release button of the serter to release the infusion set.

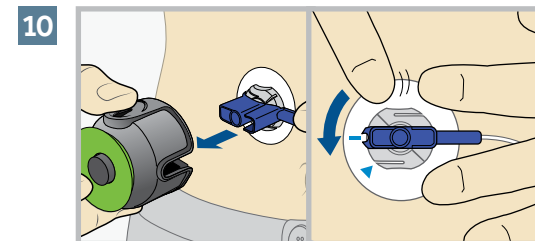

10 Pull serter away from body. Press adhesive against skin.

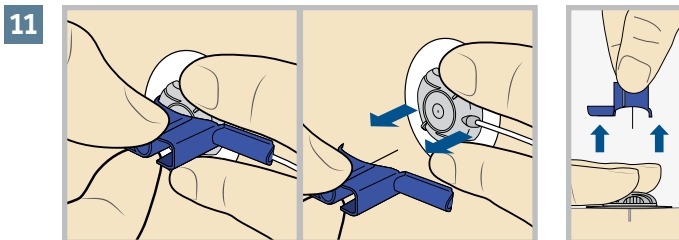

11 Hold infusion set. Pull **blue** handle straight out to remove needle.

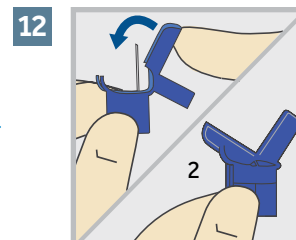


12 Fold **blue** handle until locked.

Dispose blue handle into a sharps container.

FILL CANNULA

You will now fill the cannula, the little tube under your skin, with insulin.

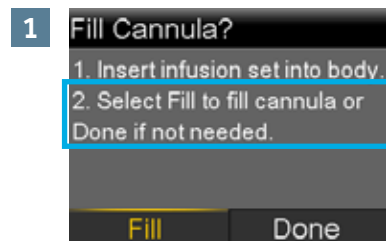

Select **Fill**.



Select **Fill amount** and enter :

- 0.300 if using 6mm cannula
- 0.500 if using 9mm cannula

Then press **0**.



NOTE: Your pump will remember the **Fill amount** that you used last. Always verify that the Fill amount is correct.

- If it is correct, press **✓** to **Fill Now** and press **0**.
- If it is incorrect, press **0**. Change to correct amount and press **0**. Select **Fill Now**.



Select **Fill Now**.



The Home screen displays the insulin as it fills the cannula.



NOTE: Select **Stop Filling** if you need to stop, for example, if you notice the **Total** amount is incorrect. This should rarely happen if you have verified the **Fill amount** on the previous screen.

Your infusion set change is now complete!